BodiLink® on Permobil Corpus

OPERATION MANUAL

BodiLink® Head Support
BodiLink® Lateral Trunk Support
BodiLink® Lateral Pelvic/Thigh Support
on the Permobil Corpus* Power Wheelchair

**Compatible models only.*

Thank you for choosing BodiLink®

CONTACT INFORMATION

The Comfort Companies, LLC
1960 South Calhoun Road
New Berlin, Wisconsin 53151 USA

Phone:
U.S.A.: 800-736-0925
Outside the USA: 262-432-0830

Fax: 406-522-8563

Customer Support:
orders.comfort@permobil.com

permobilus.com
comfortcompany.com

**WARNING**

Indicates that not following the specified procedure could lead to potentially hazardous conditions resulting in serious injury.

**CAUTION**

Indicates that not following the specified procedure could lead to potentially hazardous conditions resulting in minor to moderate injury or damage to the equipment or other property.

PLEASE NOTE

Supplier: These instructions must be given to the user of this product.

Operator (Individual or Caregiver): Before using this product, read the instructions and save for future reference.

Products may be covered by one or more U.S. and foreign patents and trademarks, including Permobil®, BodiLink® and Stretch-Air®.

Product information is changed as needed; current product information is available at comfortcompany.com.

The Comfort Companies, LLC is a part of Permobil Seating and Positioning.

© 2019 Permobil

REV 2019-11-20

TABLE OF CONTENTS

INTENDED USE - INTENDED USER	2
IMPORTANT SAFETY INFORMATION	3
PARTS DETAIL	4
BODILINK HEAD SUPPORT	HS
Installation - Permobil "A" Mounting Block	5
Installation - Permobil "B" Mounting Block	6
Installation - Corpus Support Pad	7
Hardware Adjustment	
BODILINK LATERAL TRUNK SUPPORT	LTS
Support Pad Attachment	8
Installation - Track Mount	9
Adjustment and Use	10
BODILINK LATERAL PELVIC/THIGH SUPPORT	LPTS
Installation - Support Pad to Mounting Hardware	11
Installation - Power Mount Hardware	11
Adjustment - Support Pad	12
Adjustment - Power Mount Hardware	13
Use - Hardware	14
Use - Hardware Release Lever	14
CLEANING	15
MAINTENANCE AND WARRANTY	16
Español	17
Français	33
Deutsch	49
Português	65
Italiano	81
Nederlands	97
Dansk	113
Svenska	129
Norsk	145
Suomi	161
Čeština	177
Slovenščina	193

INTENDED USE - INTENDED USER**BodiLink Head Support - Intended Use****HS**

The BodiLink® Head Support is a covered head support pad with soft, durable foam layered over high-resiliency foam, mounted on adjustable hardware. The BodiLink Head Support is intended to support and conform to the user's head and to provide comfort, positioning, and pressure redistribution. The product must be properly sized to the user and must be used with the provided cover.

BodiLink Lateral Trunk Support - Intended Use**LTS**

The BodiLink® Lateral Trunk Support is a support pad with compatible attachment hardware. The BodiLink Lateral Trunk Support is intended to provide support, positioning, and comfort.

The BodiLink Premium Lateral Trunk Support Pad is a cushioning pad molded over a moldable or rigid metal insert. It may be used with or without the cover.

BodiLink Lateral Pelvic/Thigh Support - Intended Use**LPTS**

The BodiLink® Lateral Pelvic/Thigh Support is a support pad with compatible attachment hardware. The BodiLink® Lateral Pelvic/Thigh Support is intended to provide support, positioning, and comfort.

The BodiLink® Lateral Pelvic/Thigh Support Premium Pads are constructed of molded foam with a track for adjustability. A cover is always required for use.

BodiLink Supports - Intended User**HS****LTS****LPTS**

BodiLink Supports are intended for those with positioning and support needs as deemed necessary by a qualified healthcare professional. BodiLink Supports can be positioned to accommodate user rehabilitation needs and goals while allowing for function and comfort. The sizes available and range of adjustability fit pediatric to full grown users and should be matched by the healthcare provider to suit the user needs. Adjustments or replacement components may be necessary as changes in condition occur or to accommodate developments such as growth.

BodiLink Supports - Compatibility

BodiLink Head Support and Lateral Pelvic/Thigh Support: Power Wheelchairs: Corpus, Corpus VS, Corpus HD

BodiLink Lateral Trunk Support: Power Wheelchairs: Corpus, Corpus VS

A clinician who is experienced in seating and positioning should determine whether the product is appropriate for the individual's particular seating needs. The product's intended use statement alone should not be relied upon to make that determination.

The medical products featured in this document are intended to be only one part of an overall care regimen that includes all seating and mobility equipment and therapeutic interventions. A clinician who is experienced in seating and positioning should identify this care regimen after evaluating the individual's physical needs and overall medical condition.

A clinician should also evaluate for visual, reading, and cognitive impairment to determine the need for caregiver assistance or the need for other assistive technology, such as large-print instructions, to ensure proper product use.

IMPORTANT SAFETY INFORMATION**WARNINGS: Follow These Instructions**

- User and Caregiver: Follow all instructions and safety information provided with the product and provided by the manufacturer of your wheelchair or other accessories. The product should be installed and adjusted by a qualified and competent medical equipment provider or licensed healthcare professional.
- Permobil cannot be held responsible for damages or injury due to misuse of the product or failure to follow instructions provided in this manual. Contact Customer Support if you require further assistance in understanding cautions, warnings, or instructions for using or maintaining your product.

**WARNINGS: Use As Intended:**

- This product is designed for the intended user and intended use only. Any modifications could potentially endanger the user and cause discomfort or damage.
- Never use the equipment to maneuver a wheelchair or anything that it is attached to. It is never to be used as a handle, lift, guide, or support, except as described in this manual.
- DO NOT use any BodiLink Supports to transfer in or out of the wheelchair. BodiLink Supports are not designed to support the full weight of the user.
- During the first few hours of use, a healthcare professional should be available to observe skin condition and assess the likelihood of skin breakdown. During regular use, periodically check for skin discoloration and/or irritation. If there is any discoloration to skin/soft tissue, STOP USE immediately. If the discoloration does not disappear within 30 minutes of disuse, immediately consult a healthcare professional.
- DO NOT place anything between the user and the product. Additional materials may reduce the effectiveness of the product and increase risk factors for skin breakdown and instability.

**WARNINGS: Installation, Interference and Use**

- Never use the equipment unless all components are fully functional. Refer to the "Warranty" section in this manual.
- ATTENTION! Before use: Check regularly to make sure the equipment is securely mounted to the wheelchair system and is properly adjusted, and check all hardware for wear, breakage, and tightness. Re-tighten any loose screws/bolts, referring to the instructions in this manual. DO NOT continue to use the equipment if any part cannot be tightened or is broken.
- Unless otherwise specified, fasteners should be tightened to 9.6 Nm (85 lb-in).
- ATTENTION! Always check that the gear teeth in each joint are aligned for full tooth engagement before final tightening. Failing to do so could cause wear on the teeth, leading to parts slipping out of adjustment or hardware failure. (GT Hardware only)
- Ensure that the product does not interfere with any other wheelchair components. Interference could jeopardize effectiveness and user safety and cause damage. For the Permobil Corpus wheelchair: Use of tilt or recline may result in equipment interference. After installation, slowly articulate the full range of power recline / tilt, watching closely to make sure there is no interference.
- DO NOT modify the product. Doing so will void its warranty and may lead to personal injury or alter the effectiveness of the product by increasing the risk factors for skin breakdown and/or instability.

- Installing alternate parts to your wheelchair may change the structure and function of the equipment. To ensure stability during use, assess the need for additional safety features on the wheelchair, such as anti-tip bars or other available options.
- Adding additional links to the hardware may reduce the overall strength.
- Use caution when installing, adjusting, or removing hardware, to avoid pinching or trapping fingers in openings.
- Do not expose the product or components to high heat, open flames or hot ashes. Testing or certification claims, including for flammability, may no longer apply to this device when it is combined with other products or materials. Check testing and certification claims for all products used in combination with this device.
- Metal components may become hot if exposed to high temperatures or cold if exposed to low temperatures. Take appropriate precautions, especially when the product will contact unprotected skin.
- Avoid sharp objects.
- Protect the product from blunt force shocks that may cause damage and/or breakage.

**WARNINGS: Motor Vehicle Transportation:**

- Failure to pay attention to these warnings could result in severe injury to the individual in the wheelchair or to others.
- Whenever possible, transfer out of your wheelchair installed with BodiLink Supports and into a manufacturer-installed vehicle seat, and use the vehicle's crash-tested occupant restraint system.
- Children weighing less than 48½ lb. (22 kg) must be transferred into an appropriate child-restraint system intended for use in motor vehicles.
- ALWAYS store BodiLink Supports safely during transport to avoid damage.
- If a BodiLink Support has been involved in an accident during transport, see the "Inspection" section in this manual.
- The motor vehicle safety information in this operation manual applies to BodiLink Supports a) only when correctly and securely installed, as instructed in this manual; b) only with a wheelchair that complies with the performance requirements of ISO 7176-19 and that is installed, used, and maintained according to the wheelchair manufacturer's instructions.
- The individual should be seated in an upright position and the wheelchair should face forward during transportation in a motor vehicle. BodiLink Supports have only been crash tested on a wheelchair in the forward-facing position, per the test information provided in this operation manual.
- The wheelchair with a BodiLink Support must be used with an effective wheelchair securement system and a properly positioned, crash-tested pelvic- and shoulder-belt restraint, or a Wheelchair Tiedown and Occupant Restraint System (WTORS), following the manufacturer's instructions.

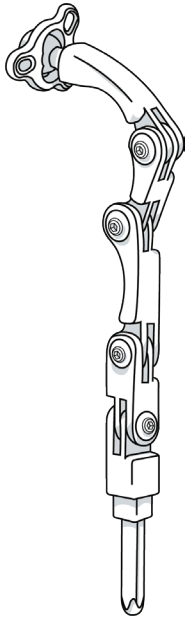
The Comfort Company seating system with a back support and accessories, including BodiLink Supports, installed on a surrogate wheelchair frame, has been dynamically tested by a third-party testing facility for use in a motor vehicle, and meets all applicable performance criteria in Section 5.1 of ANSI/RESNA WC-4:2012 Section 20 "Wheelchair Seating Systems for Use in Motor Vehicles" and in 5.1 of ISO 16840-4 when tested with a surrogate seat.

PARTS DETAIL

BodiLink Head Support (HS) Hardware


HS

Removable Hardware - Max 8


Adapter Plate

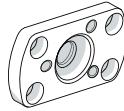For Corpus Headrest Support Pad


Tool-Less Wing Lock

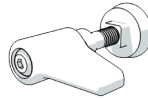Available through Special Adaptation


Mounting Blocks

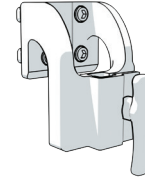Permobil "A"

(Corpus 23" - 28" H
Corpus HD 24"-29" H)

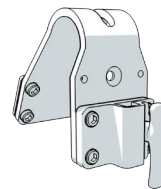
Permobil "B"

(Corpus 20" H
Corpus HD 21" H)


BodiLink Lateral Trunk Support (LTS) Hardware and Support Pads

LTS

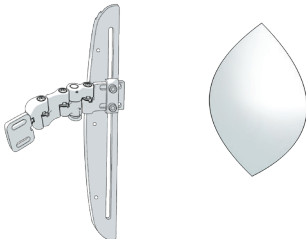Permobil Power Mount Hardware

Support Pads

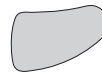L x D: 3.5" x 5.5" (8.9cm x 14cm)

L x D: 4" x 6.5" (10.2cm x 16.5cm)

BodiLink Lateral Pelvic/Thigh Support (LPTS) Hardware and Support Pads

LPTS

Power Mount Hardware - Medium

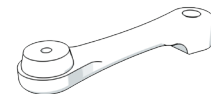
Support Pads

L x D: 3.5" x 4.5"
(8.9cm x 11.4cm)L x D: 3.5" x 8"
(8.9cm x 20.3cm)L x D: 3.5" x 12"
(8.9cm x 30.5cm)L x D: 5.5" x 6"
(14 cm x 15.2cm)

Long Extension Arm (4")

Available through Special Adaptation

HS BodiLink Head Support

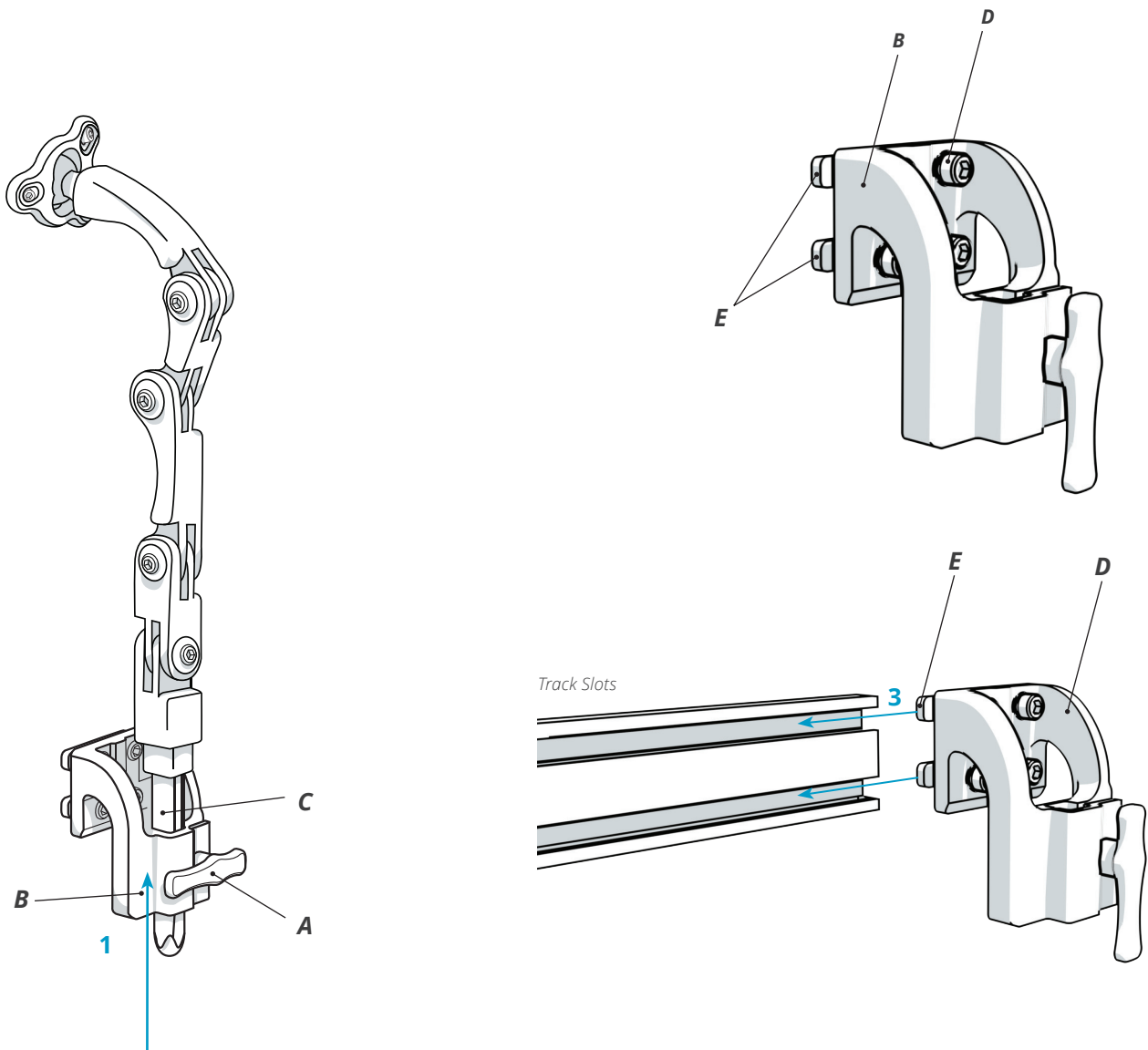
INSTALLATION - PERMOBIL "A" MOUNTING BLOCK

Corpus (23"-28" height) and Corpus HD (24"-29" height) with UniTrack System

- 1 If the head support hardware stem (C) is inserted in the mounting block (B), turn the knob (A) counterclockwise to release it. Remove the hardware stem; set aside.
- 2 Using the 5mm hex key, loosen but do not remove the screws (D) to create a gap between the double nuts (E) and the mounting block (B).
- 3 Slide the double nuts (E) into the UniTrack System slots and move into position.
- 4 Tighten the screws (D) with the 5mm hex key to secure.
- 5 Insert the hardware stem (C) into the mounting block (B) and turn the knob (A) clockwise to secure.

! IMPORTANT! Make sure that at least ½" (1.27cm) of the hardware stem is completely through the bottom of the mounting block (B) at all times. See the warning in the "Important Safety Information" section in this manual.

Note: See SUPPORT PAD installation, positioning and adjustment instructions in this manual to make final adjustments.


INSTALLATION - PERMOBIL "B" MOUNTING BLOCK**Corpus (20" height) and Corpus HD (21" height)**



- 1 If the head support hardware stem is inserted in the mounting block (B), turn the knob (A) counterclockwise to release it. Remove the hardware stem; set aside.
- 2 Using the 5mm hex key, loosen the screws (C) to release the washers (D), and the double nuts (F) from the mounting bracket (E).
- 3 Move the back support cover (not shown) allowing access from both sides to the head support holes on the back support shell (G).
- 4 Align the mounting bracket (E) with the head support holes on the back support shell (G) and the double nuts (F).
- 5 Replace the screws (C) and the washers (D) and tighten.
- 6 Replace and secure the back support cover.
- 7 Insert the hardware stem into the mounting block (B) and turn the knob (A) clockwise to secure.

! IMPORTANT! Make sure that at least ½" of the hardware stem is completely through the bottom of mounting block (B) at all times. See the warning in the "Important Safety Information" section in this manual.

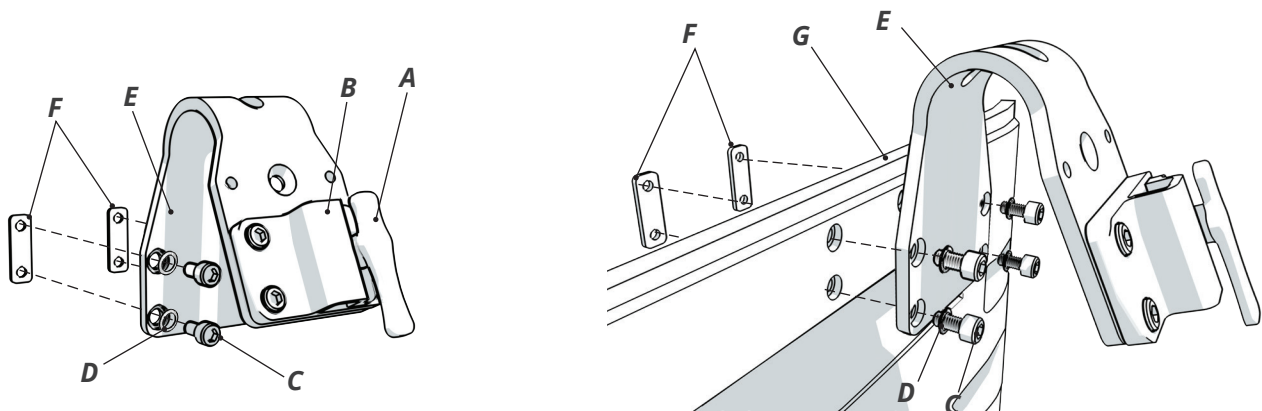Note: See SUPPORT PAD installation, positioning and adjustment instructions in this manual to make final adjustments.

**Alternate options when using the Permobil "B" Mounting Block.**

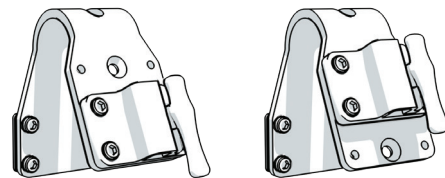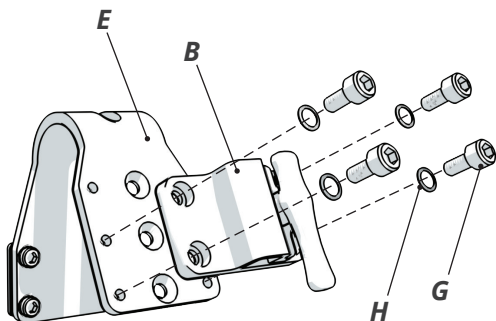- 1 Remove the screws (G) and the washers (H) on the mounting block (B).
- 2 Align the mounting block (B) with the alternate holes on the mounting bracket (E).
- 3 Replace the screws (G) and the washers (H) and tighten.

Standard and Alternate Mounting Block Placement on the Mounting Bracket

INSTALLATION - CORPUS SUPPORT PAD AND HARDWARE ADJUSTMENT

Attaching BodiLink Head Support Hardware to the Corpus Headrest Support Pad

Package contents: Corpus Adapter Plate; (4) M6 x 1 x 10mm socket head cap screws; (3) M6 x 1 x 20mm socket head cap screws; (3) safety washers / Not included: 5mm hex key / Material: aluminum

WARNING: The Corpus Adapter Plate is intended **ONLY** to be used with the Corpus Headrest Support Pad to attach to BodiLink Head Support Hardware.

1 Align the holes in the Adapter Plate (A) and Corpus Headrest Support Pad (B). Insert the four screws (C) and tighten.

2 Align the ball link and ball ring (D) with the ball receiver (E). Loosely engage the three screws (F) and the safety washers (G) with the three holes in the Adapter Plate (A).

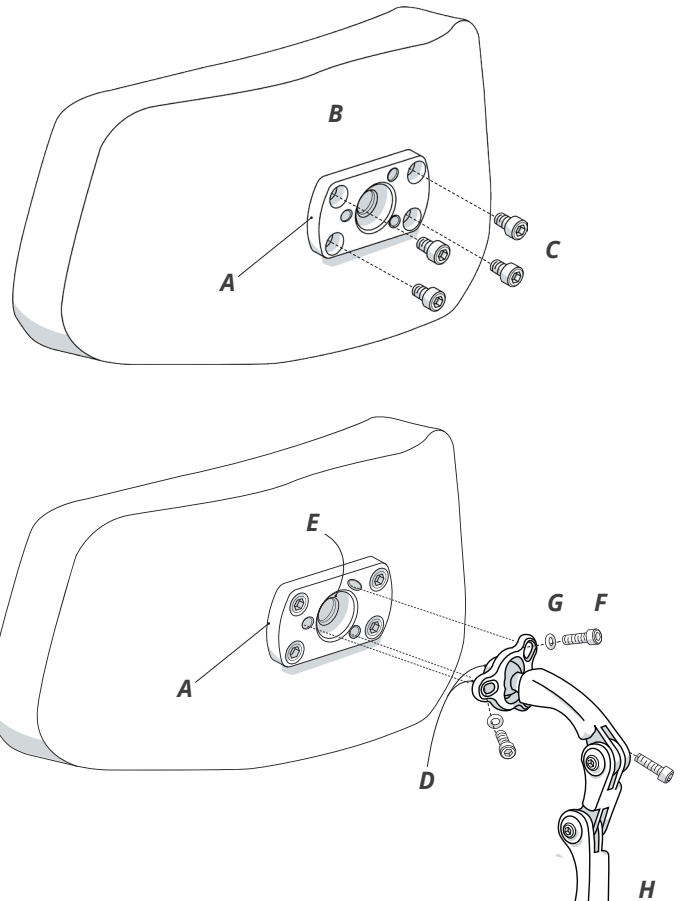
Once all three screws are engaged, tighten them securely.

3 To position the headrest: Loosen all screws at adjustment points, and position as desired.

Note: Loosening the three screws (F) around the ball ring will allow pad orientation adjustment.

(H) = BodiLink Head Support Hardware

WARNING: If the three screws (F) around the ball ring are not evenly tightened, unwanted movement may occur. After installation and adjustment is complete, and before using the headrest, make sure all the screws are evenly tightened.

Hardware Adjustment

Note: Adjustments should be made simultaneously with support pad positioning adjustments.

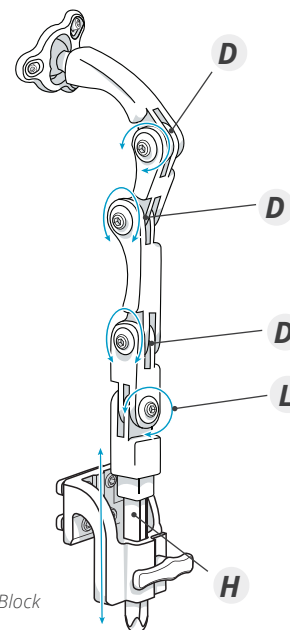
Standard Taper Inserts

- Loosen all of the hardware adjustment point screws, allowing for the full articulation of the head support to a desired position.
- Hold the hardware and the support pad in place while tightening each joint.

Tool-Less Wing Locks (not shown)

- Turn the wing lock(s) counter-clockwise to loosen the hardware adjustment points, allowing for articulation.
- Turn the wing lock(s) clockwise to tighten.

D = Depth
H = Height
L = Lateral adjustment

Shown: Permobil "A" Mounting Block

LTS BodiLink Lateral Trunk Support

SUPPORT PAD ATTACHMENT

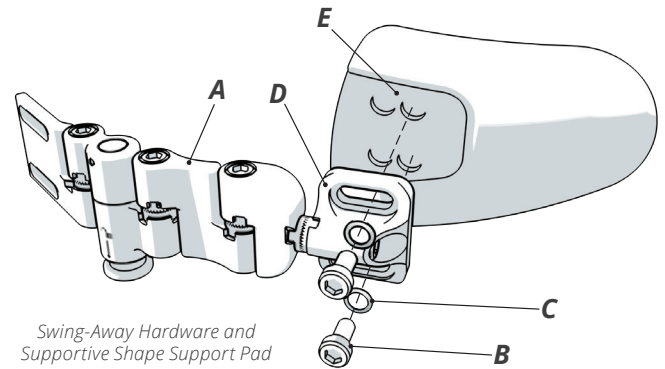
All required tools are provided.

- 1 If the support pad (E) has a cover, start by pulling back the cover flap to expose the holes in the support pad. Line up the bolts (B) and the safety washers (C) with the short slots on the hardware link (D) and the holes inside the support pad (E). The flat face of the link (D) must face the support pad.
- 2 Tighten the bolts (B) using the 4mm hex key.

Note: Follow these instructions if the support pad is already attached to the hardware assembly (A), and depth adjustment of the support pad is needed:

- 1 Pull back the cover flap. Loosen but do not remove the bolts (B) using the 4mm hex key. Slide the support pad (E) into position.
- 2 Tighten the bolts (B) with the 4mm hex key and release the cover flap.

Note: The LTS support pad cover is not shown but may be included with the support pad. See "Cleaning" for cover care.

INSTALLATION - TRACK MOUNT

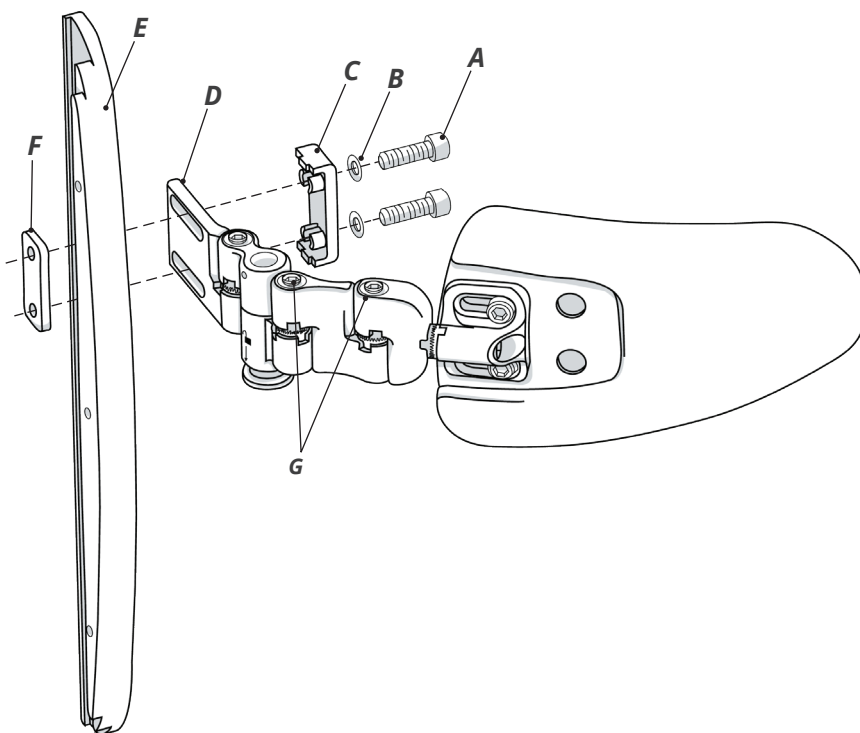
Attaching the hardware to the track.

In some cases, the LTS hardware will already be attached to the track mount bracket (E). If height, width or depth adjustments are needed, loosen but do not remove the bolts (A) using the 5mm hex key, and slide the link (D) into position. Tighten the bolts (A) to secure. This can also be done after the track mount bracket is installed on the Corpus shell.

- 1 Line up the bolts (A) with the safety washers (B), the slot-lock double washer (C), the slots on link (D), the slot on track mount bracket (E) and the slot-lock double nut (F).

Note: The bolts (G) on the LTS hardware should face up, and the curved edge of the track mount bracket (E) should be oriented towards the support pad.

- 2 Tighten the bolts (A) using the 5mm hex key. Follow these steps for each set of hardware needed. Continue to step 3.

Attaching the track to the back shell.

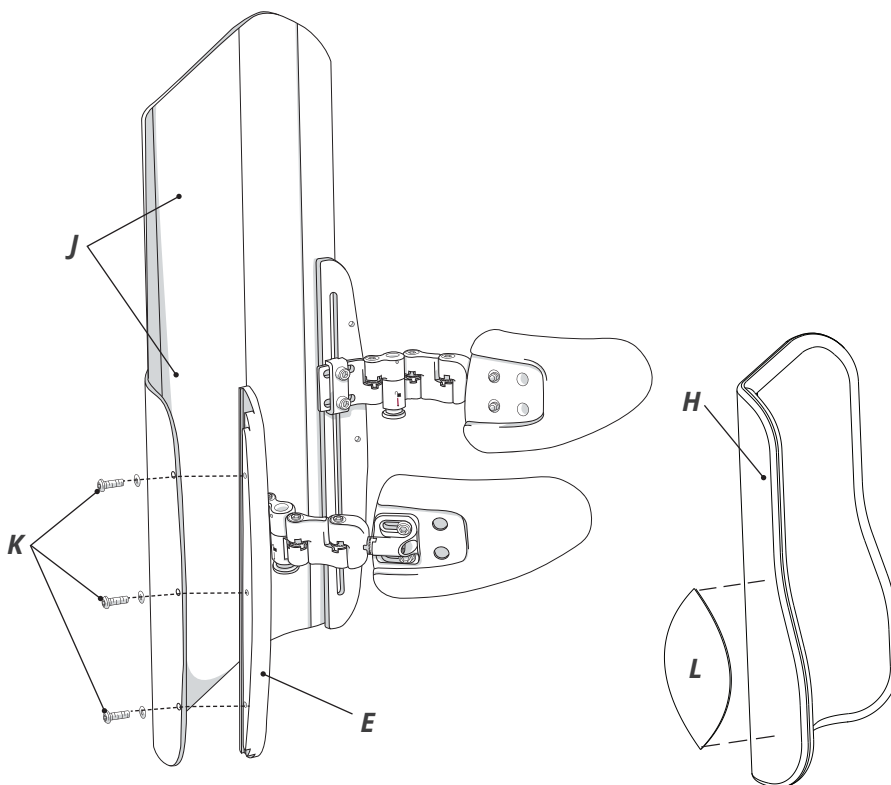
- 3 Remove the back support cushion (H) from the back shell (J), pulling forward to separate the hook and loop attachment.
- 4 Locate the existing 3 bolts along the rear-facing, outside edge of back shell (J), next to the armrest. Remove the 3 bolts to release the existing height adjustment track and the plastic lateral wedges (not shown) attached to the front face of the back shell (J). Discard the removed parts or retain for future use.

- 5 Line up the provided button head bolts (K) with the back shell (J) and the track mount bracket (E). Tighten the bolts (K), using the 3mm hex key.

Recommendation: For better access to the bottom bolt, lift the armrest.

- 6 Adhere the replacement foam support pads (L) to the back support cushion (H). The support pads (L) adhere with a hook and loop attachment and can be placed wherever the lateral support is needed.

- 7 Replace the back support cushion (H) on the back shell (J). Continue to the "Adjustment and Use" section for final hardware adjustments and support pad positioning.

ADJUSTMENT AND USE

Mounting Hardware and Adjustments.

WARNING: Adding additional links to the hardware may reduce the overall strength.

IMPORTANT! After hardware and support pad positioning, make sure that the teeth on all link joints are properly aligned and engaged before final tightening. See the warning in the "Important Safety Information" section in this manual.

H = Height

Loosen the bolt(s) and slide mounting link up and down on the back support shell, mounting track, chair cane or wing slot.

D = Depth

Adjustable in multiple ways:

- Loosen the bolts and slide the mounting link* in and out on the back support shell or wing slot. This may also affect the width.
- Loosen support pad mounting bolts under the support pad cover, and slide the support pad along the slots.
- Remove the bolts (under the cover, if applicable) and reattach the bolts, using alternate threaded holes in the support pad.

W = Width

Loosen the bolts and rotate the hardware links.

* Note: Extra hardware links are available.

A = Angle

Adjustable in multiple ways:

- Loosen the bolts (under the cover, if applicable) and tilt the support pad up and down.
- Loosen the bolt at the end of the support pad mounting link and rotate the support pad.**

** Note: Gear tooth adjustment points are spring loaded for easy adjustment. When the bolts are loosened, the spring will hold the support pad and hardware in place while other areas are tightened.

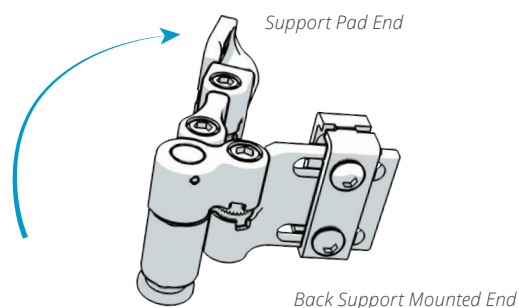
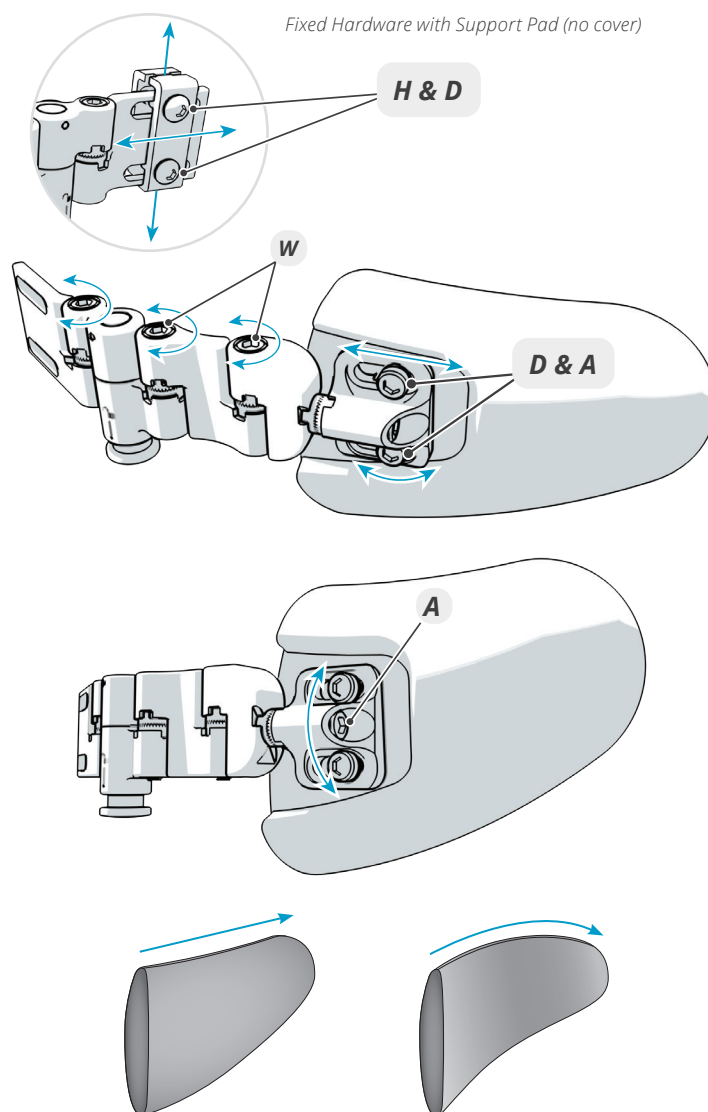
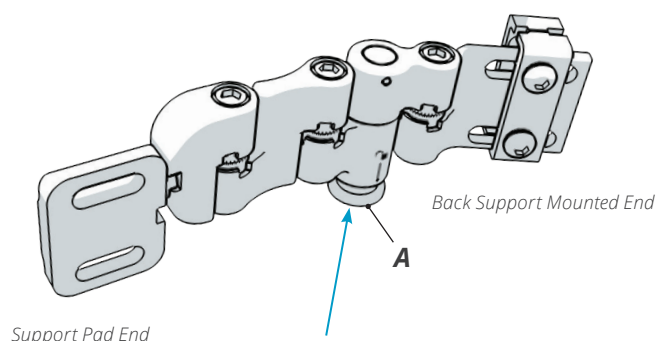
Support Pad Adjustments

The moldable support pad can be shaped by hand or by bending it over a stable surface.

Operating the Swing-Away Feature on the Mounting Hardware

To swing away the lateral support pad, press the silver release button (A), located on the underside of the assembly. The assembly will snap back into place when it is returned to its original position. Before use, make sure that the assembly is engaged.

WARNING: If the assembly is not engaged, it will not be locked into place. If not locked, it will not provide the intended support.

LPTS BodiLink Lateral Pelvic/Thigh Support

INSTALLATION - SUPPORT PAD TO MOUNTING HARDWARE

All required tools are provided.

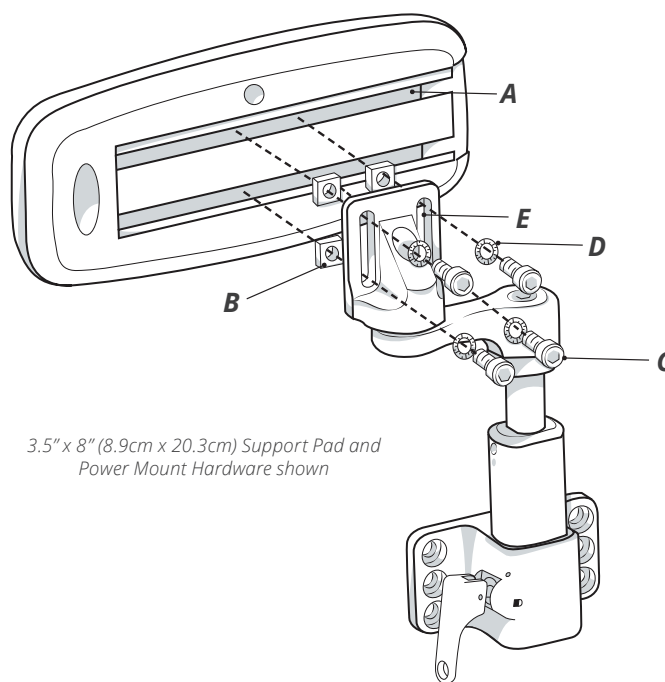
Note: In some cases, the support pad will come installed on the hardware.

- 1 Unzip and remove the cover and foam insert to gain access to the metal tracks (A) for hardware installation. Remove the plastic cap (if installed). Insert the two square nuts (B) into each of the two pad tracks (A).

- 2 Line up the bolts (C), washers (D), and slots on the mounting bracket (E) with square nuts (B) in track (A).

Recommendation: When installing a hardware assembly on a support pad, install the hardware close to the center of the support pad (depth).

- 3 Tighten the bolts (C) using the 5mm hex key. Reinstall the plastic cap (if provided).
- 4 Replace the cover by inserting the foam insert in the cover with the gripping side of the foam facing the support pad, and then add in the support pad and zip the cover closed.

3.5" x 8" (8.9cm x 20.3cm) Support Pad and Power Mount Hardware shown

INSTALLATION - POWER MOUNT HARDWARE

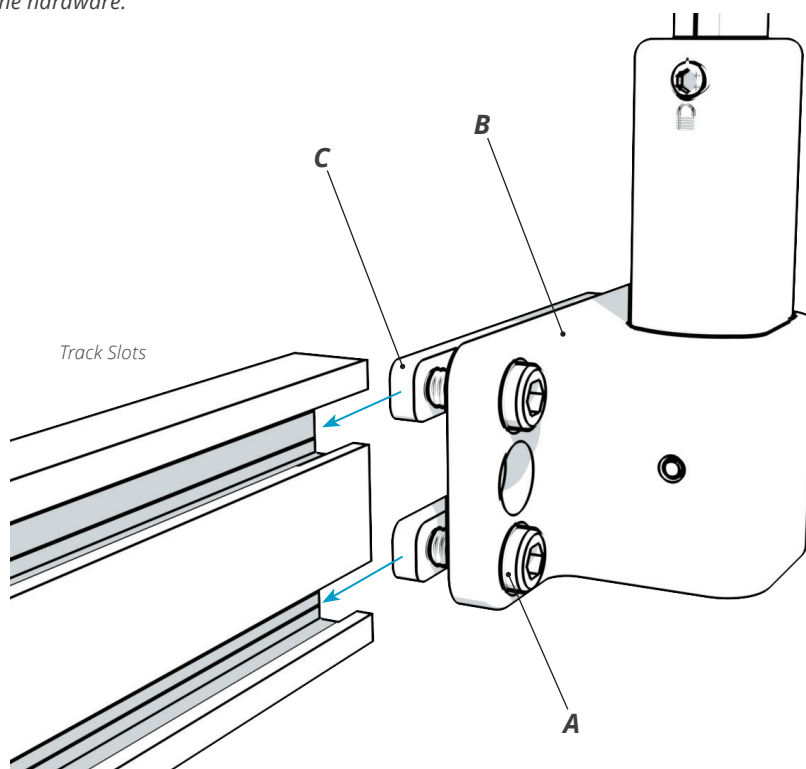
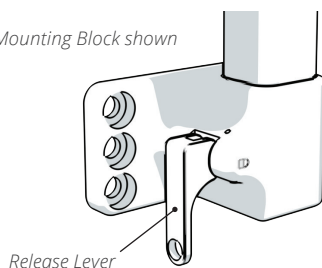
Attaching hardware to the seat frame tracks. All required tools are provided.

Note: In some cases, the support pad will come installed on the hardware.

- 1 Loosen the screws (A), connecting the hardware mounting block (B) to the double nut(s) (C) using the 5mm hex key. A gap between the mounting block and the double nut(s) is needed to slide into the seat frame track.
- 2 Attach the hardware mounting block (B) to the seat frame track by sliding the double nut(s) (C) inside of the seat frame track (requires temporary removal of a cap at the end of the track). Slide the hardware to the general location on the seat frame track for the desired support pad and hardware placement.
- 3 Do not fully tighten the screws (A). Continue with support pad and hardware ADJUSTMENT and USE instructions in this manual.

Note: For removable hardware, the release lever can be installed towards the front or rear of the wheelchair, pointing either up or down, depending on user/caregiver preferences.




Removable Mounting Block shown

ADJUSTMENT - SUPPORT PAD**Adjusting the support pad with hardware.**

- 1 To adjust the position of the support pad depth relative to the mounting hardware, unzip the cover to reveal the support pad mounting bolts. Loosen the bolts (A) using the 5mm hex key.
- 2 Move the support pad to the desired position. Retighten the bolts using the hex key. Zip the cover closed.

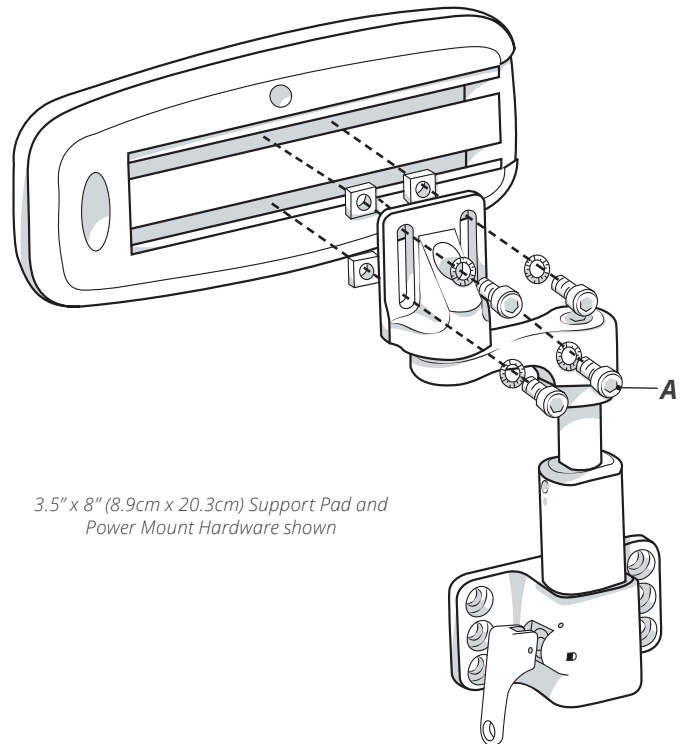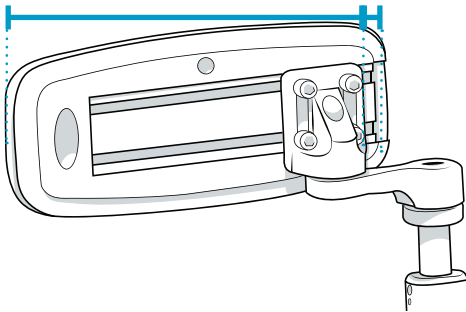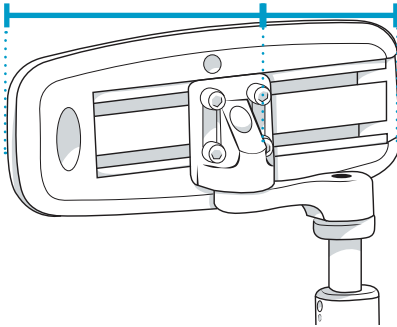Support pad depth:

*3.5" x 8" (8.9cm x 20.3cm) Support Pad and
Power Mount Hardware shown*

ADJUSTMENT - POWER MOUNT HARDWARE

Adjusting Power Mount Hardware.

W = Width/Medial Adjustment

A = Angle

H = Height

D = Depth

1 Loosen the following screws:

H Slider height set screws (each side of the post), using the 3mm hex key.

W Medial arm adjustment screw, using the 5mm hex key.
Note: If no medial positioning is needed, the medial arm can be removed and the support pad mount bracket can be fastened directly to the vertical slider post.

A Pad swivel screw, using the 5mm hex key.
 May need to be adjusted in combination with the above.

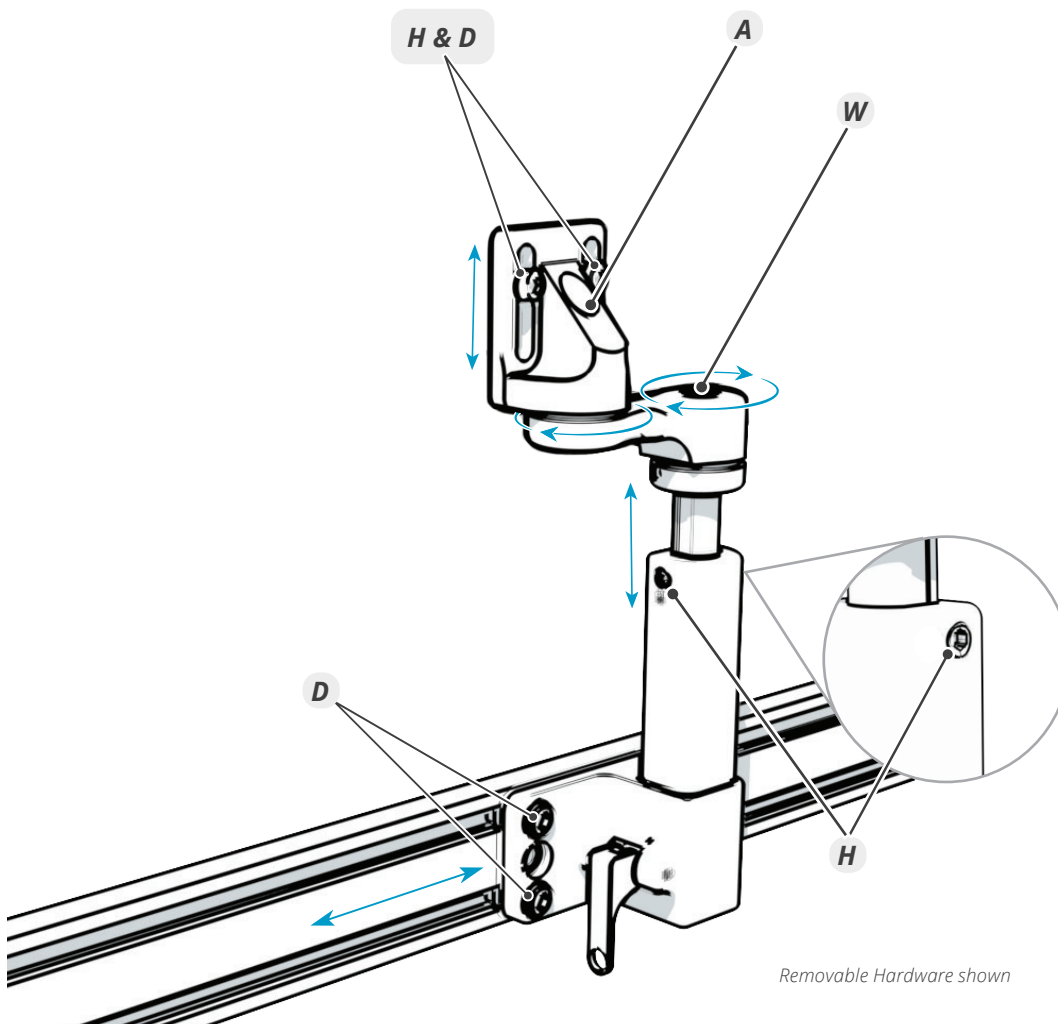
D Mounting block screws (if necessary), using the 5mm hex key.

H & D Pad mount screws (if necessary), using the 5/32" hex key.

2 Position the pad to the desired height, width location, and transverse angle.


Note: The pad and/or medial adjustment arm should not compress the seat surface.

3 Fully tighten all screws that were loosened in step 1, using the appropriate hex key.

Removable Hardware shown

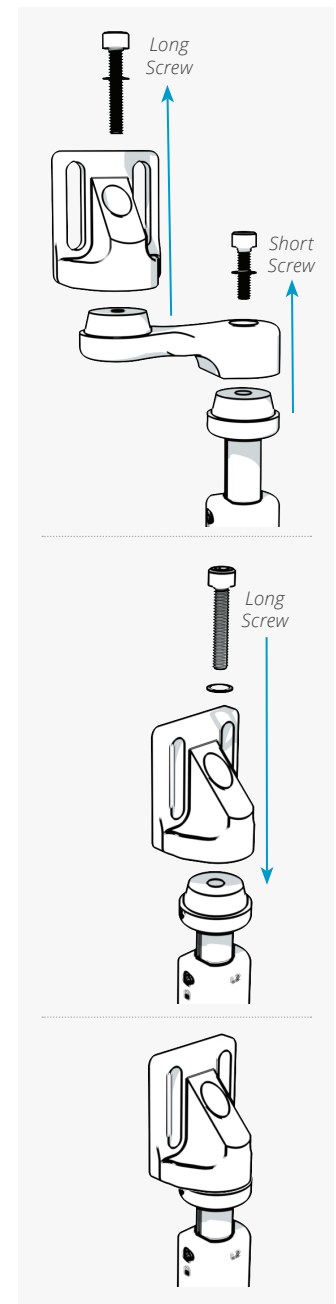Medial Arm Removal




USE - HARDWARE**Removing and reinstalling Power Mount Hardware.**

⚠ WARNING: Hardware should be removed before transferring into the wheelchair.

- 1 To remove the hardware from the mounting block, pull the release lever (A) up (or down) until the lever remains in the open position (~90°). Pull the hardware post (B) straight up and out of mounting block (C).
- 2 To re-install the hardware, insert the post (B) straight down into the mounting block receiver (C) until the post is fully seated.
Note: The release lever can be open or closed when inserting the post.
- 3 Once the hardware is fully inserted into the mounting block, ensure that the release lever is fully closed and the hardware is engaged before use.

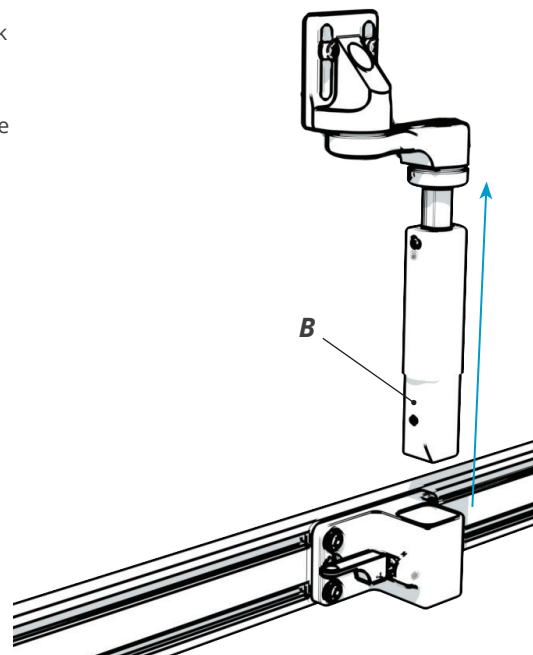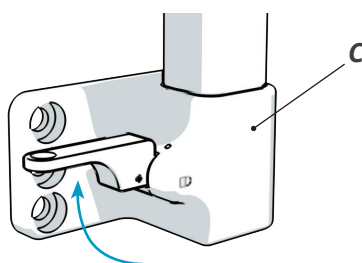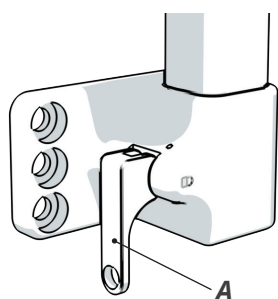Note: Release levers on the hardware have a hole near the end for an optional trigger pull lanyard. (Lanyard is not provided and should only be used at a therapist's discretion.)




**USE - HARDWARE RELEASE LEVER**

Switching the removable hardware release lever orientation. The lever location and orientation are subject to user/caregiver preference and ability to operate.

- 1 Firmly grasp the sides of the release lever base and pull straight out.
- 2 While the lever is pulled out from the mounting block/base, rotate the lever 180° and release in the new position.

Note: To switch the location of the lever from the front to the rear of the mounting base/block, switch the left and right-hand side hardware (or rotate 180° for power mount) during installation. Refer to the Hardware installation instructions in this manual.

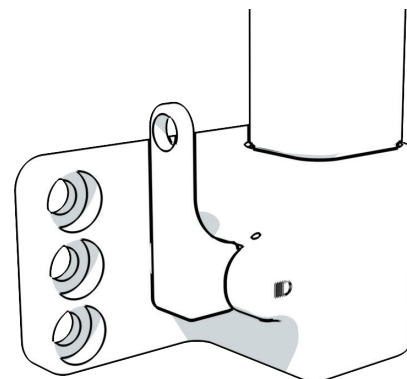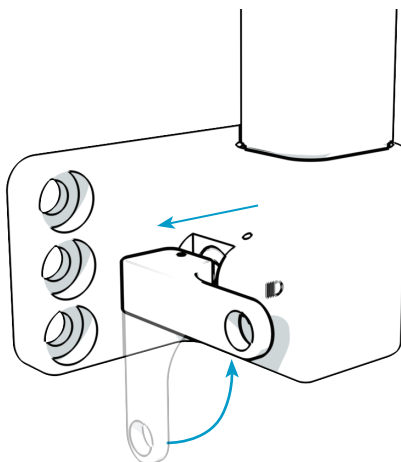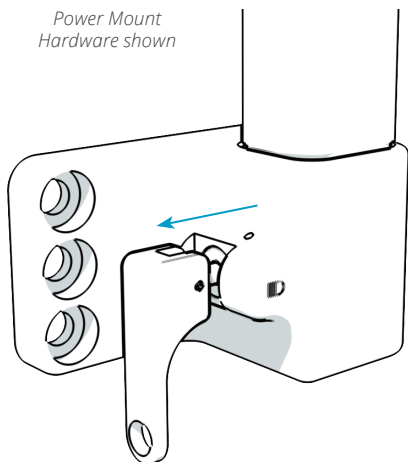Power Mount
Hardware shown

CLEANING**WARNINGS:**

- Follow all manufacturer bleach container and germicidal disinfectant safety guidelines.
- Before transferring to a new user, the product must be completely cleaned and disinfected.
- Periodically, check the cover and components for damage and wear, and replace if necessary. Some wear with regular use is expected.

CAUTION: Frequent heat drying and bleaching will cause fabric to break down. Air drying and cleaning without bleach are recommended whenever possible to help extend the life of the cover.

CLEANING AND DISINFECTING ALL FABRICS

CAUTION: Remove the cover from the foam and any inserts before cleaning or disinfecting.

Cleaning: When possible, hand wash in warm water with mild soap, rinse thoroughly with clean water, and air dry (let dry completely before use). Or machine wash using common laundry detergent and tumble dry low.

Disinfecting: Machine wash hot (without the foam insert) using common laundry detergent. Tumble dry low. Let dry completely before use.

**FOAM**

WARNINGS: DO NOT use contaminated foam with multiple users. Disinfectants are NOT effective on porous surfaces such as foam. Replacement foam components may be purchased separately.

CAUTIONS: DO NOT submerge the foam in any liquid, including water. The foam should not be exposed to light for any extended period of time. Light will cause discoloration of the foam.

Note: Discoloration will not affect the functional features. There is no need to attempt to "clean" discoloration caused by light exposure.

To clean the foam: Soiled foams should be replaced. Remove the cover. Wipe lightly with a damp cloth. DO NOT use soap. Let dry completely before replacing the cover.

**CLEANING AND DISINFECTING HARDWARE**

CAUTION: Before cleaning or disinfecting, remove the cover support pad from the hardware.

Cleaning: Hand wash in warm (40°C/104°F) water with soap and a sponge. Rinse with clean water and wipe dry with a clean cloth. Dry completely before use.

Disinfecting: Wipe gently with a cloth dampened with a household disinfectant. Rinse with clean, warm water. Wipe dry with a clean cloth. Dry completely before use.

**CARE SYMBOLS**

Wipe with a damp cloth.

Hand wash.

Dry flat.

Do not tumble dry.

Tumble dry low.

Dry flat in the shade.

Machine wash hot, normal*.

Machine wash warm, normal*.

Hand wash warm*,

* maximum temperature shown

MAINTENANCE

INSPECTION

Inspect all components of the equipment weekly.

Check and tighten all fasteners and adjustment points.

Inspect all materials for wear such as stretching, fraying, fractures, and bends. Some wear with regular use is expected.

Components that show excessive wear should be replaced. For information on replacements, refer to the "Warranty" section in this manual.

WARNING: DO NOT continue to use the product if any part is broken or cannot be tightened. Periodically, the product should be thoroughly inspected by an equipment provider.

REUSE

You must re-fit the support to the new user or after extended periods of non-use by the same user. Refer to the adjustment and use instructions in this manual.

Before transferring to a new user, the product must be completely cleaned and disinfected. Refer to the cleaning and disinfecting instructions in this manual.

STORAGE

Assembly or individual parts should not be stored in extreme hot or cold temperatures. Only store in a clean and dry place at room temperature.

Do not allow exposure to direct sunlight for prolonged periods of time. This could age the product more rapidly and decrease its effectiveness.

Do not store where there is risk of open flame or spark.

DISPOSAL

Dispose of the product and/or components in accordance with the applicable regulations in your jurisdiction.

CAUTION! DO NOT INCINERATE as a means of disposal.

WARRANTY

This product is covered by a lifetime warranty against manufacturer defects. If a manufacturing defect should occur, discontinue use immediately; note that the product will be replaced at no cost to the buyer. Wear from regular use is not considered a manufacturing defect. Replacement parts are available by contacting your supplier, distributor, or retailer.

The following will void the warranty: modifications made in attempt to interface with other products (where compatibility is limited, as described in product information); alterations made to the product; or defects caused by irregular use.

WARRANTY REPLACEMENT PROCEDURE

Warranty claims should be initiated by the original purchaser, who should contact the DME supplier, distributor or retailer from whom the product was purchased, if they are available, or any authorized Permobil DME supplier, distributor or retailer.

Whenever possible, the DME supplier, distributor or retailer with a warranty issue should contact Customer Support and provide the original purchase order number, sales order number, or invoice number. Before a warranty replacement order or quote may be processed, Customer Support must issue a Return Merchandise Authorization (RMA).

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INFORMACIÓN DE CONTACTO

The Comfort Companies, LLC
1960 South Calhoun Road
New Berlin, Wisconsin 53151 EE. UU.

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The Comfort Companies, LLC pertenece a Permobil Seating and Positioning.

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