
ROHO®

AGILITY® Max Contour Back Support

Operation Manual

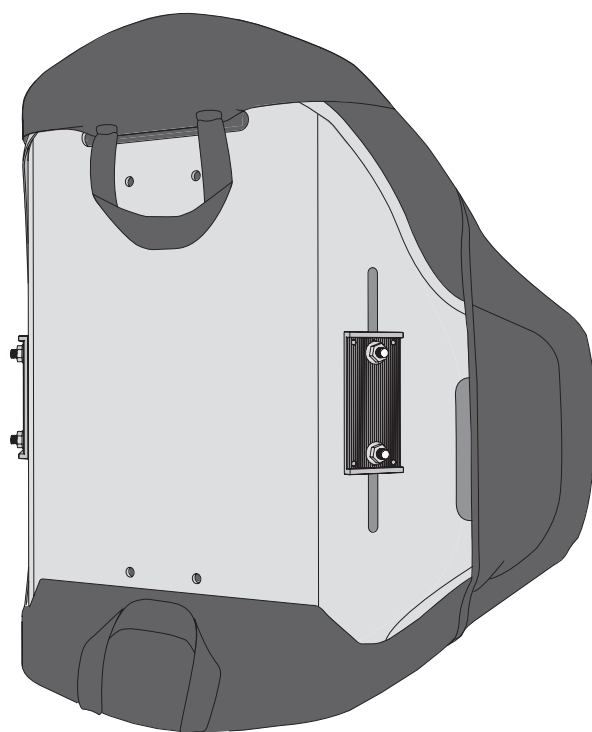permobil

s h a p e f i t t i n g t e c h n o l o g y®

Contact Information

 ROHO, Inc.
1501 South 74th St.
Belleville, IL 62223 U.S.A.

U.S.A.:
800-736-0925
Fax 888-551-3449

Customer Support:
orders.roho@permobil.com

Outside the U.S.A.:
618-277-9150
Fax 618-277-6518
intlorders.roho@permobil.com
permobil.com

 Warning

Indicates that not following the specified procedure could lead to potentially hazardous conditions resulting in serious injury.

 Caution

Indicates that not following the specified procedure could lead to potentially hazardous conditions resulting in minor to moderate injury or damage to the equipment or other property.

Incident Reporting

In the event of a serious incident or injury, contact Customer Support.

Contents

Intended Use	2
Important Safety Information	2
Parts Detail	3
Product Specifications	4
Cover Removal and Installation	5
Cover Components Installation	6
Cover Components - Air Insert Adjustment	7
AGILITY Lumbar Support Pad	
Troubleshooting	8
Storage, Transport, Disposal and Recycling	
Limited Warranty	
Cleaning and Disinfecting	9
Español	10
Français	19
Deutsch	28
Italiano	37
Nederlands	46
Dansk	55
Svenska	64
Norsk	73
Polski	82
Português	91

Symbols Glossary

Batch code

Medical Device

Manufacturer

Meets flammability

Model number

Serial number

Date of manufacture

Size (width x depth x length)

Cleaning instructions

Fits

Consult instructions for use

Weight limit

Disinfecting instructions

Quantity

Authorized representative in the European Community

Authorized representative in Switzerland

Please Note

Supplier: This manual must be given to the user of this product.
Operator (Individual or Caregiver): Before using this product, read the instructions and save for future reference.

Products may be covered by one or more U.S. and foreign patents and trademarks, including ROHO®, AGILITY®, and shape fitting technology®.

Product information is changed as needed; current product information is available at permobil.com. Contact Customer Support for previous versions of product information.

ROHO, Inc. is part of Permobil.

© 2017, 2022 Permobil

Rev. 2022-08-02

Includes instructions for ROHO AGILITY Back Shell, Cover, and Accessories. Use this operation manual with the instructions provided with the ROHO AGILITY Back Support Hardware.

Intended Use

The ROHO AGILITY Max Contour Back Support (AGILITY Max) is a replacement backrest for wheelchairs that is intended to provide trunk/pelvis support and positioning. AGILITY has a weight limit of 300 lbs. (136 kg) and must be properly sized to the individual.

 **Compatibility:** AGILITY Max is NOT intended to be used with any laterals. See "Product Specifications" for compatible attachment hardware and accessories.

A clinician who is experienced in seating and positioning should determine whether the product is appropriate for the individual's particular seating needs. The product's intended use statement alone should not be relied upon to make that determination.

The medical products featured in this document are intended to be only one part of an overall care regimen that includes all seating and mobility equipment and therapeutic interventions. A clinician who is experienced in seating and positioning should identify this care regimen after evaluating the individual's physical needs and overall medical condition.

A clinician should also evaluate for visual, reading, and cognitive impairment to determine the need for caregiver assistance or the need for other assistive technology, such as large-print instructions, to ensure proper product use.

Important Safety Information

Warnings:

- DO NOT use the product without first reading and understanding all of the provided instructions. If you are unable to perform any tasks in this manual, seek assistance by contacting a clinician, equipment provider, distributor, or Customer Support.
- Follow all instructions and safety information provided with the product and provided by the manufacturer of the wheelchair or other accessories.
- Skin/soft tissue breakdown can occur due to a number of factors, which vary by individual. Check skin frequently, at least once a day. Redness, bruising, or darker areas (when compared to normal skin) may indicate superficial or deep tissue injury and should be addressed. If there is any discoloration to skin/soft tissue, STOP USE immediately. If the discoloration does not disappear within 30 minutes after disuse, immediately consult a healthcare professional.
- The product should be installed by an equipment provider or clinician.
- The product must be sized correctly to the individual and the wheelchair.
- The AGILITY is intended to be used with an AGILITY cover, foam pad, lumbar support pad (if applicable), and either the air inserts or the foam inserts. All components MUST be compatible sizes and MUST be used as directed in this manual. Use of any hardware or components other than those approved for use with the product may result in product damage or failure and may void the warranty. Compatible components and accessories are listed in "Product Specifications" in this manual.
- As with any addition to a wheelchair, the product may affect the wheelchair's stability. Before using the product, consult with a clinician or equipment provider to determine whether the wheelchair should be modified.
- DO NOT lean over the top of an AGILITY Back Support. Doing so may cause a change in the center of gravity, and the wheelchair could tip over.
- DO NOT lift or push a wheelchair using any part of an AGILITY because it may become detached from the wheelchair.
- DO NOT expose the product to high heat, open flames or hot ashes. Testing or certification claims, including for flammability, may no longer apply to this device when it is combined with other products or materials. Check testing and certification claims for all products used in combination with this device.
- Metal components may become hot if exposed to high temperatures or cold if exposed to low temperatures. The cover assembly surface will acclimate to the surrounding temperature. Take appropriate precautions, especially when the product will contact unprotected skin.
- Altering the product may affect its performance, result in product failure, and void the warranty.
- If a headrest is necessary, follow the headrest manufacturer's instructions and safety information. Note: Before installing another manufacturer's headrest, consult with a clinician or service provider to ensure the headrest is compatible with the AGILITY.

Cautions:

- Protect the AGILITY from blunt force shocks that may cause damage and/or breakage.
- The components of an AGILITY are not intended to be used separately from an AGILITY or in conjunction with other wheelchair back supports.
- Check AGILITY components for wear or damage. DO NOT use an AGILITY if any component is not functioning properly or is broken. Periodically, the AGILITY should be thoroughly inspected by a qualified technician.
- Carry the AGILITY by the cover handle or by firmly grasping the edges of the back shell.

Warnings - Motor Vehicle Transportation Important Safety Information:


- Whenever possible, transfer out of your wheelchair installed with AGILITY and into a manufacturer-installed vehicle seat, and use the vehicle's crash-tested occupant restraint system.
- To determine whether motor vehicle third-party testing has been performed with a specific AGILITY-hardware or other accessory combination, refer to the operation manual provided with AGILITY, attachment hardware, and other accessory, or contact the equipment manufacturer.

Parts Detail

Package contents: back shell, back shell attaching hardware*, cover, foam pad and inserts, hand pump, repair kit, AGILITY Lumbar Support Pad (select models only), operation manual, limited warranty supplement, product registration card

To order replacement parts, contact the product provider or Customer Support.

*For additional product information, refer to the instructions provided with the AGILITY back support hardware. 

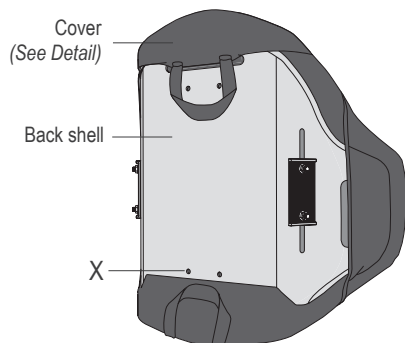
Pre-drilled holes (X) for attaching a universal headrest mount (x4).

Note: Pre-drilled holes are not available on AGILITY models that are 10 in. (25.5 cm) and 13 in. (33 cm) in length.

Back Shell shown with attaching hardware.*

Cover

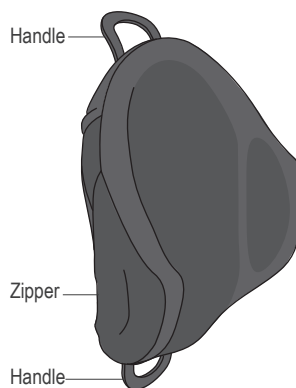Hand pump


AGILITY Lumbar Support Pad
(select models only)


Repair kit

AGILITY Max Cover Components

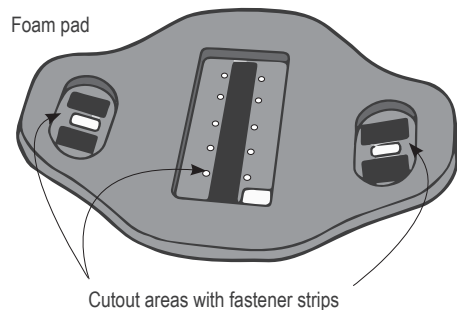
Foam insert

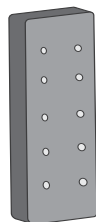Wing air insert


Back air insert

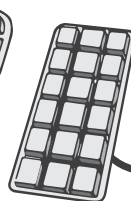Wing air insert

Valve hose and inflation valve (x3)

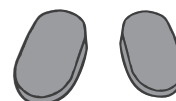Wing foam inserts

Product Specifications

All measurements are approximate.

1 in. = 2.54 cm = 25.4 mm

Adjustable Width (W) Fits Wheelchair Width Range - in. (cm):

13 to 14 (33 to 35.5)	Width (W) measurements are taken from the outside of the back cane to the outside of the back cane.
14 to 15 (35.5 to 38)	
15 to 16 (38 to 40.5)	Adjusting AGILITY to the smallest width within a wheelchair width range will result in ½ in. (1.5 cm) loss of seat depth.
16 to 17 (40.5 to 43)	
17 to 18 (43 to 45.5)	
18 to 19 (45.5 to 48.5)	
19 to 20 (48.5 to 51)	


Note: The width measurement information above applies only when using the product with AGILITY Fixed Hardware or AGILITY Quick Release Hardware. It does not apply when using the product with AGILITY Direct Mount Hardware. For further information, refer to the instructions provided with AGILITY Direct Mount Hardware.

AGILITY Length (L) - in. (cm):

16 (40.5)	Length (L) measurements are taken from the center edge to the center edge between the narrow and wide ends of the back shell.
18 (45.5)	
20 (51)	

Note: The cover assembly is not included in the length measurement.

Shown with ROHO AGILITY Back Support Quick Release Hardware

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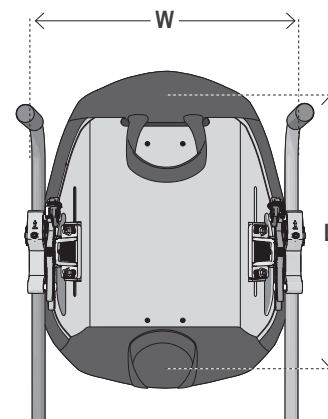
Depth Range: 6½ to 7¼ in.; 17-18 cm

AGILITY Back Supports are measured for depth (D) with covers and with air inserts that are inflated to allow for ½ in. (1.5 cm) of immersion. Measurements are taken from the center, front surface to the highest point of the the contoured edge.

Weight Range (lb.): 4.25 to 6½ lb. ; 2-3 kg

Depths and weights are based on W x L in. (cm): 14 x 16 (35.5 x 40.5); and 20 x 20 (51 x 51)

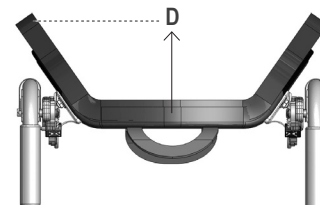
Adjustable Back Air Insert

Width (BW): 3 cells, 6.5 in. (16.5 cm)

Cell Depth (CD) when fully inflated: 1¼ in. (3 cm)

	Length (BL)		
# Cells	6	7	8
in.	12½	14½	16½
(cm)	(32)	(37)	(42)

Perimeter included in width and length measurements: ¼ in. (1 cm)

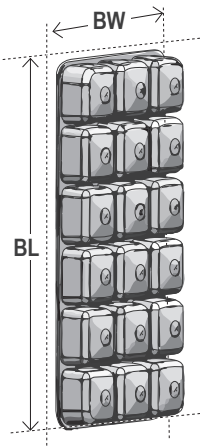
Adjustable Wing Air Inserts

Width (WW): 4¼ in. (10.5 cm)

Length (WL): 7 in. (17.5 cm)

Cell Depth (CD) when fully inflated: 1¼ in. (4.5 cm)

Perimeter included in width and length measurements: ¼ in. (1 cm)

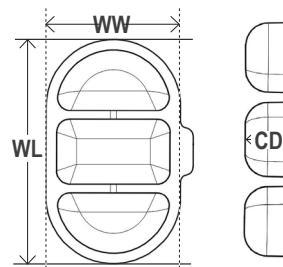Materials:

Back Shell: aluminum

Attaching Hardware: steel, aluminum, Zytel

Cover: spacer fabric front and back center panel, fluid-resistant back panels and flaps, with hook and loop fasteners to secure the cover to the back shell. Refer to cover care label for material content and for compliance to flammability requirements.

Foam Pad: polyurethane front; polyethylene back

Foam Insert: polyurethane

Air Inserts: polyether polyurethane front; urethane-coated nylon back

AGILITY Lumbar Support Pad: polyurethane

Expected product lifetime: 5 years

Attachment Hardware and Accessories Compatible with AGILITY Max:

Refer to the instructions provided with the hardware or other accessory to determine limitations or restrictions.

- ROHO AGILITY Back Support Quick Release Hardware
- ROHO AGILITY Back Support Direct Mount Hardware
- ROHO AGILITY UniTrack Rail

Not available in all markets; contact Customer Support:

- ROHO AGILITY Headrest Adapter Plate
- ROHO AGILITY Back Support Fixed Hardware
- ROHO AGILITY Lumbar Support Pad
- Therafin Headrest for ROHO AGILITY Back Support


Cover Removal and Installation

⚠ Warnings:


- DO NOT remove or install the AGILITY cover while the individual is in the wheelchair.
- DO NOT use an AGILITY without the cover, foam pad, and, if applicable, the air insert or foam insert.
- Inspect the cover for damage and replace if necessary.

To Remove the Cover



1. On the wide end of the shell: Pull down to release the cover flap.

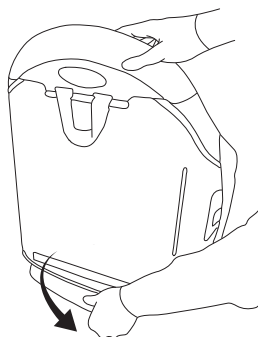2. Pull the cover away from the shell.


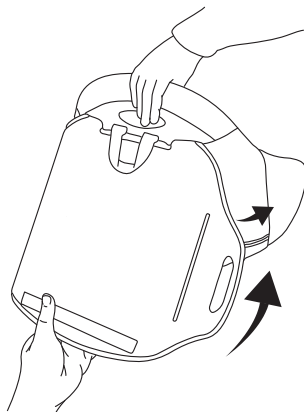
3. Lift the cover over the narrow end of the shell. Pull to release the hook fastener strip on the tab. Pull the tab through the slot.

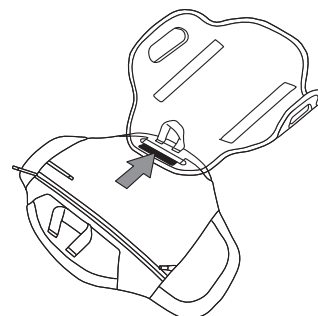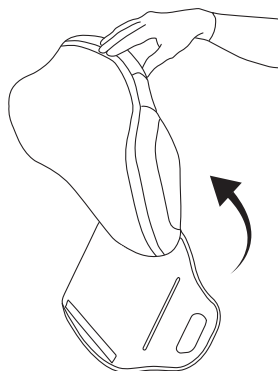
To Install the Cover


1. On the back side of the shell, push the tab, but not the handle, through the slot at the narrow end. Press to attach the tab to the shell. Lift the cover over the narrow end of the shell.

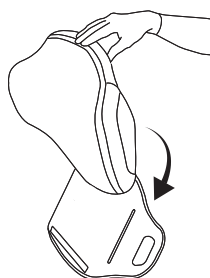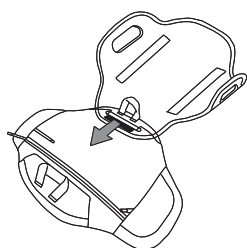
2. Make sure the zipper faces the shell. Slide the sides of the back shell into the side pockets of the cover. Press to attach the cover to the front side of the shell.

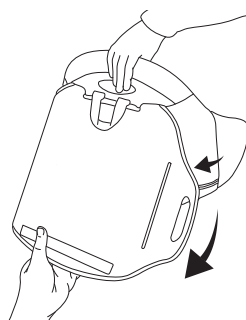Note: If the air insert will need to be adjusted after installing the AGILITY on the wheelchair, see "Cover Components Installation": "Installing the Cover on the Back Shell - Valve Hose Placement".

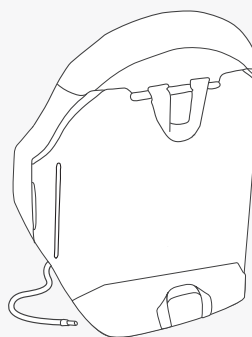3. Pull the wide end of the cover over the wide edge of the shell. Press to attach the cover flap to the shell.

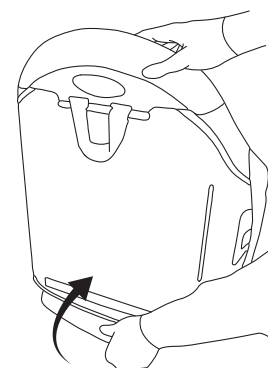To Remove the Foam Pad, Air Inserts or Foam Inserts From the Cover

1. Unzip the cover.
2. Remove the foam pad from inside the cover with the inserts still attached.
3. Detach the inserts from the fastener strip in the cutout areas of the foam pad.
4. For air inserts: Pull the valve hose back through the holes in the foam pad.


Cover Components Installation

Warning: The air insert and the foam insert are not intended to be used together.


Warning: DO NOT install or remove AGILITY cover or components while the individual is in the wheelchair.

To Install the Foam Pad with Air Inserts

- Slide the valve hose through the hole in the foam pad. The air cells should face out.

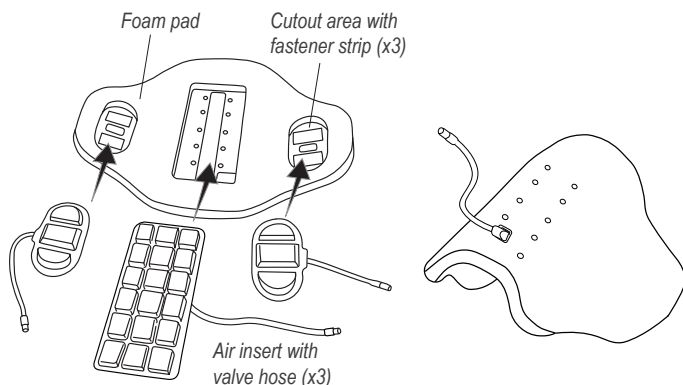- For each air insert, tuck the edges of the air insert under all of the edges of the cutout area. Press to secure the air insert to the fastener strip in the cutout area.
- Unzip the cover. Place the foam pad and the air inserts into the cover.
 - The air inserts should face the front of the cover.
 - The wide and narrow ends of the foam pad should be aligned with the wide and narrow ends of the cover.
- To complete installation of the cover on the back shell, follow the instructions below: "Option A" or "Option B" and "Wing Air Inserts".

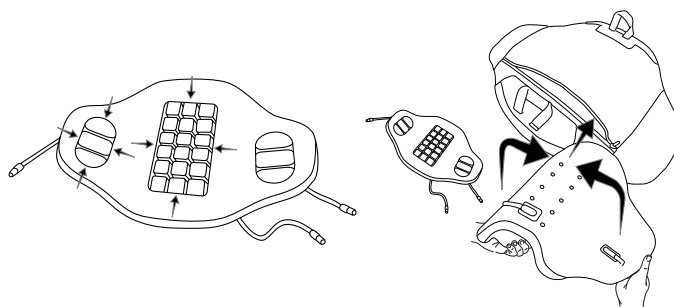
Installing the Cover on the Back Shell - Valve Hose Placement

Air Insert - Option A:

If no air insert adjustment is needed after the cover is installed on the back shell:

- Place the valve hose inside the zipper.
- Zip the cover closed.
- Install the cover on the back shell.
- See "Cover Removal and Installation".

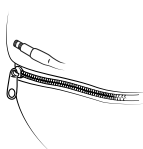Air Insert - Option B:


If the air insert will need to be adjusted after installing the AGILITY on the wheelchair:

- Leave the zipper partially open.
- The valve hose for the air insert should be routed down the INSIDE of the back shell, and then around the bottom of the back shell, so that it is accessible on the OUTSIDE of the back shell.
- Tuck the valve hose into the side pocket.
- Install the cover on the back shell. See "Cover Removal and Installation".

Wing air inserts:


Tuck the valve hoses into the side pockets, as shown.

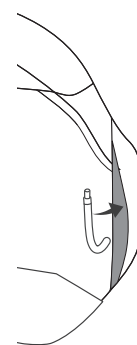To Install the Foam Pad with Foam Inserts

- Align the foam inserts with the cutout areas. Press to secure the fastener strips.
- Unzip the cover and place the foam pad, with the foam inserts attached, into the cover.

Note: The side of the foam pad with the foam inserts should face the front of the cover. The wide and narrow ends of the foam pad should be aligned with the wide and narrow ends of the cover.
- Zip the cover closed.
- Install the cover on the back shell. See "Cover Removal and Installation".

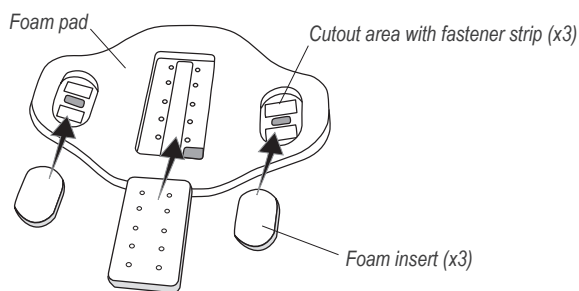Cover Components - Air Insert Adjustment

⚠ Warnings:

- Check inflation frequently, at least once a day!
- The air insert must be used with the air cells facing the individual.
- The individual's comfort level may be reduced if an air insert is not positioned correctly, is not properly sized to the foam pad, or is under-inflated or over-inflated.
- If an air insert does not appear to be holding air, see "Troubleshooting". If the air insert is still not holding air, replace the air insert with the foam insert, and contact an equipment provider or Customer Support.
- Changes in altitude may require adjustment to an air insert. Check the inflation level of any air inserts when changing altitude.

⚠ Cautions:

- DO NOT use the valve hose as a handle to carry or pull an air insert.
- DO NOT use an air insert as a water flotation device (e.g., a Life Preserver).
- If an air insert has been in temperatures lower than 32° F (0° C) and exhibits unusual stiffness, allow the cushion to warm to room temperature. Repeat adjustment instructions before using.
- Keep air cells away from sharp objects.
- DO NOT use an air insert as a seat cushion. DO NOT use an air insert and a foam insert together.
- DO NOT use a pump or repair kit other than those approved for use with an air insert.

To Adjust the Air Inserts



Notes: Before you adjust the air inserts, make sure the cover with the foam pad and air inserts is correctly installed on the back shell. Make sure that the AGILITY is correctly installed and adjusted on the wheelchair.

Follow these instructions to adjust the air insert and the two wing air inserts.

1. Slide the hand pump nozzle over the inflation valve. Open the inflation valve, at least one full rotation. Inflate the air insert until it feels firm.
2. Close the inflation valve until it comes to a stop. Remove the hand pump. Make sure the inflation valve is closed.

⚠ Warning: Safely transfer the individual into the wheelchair.

3. Adjust the air insert, if desired:
 - a. The individual may be seated in the wheelchair in their normal position.
 - b. Open the inflation valve and remove air until the individual is comfortable. Close the inflation valve. Make sure the inflation valve is closed. Follow the valve hose placement instructions.

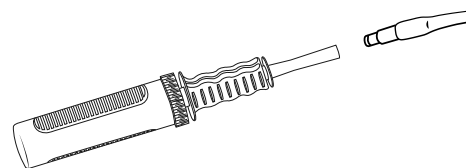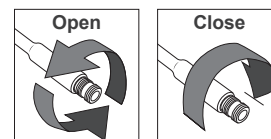
Valve Hose Placement

With the cover installed on the back shell, place the valve hoses in the cover:

Air insert:


Note: The valve hose should be accessible, as shown. If the valve hose is not accessible, see "Cover Components Installation": "Installing the Cover on the Back Shell - Valve Hose Placement".

Tuck the valve hose into the side pocket.

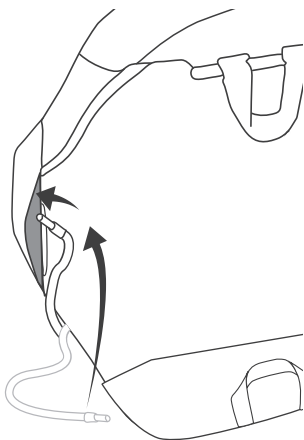Wing air inserts:

Tuck the valve hoses into the right and left side pockets.

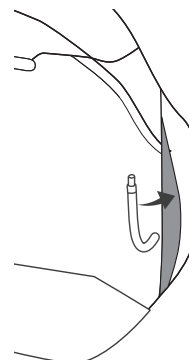AGILITY Lumbar Support Pad*

*Not included with all models.

Notes: The AGILITY Lumbar Support Pad provides added support to the lumbar region.

Package contents include hook and loop fastener strips. The hook side of the fastener strip is pre-attached to the back of the Lumbar Support Pad.

Warning: DO NOT install while the individual is in the wheelchair.

To Install

1. Detach the cover flap from the back shell and unzip the cover. Place the lumbar support pad inside the cover, at the curve of the spine, or where it is most comfortable.

Note: The Lumbar Support may be placed either in front of or in back of the foam pad. Placement behind the foam pad is recommended for models with the air insert. Placement in front of the air insert may reduce the performance of the air insert.

2. Remove the backing from the loop side of the fastener strip, and secure the lumbar support pad in the desired position.

Note: If the loop fastener strip is attached to the inside of the cover, the fastener strip may need to be replaced after multiple cover washings.

3. Zip the cover closed and re-attach the cover flap on the back shell.

To Remove

1. Detach the cover flap from the back shell and unzip the cover.
2. Carefully detach the hook and loop fastener strips from each other.

Troubleshooting

For additional assistance, contact your equipment provider, distributor, or Customer Support.

The air insert is not holding air.	Confirm that the inflation valve is completely closed. Inspect for inflation valve and hose damage. Look for holes in the air insert. If very small holes or no holes are visible, follow the instructions in the Repair Kit provided with the product. For damage to the inflation valve or for large holes or leaks in the air insert, see the Limited Warranty supplement, or contact Customer Support.
Not able to inflate or deflate an air insert.	Confirm that the inflation valve is open. Inspect the inflation valve and hose for damage.
The air insert does not feel comfortable.	Add or remove air to adjust to your comfort level. (See "Cover Components - Air Insert Adjustment"). Make sure air cells face the front of the wheelchair. Make sure the cover is properly positioned (see "Cover Removal and Installation" and "Cover Components Installation").
The air insert slides inside the cover.	Make sure the air insert is properly positioned and the fabric back is secured to the fastener strip in the cutout area of the foam pad (see "Cover Components Installation"). The perimeter of the air insert should be tucked in between the foam layers around the edge of the cutout area.
A component is damaged.	Contact Customer Support for a replacement.

Storage, Transport, Disposal and Recycling

Storage or Transport: Clean and disinfect the product. Open the inflation valves. Deflate the product. Place in a container that protects the product from moisture, contaminants and damage. After storing or transporting, inspect the product to ensure it was not damaged (see "Troubleshooting"), and follow instructions to prepare the product for use.

Disposal: The components of the products in this manual are associated with no known environmental hazards when properly used and when disposed of in accordance with all local/regional regulations. At the end of its useful life, treat as healthcare waste and dispose of according to local healthcare-waste guidelines. Incineration must be performed by a qualified, licensed waste management facility.

Recycling: Contact a local recycling agency to determine recycling options for the product.

Limited Warranty

Limited Warranty Term, from the date the product was originally purchased: back shell: 60 months; Cover, Foam Pad, Foam Insert, Lumbar Support Pad and Air Insert: 18 months. See also the Limited Warranty supplement provided with the product, or contact Customer Support.

Cleaning and Disinfecting

Warnings:

- Cleaning and disinfecting are separate processes. Cleaning must precede disinfection. Before use by a different individual: clean, disinfect, and check product for proper functioning.
- Clean the product regularly and when it becomes soiled.
- If instructed to use bleach, follow the safety guidelines on the bleach container.

Cautions:

- DO NOT use abrasives such as steel wool or scouring pads on components of an AGILITY.
- DO NOT use cleaning products containing petroleum or organic solvents [i.e. acetone, toluene, Methyl Ethyl Ketone (MEK), naphtha, dry cleaning fluids, adhesive removers] on components of an AGILITY.

Recommendation: When washing by hand, use room-temperature water.

Back Shell

To Clean: Hand wash with a sponge, using soap and water. Rinse with clean water. Wipe dry with a clean cloth and then allow to air dry completely.

To Disinfect: Wipe gently with a cloth dampened with a household disinfectant. Rinse with clean water. Wipe dry with a clean cloth and then allow to air dry completely.

Cover

Notes: Before cleaning and disinfecting, remove the foam and air components from the cover. Secure the hook and loop fastener strip on the handle at the wide end of the cover. Turn the cover inside out. Zip the cover closed to avoid snags.

To Clean: Machine wash cold, with mild detergent, gentle cycle, or wipe clean with neutral detergent and water. Rinse thoroughly. Tumble dry low or air dry.

To Disinfect: Hand wash, using quaternary ammonium disinfectant (follow product label instructions). Rinse thoroughly with clean water. Or machine wash in hot (60°C) water. Tumble dry low or air dry. DO NOT use bleach.

Foam pad, foam inserts and lumbar support pad

Note: Remove the foam components from the cover prior to cleaning.

To Clean: Wipe lightly with a damp cloth. DO NOT use soap. DO NOT submerge in water. Wipe off with a clean cloth. Completely air dry before replacing in cover.

Warning: Disinfectants are NOT effective on porous surfaces such as foam. Replace soiled or contaminated foam components with authorized ROHO foam components only. Replacement foam components may be purchased separately. Soiled or contaminated foam components MUST NOT be used by different individuals.

Air Inserts

Note: Remove the air inserts from the cover and foam pad prior to cleaning and disinfecting.

Cautions: Make sure the inflation valves are closed. DO NOT allow water or cleaning solution to enter the air insert. DO NOT machine wash or dry. DO NOT expose to direct sunlight. DO NOT autoclave. Rinse thoroughly. Detergent residue may cause air cells to stick together.

To Clean: Close valves. Place in a large sink. Hand wash, using mild liquid hand soap, hand dishwashing liquid, laundry detergent, or a multipurpose detergent (follow product label instructions.) Gently scrub all cushion surfaces, using a soft, plastic bristle brush, a sponge, or a cloth. Rinse with clean water. Air dry.

To Disinfect: Follow the cleaning instructions, using 1 part household liquid bleach per 9 parts water. Keep the product wet with the bleach solution for 10 minutes. Rinse thoroughly with clean water. Air dry.

Note: Deflated air cells may be easier to thoroughly clean. To deflate the air cells: with the inflation valve open, roll up and/or compress the cushion, and then completely close the inflation valve.

Care Symbols

Hand wash.	Machine wash cold, gentle.	Machine wash hot, normal, at the temperature shown.	Do not machine wash.
Dry flat.	Tumble dry low.	Do not bleach.	Do not tumble dry.
Wipe with a damp cloth.	Dry flat in the shade.	Bleach (1 part bleach : 9 parts water).	Line dry.

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Permobil LTD