Operation Manual

Compass® Power Mount Hardware

permobil

Contact Information

The Comfort Companies, LLC
1960 South Calhoun Road
New Berlin, Wisconsin 53151 USA

Phone:
U.S.A.: +1-800-736-0925
Outside the USA: +1-618-277-9150

Fax: +1-406-522-8563

Customer Support:
orders.comfort@permobil.com
permobil.com

WARNING

Indicates that not following the specified procedure could lead to potentially hazardous conditions resulting in serious injury.

CAUTION

Indicates that not following the specified procedure could lead to potentially hazardous conditions resulting in minor to moderate injury or damage to the equipment or other property.

Incident Reporting

In the event of a serious incident or injury, contact Customer Support.

Table of Contents

Intended Use	2
Important Safety Information	2-3
Parts Detail	4
Compatibility and Height	5
Seat Depth	6
Installation - Pre-Assembled	7
Installation - Mounting Bracket and Track Alignment	8
Installation - Quantum, Permobil, Rovi and AVID	9
Installation - Quickie	10
Installation - Back Support Assembly	11
Assembly - Comfort Relief Shell Spacers	12
Assembly - Comfort Relief Shell Bracket and Track Spacer	13
Assembly - Comfort Unite and Comfort Embrace Shell Bracket	14
Assembly - Tru-Shape Shell Bracket	15
Cleaning and Alignment - Tracks	16
Cleaning and Disinfecting	17
Maintenance	
Warranty	
Español	18
Français	35
Português	52

Symbols Glossary

Manufacturer

Model number

Consult instructions for use

Size (width x length x depth)

Unique Device Identifier

Quantity

Medical Device

Please Note

Supplier: These instructions must be given to the user of this product.

Operator (Individual or Caregiver): Before using this product, read the instructions and save for future reference.

Products may be covered by one or more U.S. and foreign patents and trademarks, including BodiLink®, Compass®, Corpus® and Tru-Shape®.

Quickie® is a registered trademark of Sunrise Medical, LLC.

BOA® is a registered trademark of Boa Technology Inc.

Quantum®, Quantum Rival®, Q6 Edge® and TRU-BALANCE® are registered trademarks of Pride Mobility Products Corporation.

Product information is changed as needed; current product information is available at permobil.com. Contact Customer Support for previous versions of product information.

The Comfort Companies, LLC is a part of Permobil.

© 2019, 2025 Permobil

Rev. 2025-01-29

Intended Use

Compass® Power Mount Hardware is intended to attach compatible Comfort back supports to compatible power wheelchairs and has a user weight limit of 450 lbs (204 kg). Compatible Comfort back supports are: Comfort Relief, Comfort Unite, Comfort Embrace and Tru-Shape.

⚠ Compatibility: The Compass Power Mount Hardware is NOT intended for use with other manufacturers' back supports or for attachment to wheelchairs other than those specified in this manual.

The medical products featured in this document are intended to be only one part of an overall care regimen that includes all seating and mobility equipment and therapeutic interventions. A clinician who is experienced in seating and positioning should identify this care regimen after evaluating the individual's physical needs and overall medical condition.

A clinician who is experienced in seating and positioning should determine whether the product is appropriate for the individual's particular seating needs. The product's intended use statement alone should not be relied upon to make that determination.

Important Safety Information

⚠ WARNINGS:

- This product is designed for the intended user and intended use only. Any modifications could potentially endanger the user and cause discomfort or damage.
- Ensure that the product does not interfere with any other wheelchair components. Interference could jeopardize effectiveness and user safety as well as cause damage.
- Avoid sharp objects.
- DO NOT modify this product. Doing so will void its warranty and may lead to personal injury or alter the effectiveness of the product by increasing the risk factors for skin breakdown and/or instability.
- Follow all instructions and safety information provided with the product and provided by the manufacturer of the wheelchair or other accessories.
- Understanding how to care for the product is the key to safe and proper use. Share this information with all caretakers to help them meet healthcare needs and ensure safety.
- Permobil cannot be held responsible for damages or injury due to misuse of the product or failure to follow instructions provided. Please contact us if further assistance is required in understanding cautions, warnings, or instructions for using or maintaining the product.
- Installing alternate parts to the wheelchair may change the intended structure and function of the equipment. Be sure to assess the need for additional safety features for your wheelchair, such as anti-tip bars or other available options in order to maintain stability during use. Use this manual in conjunction with the manual that came with the wheelchair to ensure safety guidelines are being met.
- It is necessary to fully articulate power recline features to ensure there is no possible equipment interference after installation. Do so slowly while watching closely.
- Never use the equipment to maneuver a wheelchair or anything that it is attached to. It is never to be used as a handle, lift, guide, or support for anything other than the intended use.

⚠ WARNINGS:

- Never use the equipment unless all components are fully functional. For information on assistance with any issues, refer to the "Warranty" section.
- To prevent injury, damages, or voided warranty DO NOT make any alterations to the equipment. Modifications in attempt to interface with any other products than those listed in this manual potentially compromise the safety of the user and voids the warranty.
- When possible, DO NOT use during transport. Repeated or excessive force (braking or crash) could potentially change adjustment or cause extra wear, shortening the life of the product. DO NOT use the equipment if it has been in a crash.
- Caregivers: It is necessary to learn from a healthcare provider how to properly adjust the equipment for the user's specific positioning needs and goals.
- DO NOT place anything in between the user and the product. Additional materials may reduce the effectiveness of the product and increase risk factors for skin breakdown and instability.
- ATTENTION! Make sure that the equipment is securely mounted to a wheelchair system and is adjusted properly each time it is in use. Check all hardware for wear, breakage, and tightness weekly. Re-tighten loose screws/bolts following the instructions in this manual.
- DO NOT continue to use the equipment if any part cannot be tightened or is broken.
- IMPORTANT! TIGHTEN SECURELY BEFORE USE: For all instructions in this manual, all fasteners should be tightened securely, as instructed. If a torque wrench is available, tighten to 9.6 Nm (85 lb-in).
- Keep small hardware components out of the reach of small children and any individuals who have a tendency to place inedible objects in their mouths. Swallowing or inhalation may lead to serious injury or death. IMMEDIATELY seek emergency medical assistance.

Important Safety Information, continued**⚠ CAUTIONS:**

- Check packaging and inspect all parts for damage before use. DO NOT USE if damaged. Check regularly for loose components and tighten as necessary. For information on assistance with damaged parts, refer to the “Warranty” section.
- Use caution when installing, adjusting or removing hardware to avoid pinching or trapping fingers in openings.

⚠ WARNINGS - Motor Vehicle Transportation:

- Failure to pay attention to these warnings could result in severe injury to the individual in the wheelchair or to others.
- Whenever possible, transfer out of the wheelchair and into a manufacturer-installed vehicle seat, and use the vehicle's crash-tested occupant restraint system.
- The back support must be correctly and securely installed, used, and maintained according to the instructions and warnings provided by the wheelchair and hardware manufacturers.
- Install the back support on a wheelchair so the individual is seated in an upright position and the wheelchair is facing forward during transportation in a motor vehicle. The back support has only been crash tested on a wheelchair in the forward-facing position per the test information provided in this operation manual.
- The back support must only be used with wheelchair frames that are equipped with securement points specified in ANSI/RESNA WC-4:2017, Section 19 (WC19) or ISO 7176-19, and that have been successfully crash tested using the methods described in Annex A of WC19 or ISO 7176-19 (respectively).
- The wheelchair must be used with an effective wheelchair securement system and a properly positioned, crash-tested pelvic and shoulder-belt restraint, or Wheelchair Tiedown and Occupant Restraint System (WTORS), following the manufacturer's instructions.
- The top of the back support should be positioned as close as possible to the top of the individual's shoulder to provide good support in the event of a crash. The distance between the top of the individual's shoulder and the top of the wheelchair back should not exceed 6½ in. (16.5 cm).
- ALWAYS safely store during transport to avoid damage.
- If the product has been involved in an accident during transport, see the “Maintenance” section.
- For more information about safely transporting wheelchair users, refer to: <https://wc-transportation-safety.umtri.umich.edu/ride-safe-brochure/>.


Back support-hardware combinations have been dynamically tested by a third-party testing facility for use in a motor vehicle and meets all applicable criteria for wheelchair seating systems proposed in Section 5.1 of ANSI/RESNA WC-4:2017 Section 20 “Wheelchair Seating Systems for Use in Motor Vehicles” or ISO 16840-4(2009): Wheelchair Seating - Part 4: Seating Systems for Use in Motor Vehicles.

Test results apply to seating systems intended for use by adults and children with a total body mass greater than 51 lbs (23 kg, or approximately 6 years of age) who are transported facing forward in their wheelchairs in all types of motor vehicles.


To determine whether motor vehicle third-party testing has been performed with a specific back support hardware and other accessory combination, refer to the operation manual provided with back support, attachment hardware, and other accessory, or contact the equipment manufacturer.

Parts Detail

Mounting Bracket

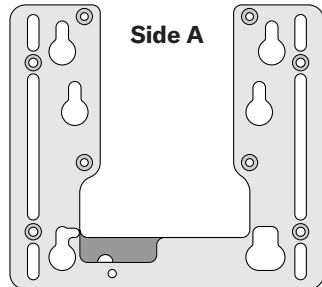
Side A

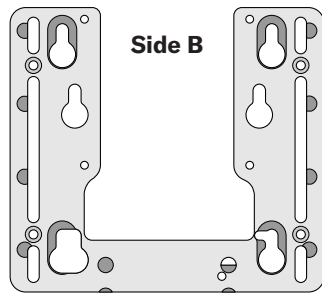Side B

**Shell Bracket
Comfort Embrace**
12 in. (30.5 cm) L -
22 in. (55.9 cm) L

L = Length

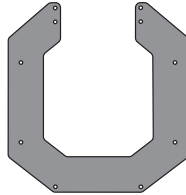
**Shell Bracket
Comfort Unite 16 in. (10.6 cm) L -
22 in. (55.9 cm) L**
Tru-Shape 22 in. (55.9 cm) H -
30 in. (76 cm) H


H = Height


**Shell Bracket
Comfort Unite**
12 in. (30.5 cm) L -
14 in. (35.6 cm) L
Tru-Shape 18 in. (45.7
cm) H - 21 in. (53 cm) H

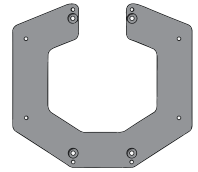
**Shell Bracket
Comfort Relief**
12 in. (30.5 cm) L

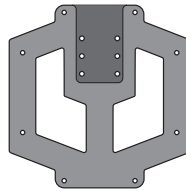**Shell Bracket
Comfort Relief**
14 in. (35.6 cm) L

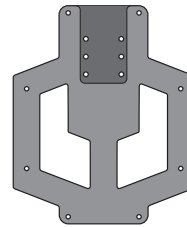
**Shell Bracket
Comfort Relief**
16 in. (40.6 cm) L

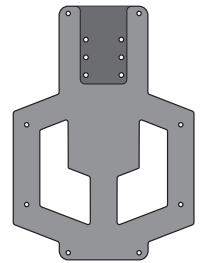
Male Track


Female Track

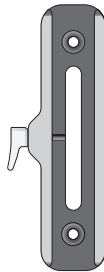
**Shell Bracket
Comfort Relief**
18 in. (45.7 cm) L

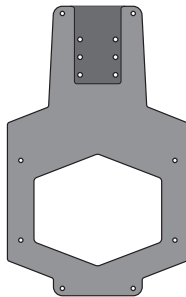
**Shell Bracket
Comfort Relief**
20 in. (50.8 cm) L

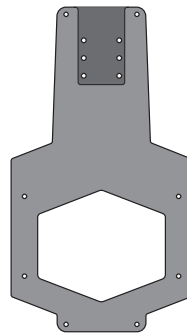**Shell Bracket
Comfort Relief**
22 in. (55.9 cm) L

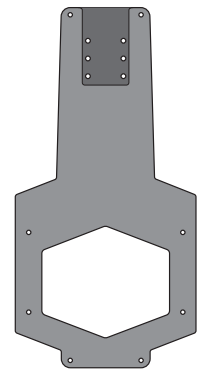Spacer
Comfort Relief 2 in. (5.1 cm) Shell (top)
Comfort Relief 1 in. (2.5 cm) Shell (bottom)

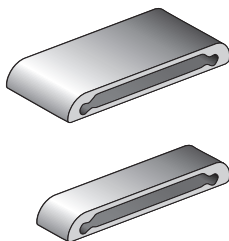Spacer
Comfort Unite
Comfort Embrace
Tru-Shape

Track Spacer
Comfort Relief

Wide Adapter Plate
Quickie

1 in. = 2.54 cm = 25.4 mm

Compatibility and Height

Compass Power Mount Hardware is compatible with the wheelchairs listed below, using only the seating systems specified. The power bases listed are a reference for what is available with the compatible seating systems.

POWER BASES AVAILABLE	Permobil with Corpus 3G Seating System	Quickie with 3.7 Power Recline Seating System	Quantum with Tru-Balance 3 and 4 Seating Systems	Rovi with Motion Concepts Seating System	AVID
	M400 Corpus	QM-710™ SC	R-TRAK	<i>Note: A Recessed Planar Interface Plate must be installed on the wheelchair.</i>	AXCEL Narrow
	M300 Corpus	QM-710 MPC	Edge 3 Stretto		Ultra w/Power Seating
	C350 Corpus	QM-715™ HD MPC	Edge 3		Vision Ultra
	F3 / M3 Corpus	QM-720™ MPD	Q6 Edge 2.0 and Q6 Edge HD	X3	Vector HD
	F5 / M5 Corpus	PULSE™ 6 SC	Q1450		Velocity
		PULSE™ 6 MPC	J4, J4E and J4HD		
			J6		

Each wheelchair manufacturer has a suggested “back gap” for setting the height of the back support on the wheelchair. This allows for the power wheelchair components to articulate as intended. Please see the height guide charts below*.

Note: This section does not apply to the Tru-Shape, for which the back height is measured instead of the “back gap”.

POWER MOUNT HEIGHT GUIDE		Quantum <u>without</u> recline		Quantum <u>with</u> recline	
		Recommended CC Back Support Length		Recommended CC Back Support Length	
BACK HEIGHT	30 in. (76.2 cm)	BACK GAP = 3 in. (7.6 cm)	22 L	BACK GAP = 7 in. (17.8 cm)	ADJUSTMENT = -3 in. (-7.6 cm), +5 in. (+12.7 cm)
	28 in. (71.2 cm)				
	26 in. (66.0 cm)				
	24 in. (60.9 cm)				
	22 in. (55.9 cm)				
	20 in. (50.8 cm)				
	18 in. (45.7 cm)				
	16 in. (40.6 cm)				

*Rovi suggested back gap is 4 in. (10.2 cm) for tilt systems and 6 in. (15.2 cm) for recline systems. Recommended Comfort Company (CC) Back Support Length and Adjustment information depends on where the holes are drilled in the Recessed Planar Interface Plate during installation.


Contact Merits for information on suggested back gap for AVID.

BACK HEIGHT = Overall dimension from the seat base to the top edge of the back support.

BACK GAP = Manufacturer recommended dimension from the seat base to the bottom edge of the back support. Recommended gap is identical for all back supports.

ADJUSTMENT = Available height adjustment range after recommended back gap. (½ in. increments)

POWER MOUNT HEIGHT GUIDE		Permobil		Quickie	
		Recommended CC Back Support Length		Recommended CC Back Support Length	
BACK HEIGHT	30 in. (76.2 cm)	BACK GAP = 6 in. (15.2 cm)	22 L	BACK GAP = 7 in. (17.8 cm)	ADJUSTMENT = -2 in. (-5.1 cm), +6 in. (+15.2 cm)
	28 in. (71.2 cm)				
	26 in. (66.0 cm)				
	24 in. (60.9 cm)				
	22 in. (55.9 cm)				
	20 in. (50.8 cm)				
	18 in. (45.7 cm)				
	16 in. (40.6 cm)				

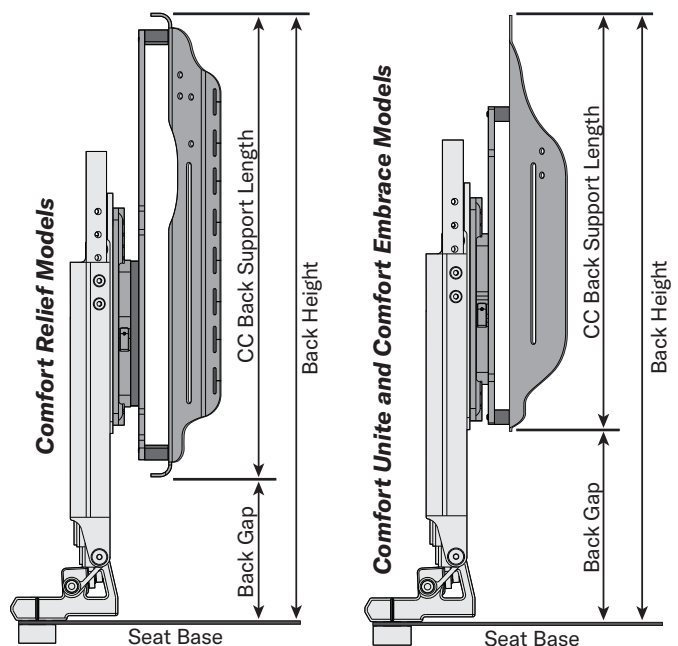Seat Depth

Comfort Relief

SEAT DEPTH USED = the distance from the base mounting plate to the rear face of the back support shell. This is the distance required for the *Compass Power Mount Hardware* assembly, including spacers. Back support foam thickness varies by model and size and therefore is not included in seat depth used.

Comfort Relief requires the use of at least a 1 in. spacer to connect the back support shell and shell bracket. 2 in. of seat depth is used in this configuration which is standard when shipped with a back support. This allows for the user to still immerse through the shell of the Comfort Relief back support effectively reducing the amount of seat depth used from that listed.



Extra spacers may be used to increase the available immersion. The possible shell spacer configurations include 1 in. spacers only, 2 in. spacers only, or 1 in. and 2 in. spacers combined.

When using a head support with Comfort Relief, an extra track spacer may be needed for clearance. This spacer requires an additional .356 in. of seat depth.





Comfort Unite, Comfort Embrace, and Tru-Shape

Comfort Unite and **Comfort Embrace** back supports require the use of four spacers to connect the back support shell to the hardware shell bracket. 1.5 in. of seat depth is used.

See assembly section if attaching Compass Power Mount Hardware to a Comfort Unite or Comfort Embrace back support.

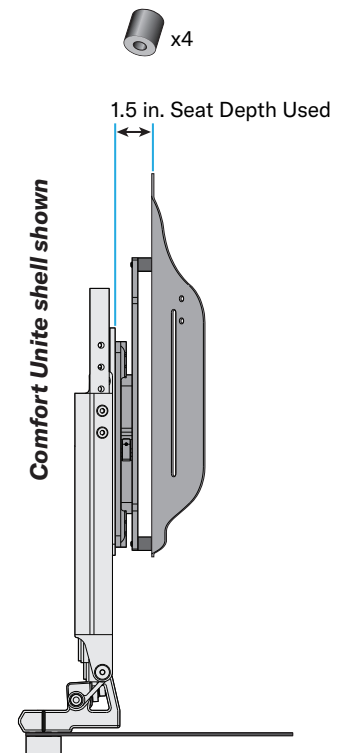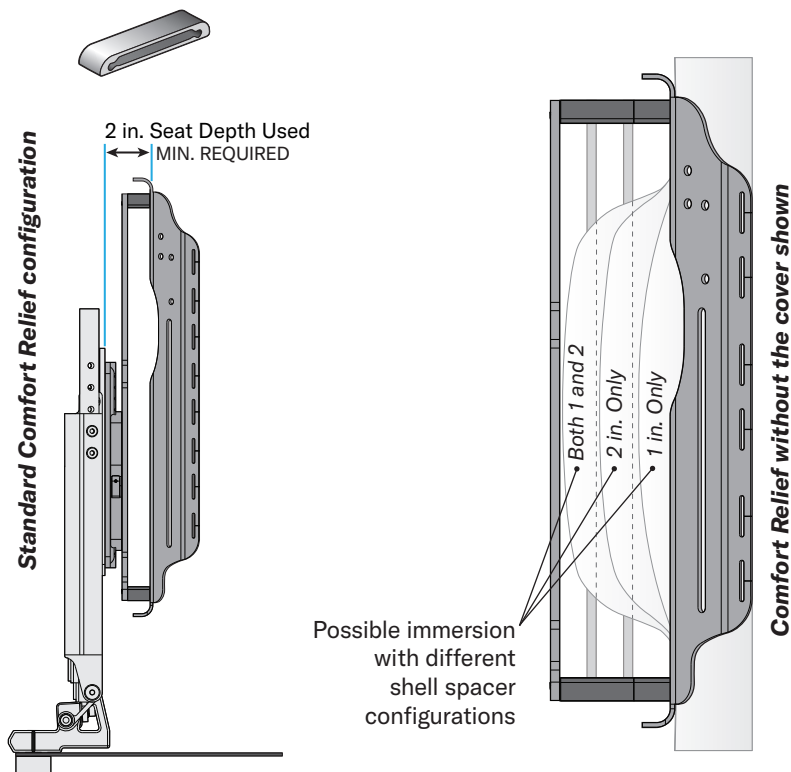Comfort Relief Seat Depth Used (SDU) for all spacer configurations

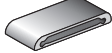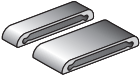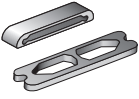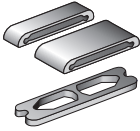2 in. (SDU)	3 in. (SDU)	4 in. (SDU)	See assembly section when using shell spacers.
			
1 in. Spacer	2 in. Spacer	1 + 2 in. Spacers	See assembly section when using track spacers.
2.356 in. (SDU)	3.356 in. (SDU)	4.356 in. (SDU)	
			
1 in. Spacer + Track Spacer	2 in. Spacer + Track Spacer	1 and 2 in. Spacers + Track Spacer	

CONVERSION CHART		
.356 in.	=	0.9 cm
1 in.	=	2.5 cm
1.5 in.	=	3.6 cm
2 in.	=	5.1 cm
2.356 in.	=	6.0 cm
3 in.	=	7.6 cm
3.356 in.	=	8.5 cm
4 in.	=	10.2 cm
4.356 in.	=	11.1 cm
5 in.	=	12.7 cm
6 in.	=	15.2 cm
7 in.	=	17.8 cm

Installation - Pre-Assembled

Use the following instructions if Compass Power Mount Hardware is already assembled and attached to the back support.

Compass Power Mount Hardware is usually provided pre-assembled, attached to the back support, and ready for use with the wheelchair designated during the ordering process. If the assembly is not ready for installation due to any of the reasons shown below, follow the instructions specified before continuing with the installation process.


Issue	Instructions
The hardware assembly is not attached to a back support.	Assembly - Comfort Relief Shell Spacers
	Assembly - Comfort Unite and Comfort Embrace Shell Bracket
The hardware is not configured for the power wheelchair.	Installation - Mounting Bracket and Track Alignment
The hardware assembly requires posterior depth adjustment.	Assembly - Comfort Relief Shell Spacers
	Assembly - Comfort Relief Shell Bracket and Track Spacer
Parts are being replaced.	Installation - Mounting Bracket and Track Alignment
	Assembly - Comfort Relief Shell Bracket and Track Spacer
	Assembly - Comfort Unite and Comfort Embrace Shell Bracket
The hardware has binding issues.	Cleaning and Alignment - Tracks

IMPORTANT! The back support shell is shown without a cover. The cover should be on the back support during installation.

⚠ WARNING: DROPPING RISK: To avoid risk of personal injury or damage to the product, make sure to support the mounting bracket while releasing, as it may slide easily, and will drop if unsupported.

- 1 Flip the release levers (A) out on both sides of the hardware assembly to release the mounting bracket (B).
- 2 Slide the mounting bracket (B) down and set the back support (C) with the remaining hardware aside. See the next page to verify orientation of the mounting bracket (B) and to make sure it is set up for installation on the desired power wheelchair.

Comfort Unite shell shown

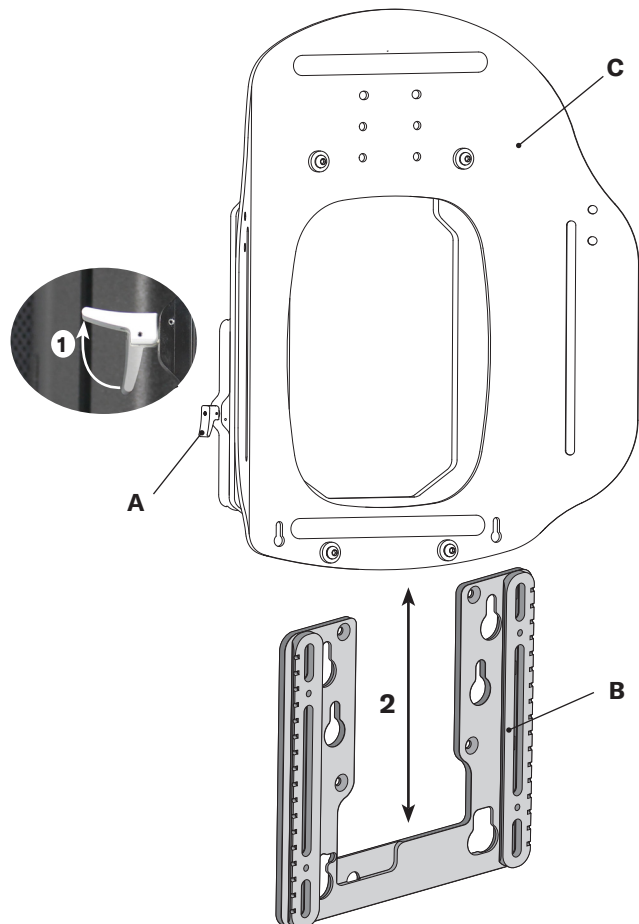Installation - Mounting Bracket and Track Alignment


Verify and/or switch the mounting bracket orientation; or replace the mounting bracket and male track parts.

If the mounting bracket is already set up for the wheelchair, see the Installation instructions on the next two pages.

If the mounting bracket orientation needs to be switched, follow steps a - f on this page.

The mounting bracket should be set up for the power wheelchair that was selected when the product was ordered.

To determine whether to use **SIDE A** or **SIDE B** of the mounting bracket, look for the location of the male track pieces that have notches, as shown here.

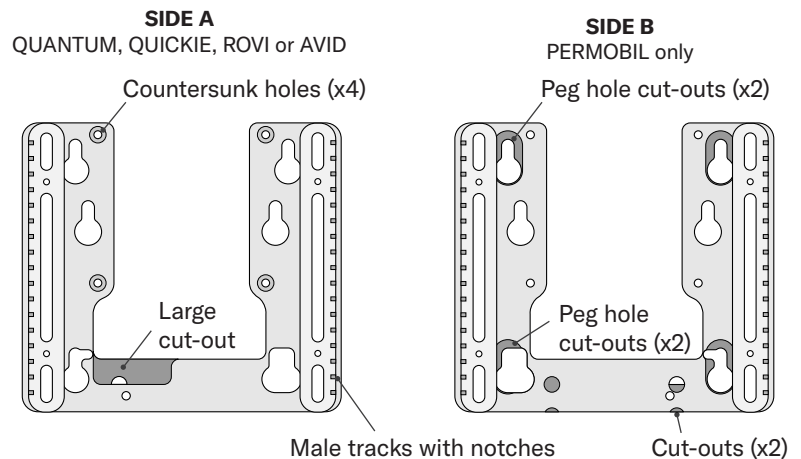
To adjust the mounting bracket:

- From the back side of the mounting bracket (B), unscrew the four flat head screws (E) with the hex key.
- Move the male tracks (D) to the opposite sides of the mounting bracket.

Note: The notches on the male tracks must be visible and run along the outer edge of the mounting bracket.

- Reassemble the mounting bracket (B) and the male tracks (D) using the screws (E).

IMPORTANT! The screws should be evenly tightened for the best alignment.

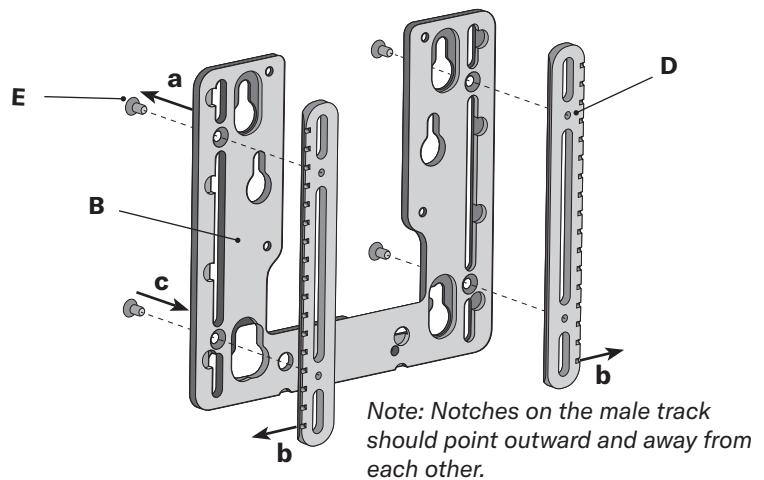
To check track alignment:

- With both release levers (A) in the open position, slide the male tracks into the female tracks of the hardware assembly. It should slide up and down easily.
- If the hardware assembly does not slide easily, locate the eight screws on the rear face of the track assembly. Gently loosen all eight screws.
- Tighten screws 1 through 8 in the order shown.

Repeat steps until the assembly slides easily, and then continue with the installation instructions on the next two pages.


If needed, see the instructions on cleaning and lubricating track pieces, see "Cleaning and Alignment - Tracks".

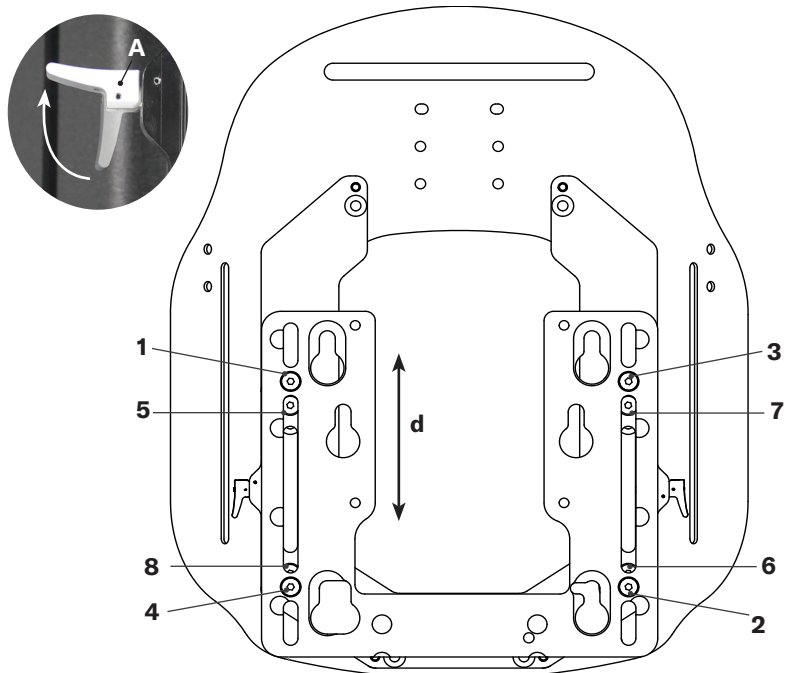Installation - Quantum, Permobil, Rovi and AVID

Installing the base mounting plate to the different brands of wheelchairs.


Quantum (Tools are not required for installation.)

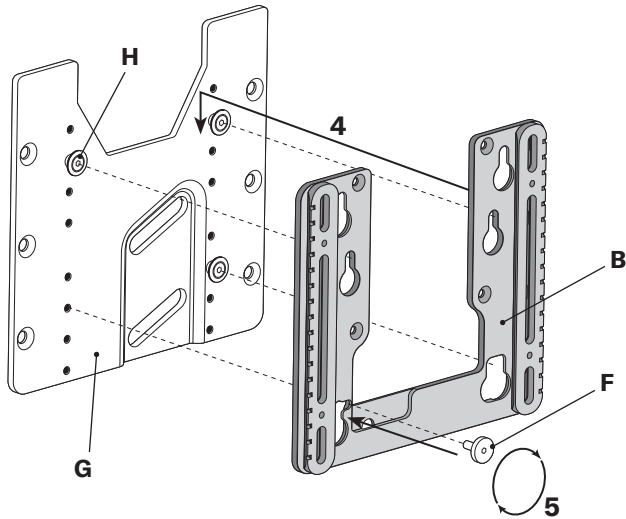- 3 Remove the threaded knob (F) from the base mounting plate (G).
- 4 Line up the three pegs (H) with the mounting bracket (B) and slide it down on the pegs.
- 5 Insert the threaded knob (F) and tighten to secure.

Continue with "Installation - Back Support Assembly," step 7.


Permobil (Tools are not required for installation.)

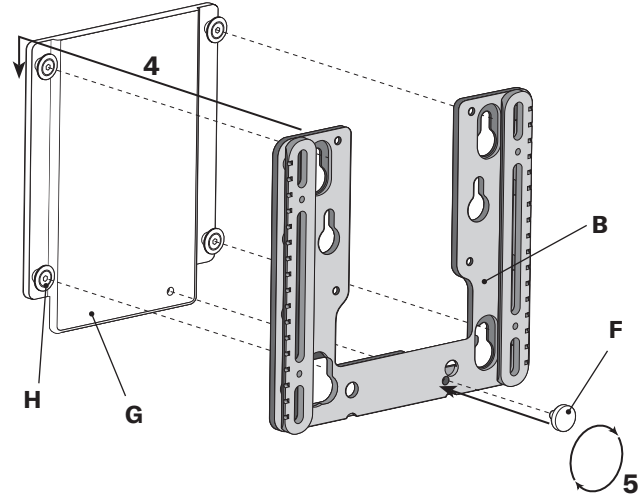- 3 Remove the threaded knob (F) from the base mounting plate (G).
- 4 Line up the four pegs (H) with the mounting bracket (B) and slide it down on the pegs.
- 5 Insert the threaded knob (F) and tighten to secure.

Continue with "Installation - Back Support Assembly," step 7.

Rovi (10mm wrench and drill with ¼ in. drill bit are required for installation. Hex key is provided.)

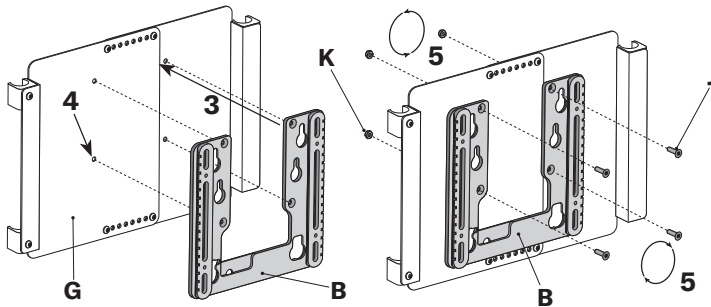
Note: A Recessed Planar Interface Plate must already be installed on the wheelchair and will be referred to as a "base mounting plate" in these instructions.

If needed, refer to the instructions provided with wheelchair for installation of the Recessed Planar Interface Plate.

- 3 Line up the back support against the base mounting plate (G) where desired, and mark the location of the four countersunk holes on the mounting bracket (B).
- 4 Ensure that the holes are level, and then drill four ¼ in. holes in the base mounting plate (G) where marked.
- 5 Use the provided longer flat head screws (J) and the lock nuts (K), and fasten securely with the hex key and the 10mm wrench.

Continue with "Installation - Back Support Assembly," step 7.

AVID (The provided hex key is required for installation.)

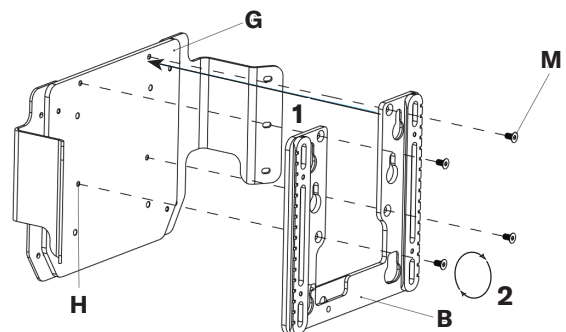
- 3 Align the four countersunk holes on mounting bracket (B) with threaded holes (H) on base mounting plate (G) as shown.
- 4 Use the provided 12mm flat head screws (M) and fasten with the hex key.

Continue with "Installation - Back Support Assembly," step 7.

Installation - Quickie


Installing the base mounting plate to the wheelchair.

Quickie - (The provided hex key is required for installation.)

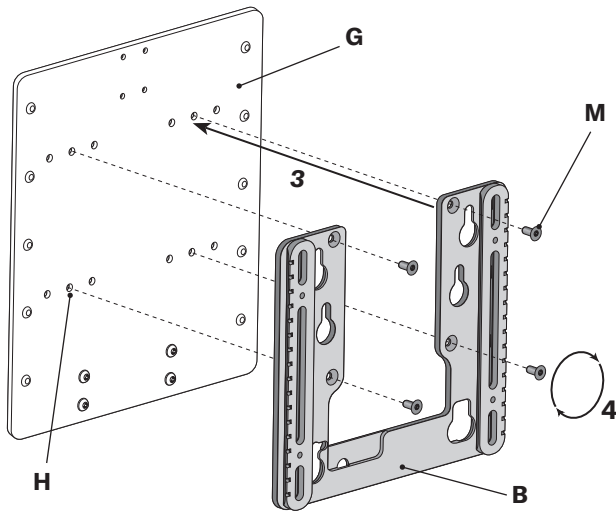- 3** Align the four countersunk holes on the mounting bracket (B) with the threaded holes (H) on the base mounting plate (G), as shown.

Note: Only the center threaded holes (H) on the base mounting plate (G) will correctly align with the mounting bracket (B).


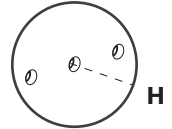
- 4** Use the provided 12mm flat head screws (M), and fasten securely with the hex key.

Continue with "Installation - Back Support Assembly," step 7.

Quickie 20W and up - (The provided hex key is required for installation.)

**Note: Quickie Wide Adapter Plates use an additional ¼ in. (6.4 mm) of seat depth.*

Contact Customer Support for Quickie Wide Adapter Plate and fasteners, if needed.

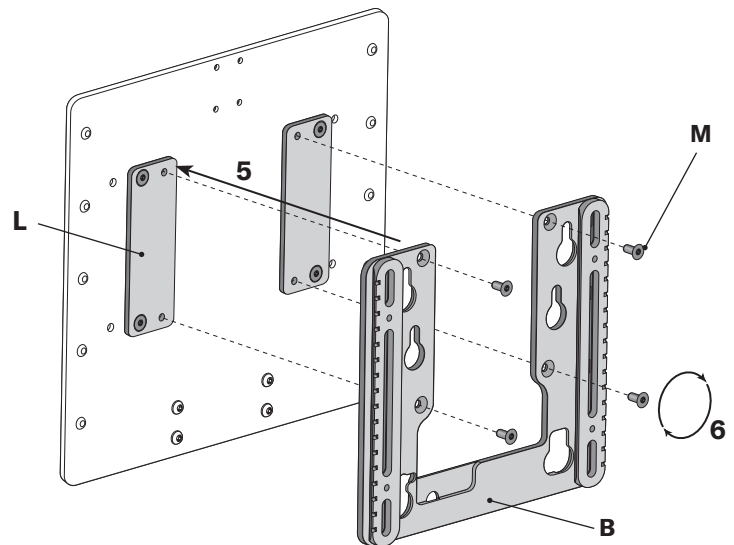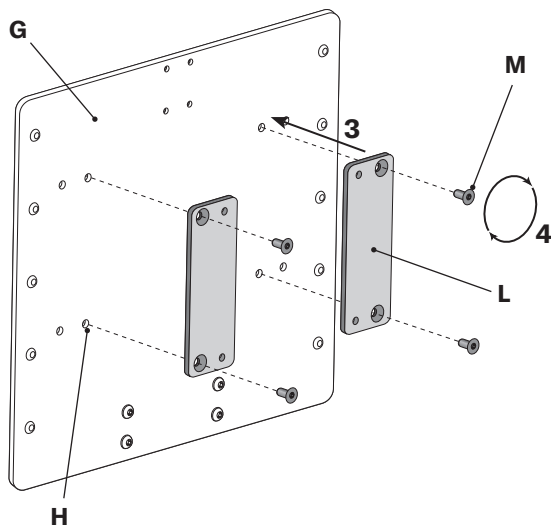- 3** Align the two countersunk holes on each of the adapter plates (L) with the inside threaded holes (H) on the base mounting plate (G), as shown.

- 4** Attach, using the four 12mm flat head screws (M) and the hex key.

Note: The threaded holes on the adapter plates should face in toward each other.

- 5** Align the four countersunk holes on the mounting bracket (B) with the threaded holes on the adapter plates (L), as shown.

- 6** Use the four additional 12mm flat head screws (M), and fasten securely with the hex key.

Continue with "Installation - Back Support Assembly," step 7.

Installation - Back Support Assembly

⚠ WARNING: DO NOT allow wires or other components to run through the support/hardware system at any time.

- 7** From behind the wheelchair, hold the back support and hardware assembly on each side at the release levers (A) with the release levers in the open position, as shown at right.

Note: Holding the back support this way allows it to be lowered evenly for better alignment.

- 8** Line up the female tracks (P) with the male tracks (not shown) on the mounting bracket installed on the wheelchair, and slide the assembly down.
- 9** Once the back support is close to the desired position, lower the release levers until they click and are completely engaged.

IMPORTANT! If a lever is not fully engaged (shown in the center image below), move the back support assembly slightly up or down until the lever clicks into place and is completely engaged (shown in the right-hand image below).


Not fully engaged

Completely engaged

Comfort Embrace without the cover shown (rear view)

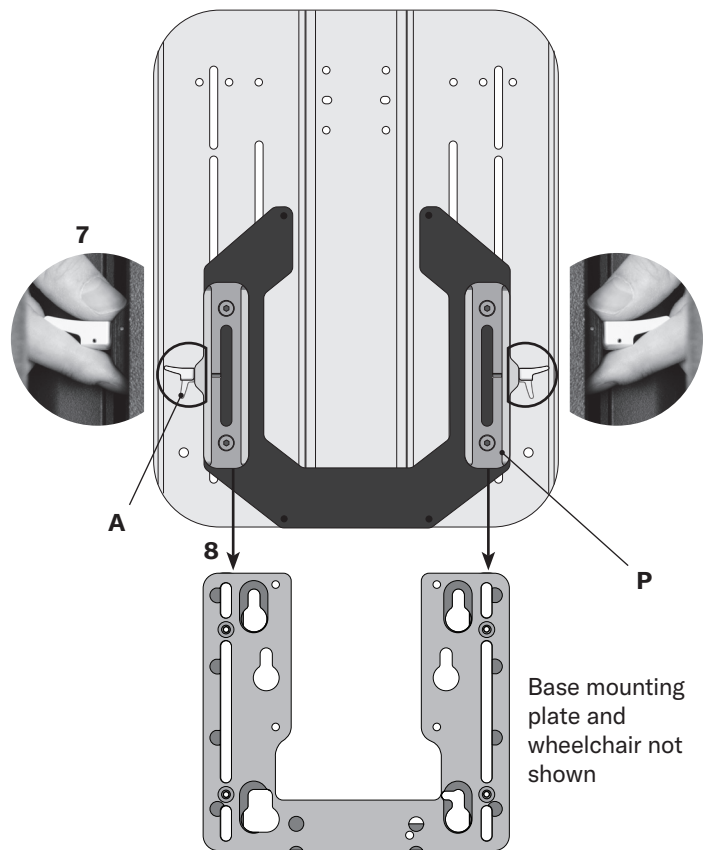If the tracks appear to be misaligned or if binding occurs, see “Cleaning and Alignment - Tracks”.

Setting the back support height.

Note: Each wheelchair manufacturer has a suggested “back gap” for setting the height of the back support on the wheelchair. Refer to the height recommendations and then evaluate the patient needs to set a final back support height.

Note: The patient may be seated during step 10.

- 10** Flip the release levers open and then slide the back support up or down to the desired height. Close the levers.


Before continuing, check the full range of articulation on the wheelchair. See the warning below.

⚠ WARNING: To prevent injury or damage, the full range of articulation on the wheelchair **MUST** be checked for any possible interference or unwanted movement of seating components.

- 11** To lock the release levers in the closed position: Using the hex key, turn the set screw on the release levers until it makes contact.

CAUTION: DO NOT remove the set screw. When the release levers are not locked, make sure the set screw is flush with the release lever surface.

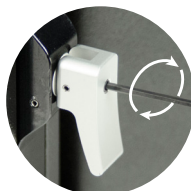Recommended Back Gap




The Back Gap is the manufacturer-recommended dimension from the seat base to the bottom edge of the back support.

For more information, see “Compatibility and Height”.

Wheelchair Model	Recommended Back Gap
Quantum without recline	3 in. (7.6 cm)
Quantum with recline	7 in. (17.8 cm)
Quickie	7 in. (17.8 cm)
Permobil	6 in. (15.2 cm)
Rovi Recline System	6 in. (15.2 cm)
Rovi Tilt System	4 in. (10.2 cm)
AVID	Contact Merits for dimensions

Assembly - Comfort Relief Shell Spacers

SDU = Seat Depth Used

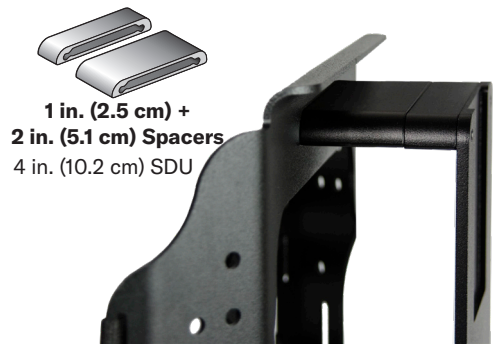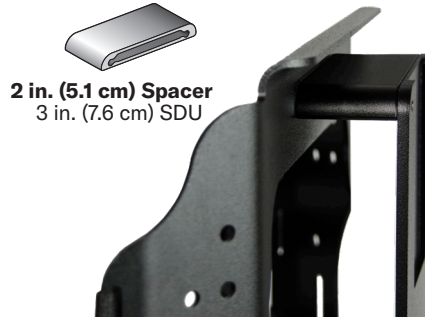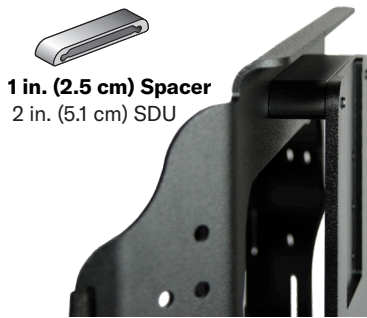Standard Comfort Relief configuration

Attaching the shell bracket and spacers to the back support:


1 Remove the cover from the back support and locate the four holes on the shell (D) for attaching the power mount hardware. If already assembled, unscrew the four bolts (A) to use as alternate spacers.

2 Line up the bolt (A), the safety washer (B), the flat washer (C), the back support shell (D), and the spacer (E) with the shell bracket (F).

Note: The bolt (A) and the spacer (E) will vary, depending on the desired configuration. Use the longer bolts with deeper spacers.

3 Lightly fasten and continue with the other three holes using the hex key. After all points are attached, tighten securely.


4 Replace the cover on the back support.

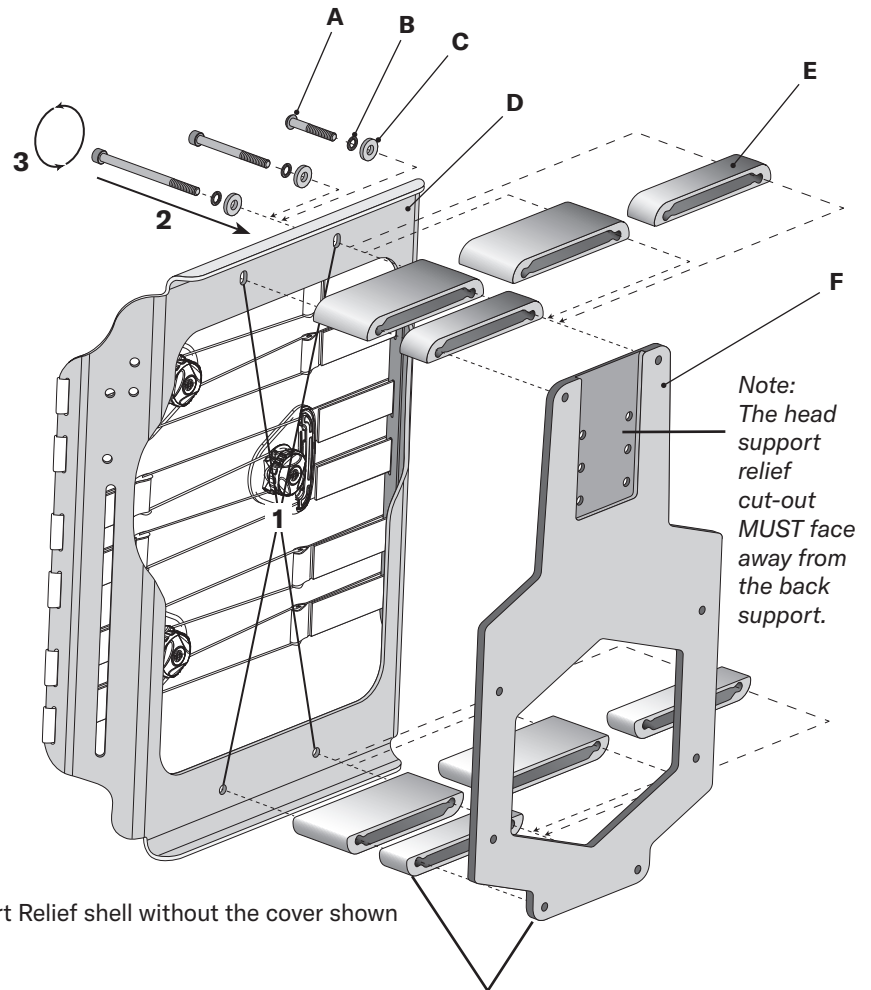
Comfort Relief shell without the cover shown

40 mm button head bolt (A)
1 in. (2.5 cm) shell spacer (E)

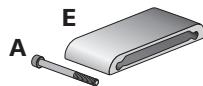65 mm socket head bolt (A)
2 in. (5.1 cm) shell spacer (E)

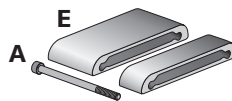90 mm socket head bolt (A)
1 in. (2.5 cm) and 2 in. (5.1 cm)
shell spacers (E)


Assembly - Comfort Relief Shell Bracket and Track Spacer

Track spacers for Compass Power Mount Hardware are to be used on Comfort Relief only when using a head support and when additional clearance is needed. Refer to the instructions provided with the head support for installation of the head support and hardware.

Attaching the female tracks to the shell bracket:

- 1 Line up the holes on the shell bracket (A) with the female track (B). Make sure that the release levers (C) point out and away from each other.
- 2 Fasten securely, using the screws (D). Make sure that the screws are evenly tightened for correct alignment.


WARNING: After installation is complete, check for alignment with the mounting bracket. See "Cleaning and Alignment - Tracks".

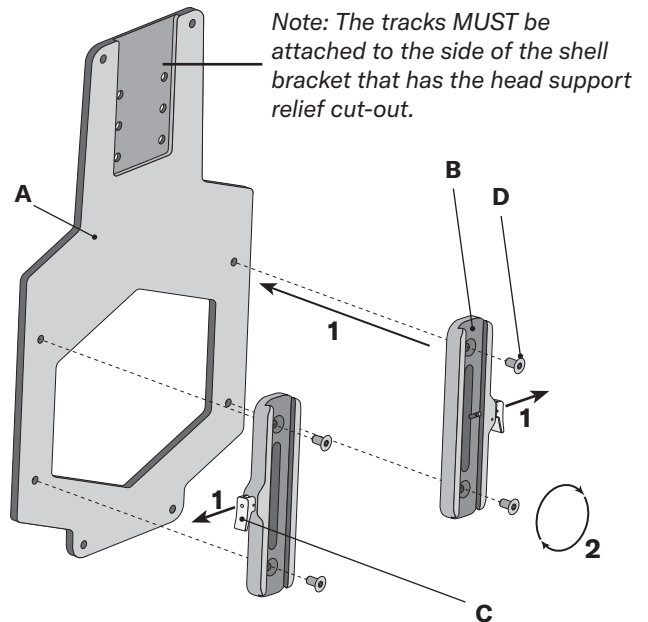
*Adding track spacers to the shell bracket:

- 1 Line up the holes on the shell bracket (A) with the track spacer (E) and the female track (B). Make sure that the release levers (C) point out and away from each other.
- 2 Fasten securely, using the alternate flat head screws (F) provided. Make sure that the screws are evenly tightened for correct alignment.

WARNING: After installation is complete, check for alignment with the mounting bracket. See "Cleaning and Alignment - Tracks".

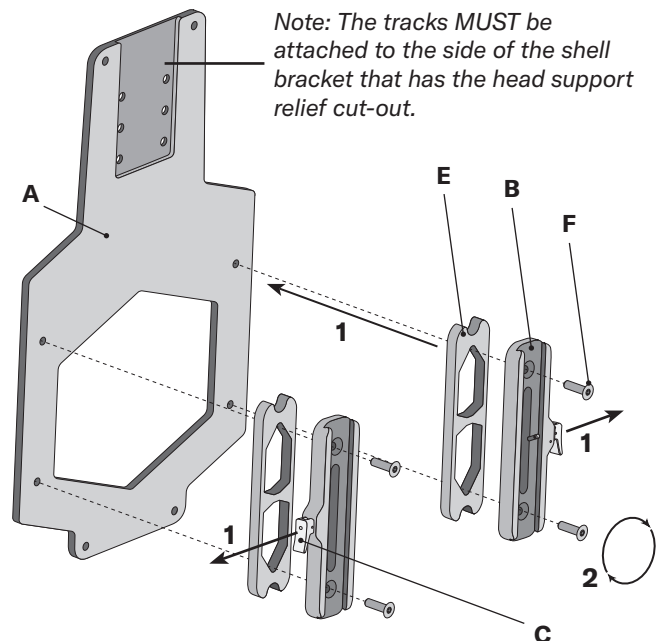Comfort Relief with BodiLink Head Support

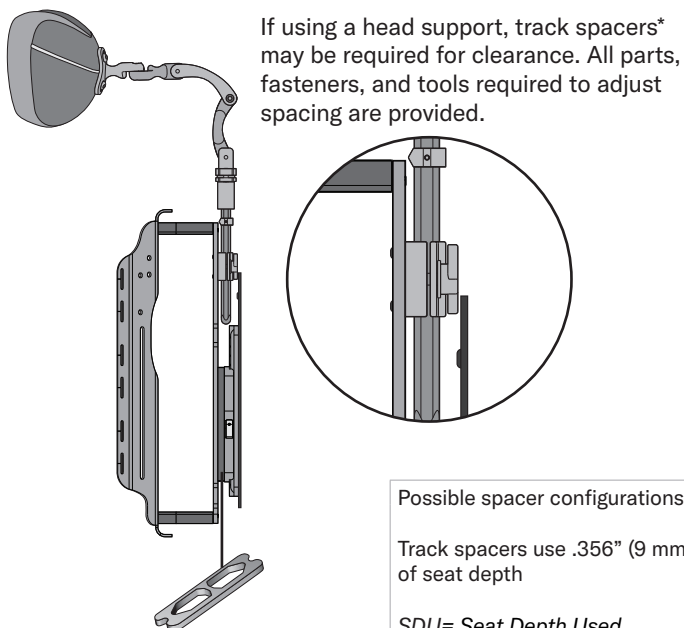Possible spacer configurations:

Track spacers use .356" (9 mm) of seat depth


SDU= Seat Depth Used

2.356 in. (6.0 cm) SDU


1 in. (2.5 cm) Spacer
+ Track Spacer

3.356 in. (8.5 cm) SDU

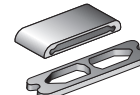2 in. (5.1 cm) Spacer
+ Track Spacer

4.356 in. (11.0 cm) SDU

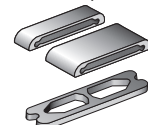1 in. (2.5 cm) and
2 in. (5.1 cm) Spacers
+ Track Spacer


Assembly - Comfort Unite and Comfort Embrace Shell Bracket

Follow these instructions only if the hardware is not already assembled and attached to the back support.

Attaching the female tracks to the shell bracket:



- 1 Line up the holes on the shell bracket (A) with the female track (B). Make sure that the release levers (C) point out and away from each other.
- 2 Fasten securely, using the screws (D). Make sure that the screws are evenly tightened for correct alignment.

WARNING: After installation is complete, check for alignment with the mounting bracket. See “Cleaning and Alignment - Tracks”.

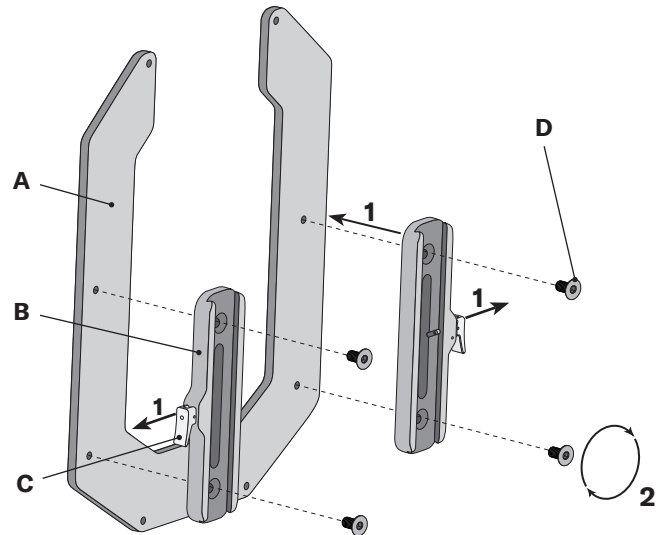Attaching the shell bracket to the back support:

- 1 Remove the cover and locate the four holes on the back support shell (E) for attaching the power mount hardware.
- Note: On Comfort Unite and Comfort Embrace shells, the two bottom holes are round and the two top holes are oblong.*
- 2 Line up the bolt (F), the safety washer (G), the flat washer (H), the back support shell (E), and the spacer (J) with the shell bracket (A).
 - 3 Lightly fasten, using the hex key, and continue with the other three holes. After attaching at all points, tighten securely.

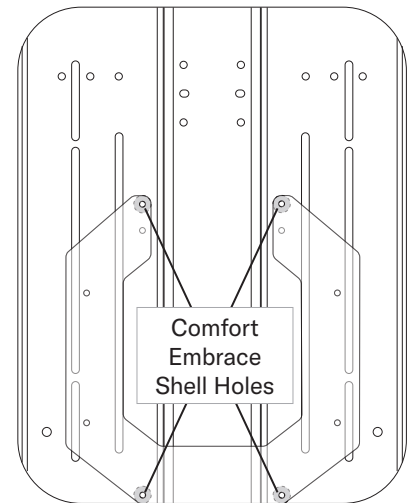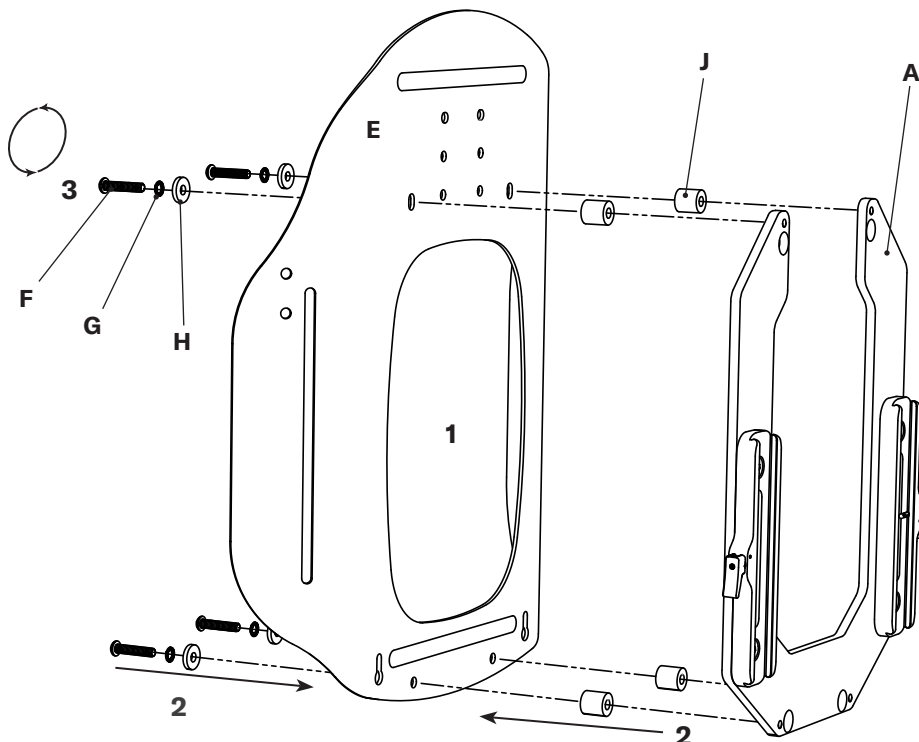Rear view Comfort Embrace shell

Note: For 12 in. (30.5 cm) L Comfort Embrace Shells (not shown) and Shell Brackets use the lower set of holes at the top of the shell and shell bracket for attachment to avoid the cover not fitting correctly.

Comfort Unite shell without the cover shown

Assembly - Tru-Shape Shell Bracket

The hardware may already be attached to the back support. Follow these instructions if the hardware is not already attached to the back support or if reattachment is required.


Notes about installation:

- The Tru-Shape back support wood plate includes multiple sets of T-nut patterns. It may be necessary to use different T-nuts, depending on the back height of the Tru-Shape (especially on Quickie chairs). After installation, confirm that sufficient back adjustment in the track assembly is available.
- The shell bracket for Tru-Shape is available in two sizes. The small shell bracket uses four screws; the large shell bracket uses six screws. The required screws are provided.

- 1** Remove the female tracks (B): On each side, remove the two screws (A) that attach the female track (B) to the shell bracket (C).
- 2** Attach the shell bracket to the back support wood plate (D):
 - a. Locate the countersunk holes (four or six) on the shell bracket (C).
 - b. Line up the screw (E), the shell bracket (C), and the spacer (F) with the back support wood plate (D).
 - c. Use the hex key to lightly fasten all screws. After all screws are engaged, tighten securely.
- 3** Attach the female tracks to the shell bracket (C):
 - a. On both sides: Line up the holes on the shell bracket (C) with the female track (B). Make sure that the release levers (G) point out and away from each other.
 - b. Using the screws (A), attach the female tracks to the shell bracket. Make sure that the screws are evenly tightened for correct alignment.

WARNING: After installation is complete, check for alignment with the mounting bracket. See “Cleaning and Alignment - Tracks”.

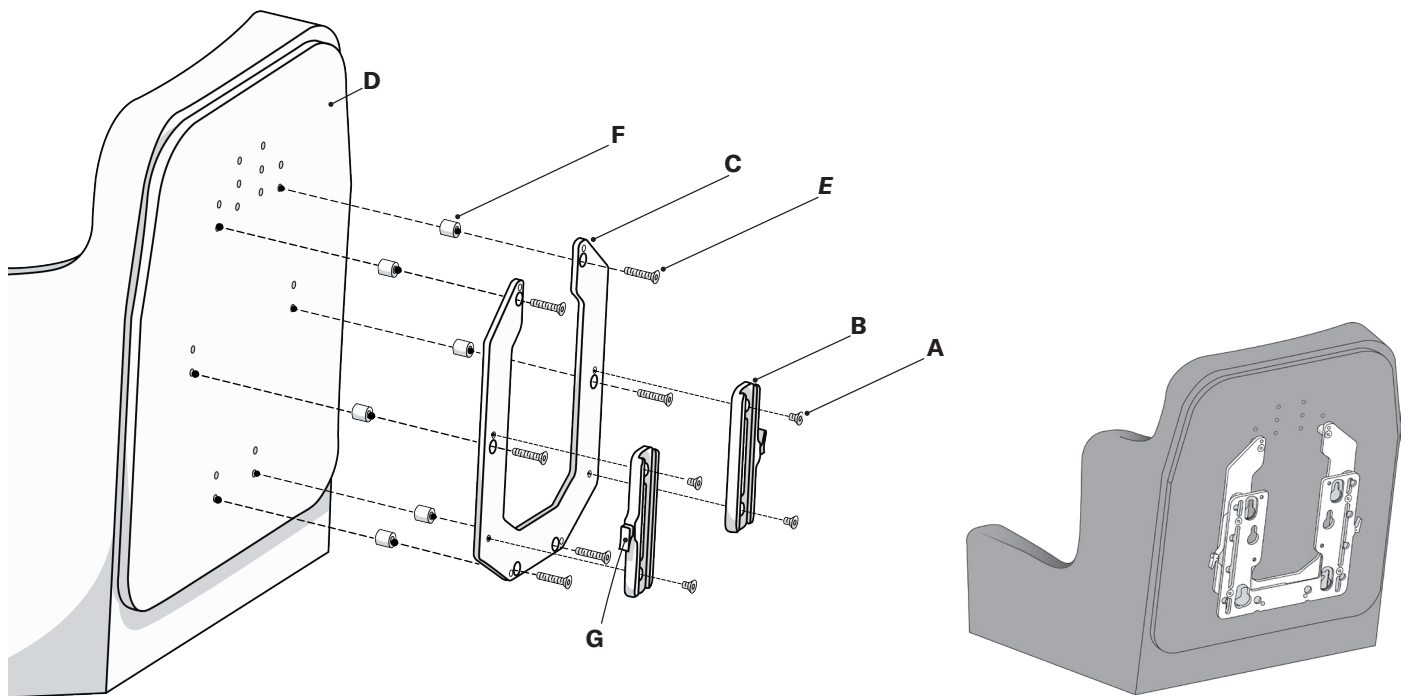Removing Spacers for Additional Seat Depth:


If additional seat depth is needed, attach the hardware to the Tru-Shape Back without using the spacers. Use the $\frac{3}{4}$ in. flathead screws instead of the $1\frac{1}{2}$ in. flathead screws.

WARNING: Removing the spacers from the attachment assembly may result in interference, which may jeopardize effectiveness and user safety and cause damage. Before use, check for component / wheelchair interference.

$\frac{1}{4}$ -20 x $\frac{3}{4}$ in.

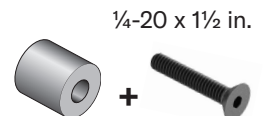$\frac{1}{4}$ -20 x $1\frac{1}{2}$ in.

Spacers add $\frac{5}{8}$ in. seat depth.


Cleaning and Alignment - Tracks

Frequent sliding of the track pieces may cause wear and binding and may require cleaning and lubrication. Follow the track alignment instructions when track pieces are assembled, altered, or appear to be binding.

Cleaning and lubricating tracks:


Note: The following may be completed with the hardware mounting bracket installed on the wheelchair.

- 1 Open the release levers (A) and slide the back support (B) apart from the mounting bracket (C), as shown below.
- 2 Use a clean, dry, cloth to wipe away debris along the outside of the male tracks (D) and the inside female tracks (E).
- 3 Using a cloth or cotton swab, apply a light lubricant to the inside edges of the female tracks (E).
- 4 Slide the mounting bracket (C) and the back support assembly (B) together, and run both tracks together the full length, allowing the lubricant to spread.
- 5 Reset the back support height and close the release levers (A).

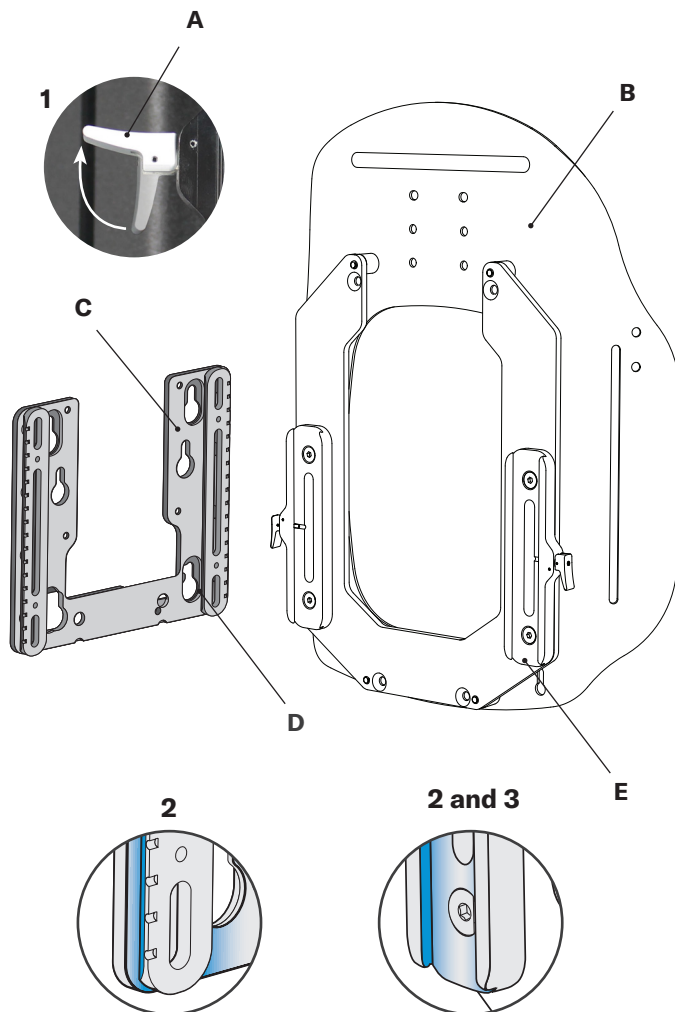Track alignment:

- 1 Start with the back support and hardware assembly removed from the wheelchair.
- 2 From the rear face of the hardware, gently loosen the eight screws in the track assembly.
- 3 For the best results, tighten the screws 1 through 8 in the order shown below.
- 4 To check for track alignment: With the release levers in the open position, slide the tracks up and down.
- 5 Repeat steps 2-4 if needed until the assembly slides easily.
- 6 Install the hardware on the wheelchair. See "Installation - Pre-assembled".

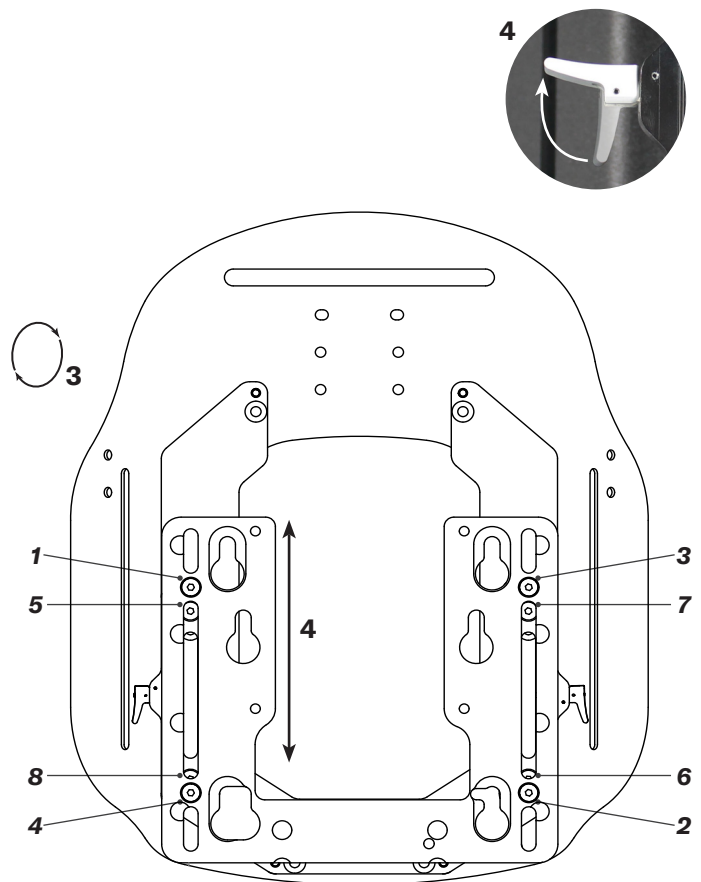Comfort Unite shell without the cover shown

Cleaning and Disinfecting


WARNINGS:

- Follow all manufacturer bleach container and germicidal disinfectant safety guidelines.
- Before transferring to a new user, the product must be completely cleaned and disinfected.
- Periodically, check the components for damage and wear, and replace if necessary. Some wear with regular use is expected.

CAUTION: Before cleaning or disinfecting, remove the cover support pad from the hardware.

Cleaning: Hand wash in warm (40°C/104°F) water with soap and a sponge. Rinse with clean water and wipe dry with a clean cloth. Dry completely before use.

Disinfecting: Wipe gently with a cloth dampened with a household disinfectant. Rinse with clean, warm water. Wipe dry with a clean cloth. Dry completely before use.

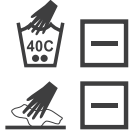Care Symbols

 Wipe with a damp cloth

 Hand wash warm

 Dry flat

Maintenance

Inspection

Inspect all components of the equipment weekly.

Check and tighten all fasteners and adjustment points.

Inspect all materials for wear such as stretching, fraying, fractures, and bends. Some wear with regular use is expected.

Components that show excessive wear should be replaced. For information on replacements, refer to the “Warranty” section.

Reuse

Re-fit the support to the new user or after extended periods of non-use by the same user. Refer to the adjustment and use instructions.

Before transferring to a new user, the product must be completely cleaned and disinfected. Refer to the cleaning and disinfecting instructions.

 **WARNING:** The Tru-Shape back support is designed for use by one person only. Re-use is NOT permitted.

Storage

Assembly or individual parts should not be stored in extreme hot or cold temperatures. Only store in a clean and dry place at room temperature.

Do not allow exposure to direct sunlight for prolonged periods of time. This could age the product more rapidly and decrease its effectiveness.

Do not store where there is risk of open flame or spark.

Disposal

Dispose of the product and/or components in accordance with the applicable regulations in your jurisdiction.

 **CAUTION!** DO NOT INCINERATE as a means of disposal.

Warranty

This product is covered by a 5 year warranty against manufacturer defects.

If a manufacturing defect should occur, discontinue use immediately. Wear from regular use, misuse, damage by an accident or damage from “an act of Nature” (e.g., flood, tornado, earthquake and/or fire) is not considered a manufacturing defect. Replacement parts are available by contacting your supplier, distributor, or retailer.

The following will void the warranty: modifications made in attempt to interface with other products (where compatibility is limited, as described in product information); alterations made to the product; or defects caused by irregular use.

Warranty Replacement Procedure

Warranty claims should be initiated by the original purchaser, who should contact the Durable Medical Equipment (DME) supplier, distributor or retailer from whom the product was purchased, if they are available, or any authorized Permobil DME supplier, distributor or retailer.

Whenever possible, the DME supplier, distributor or retailer with a warranty issue should contact Customer Support and provide the original purchase order number, sales order number, or invoice number.

Before a warranty replacement order or quote may be processed, Customer Support must issue a Return Merchandise Authorization (RMA).

Replacement Parts

Contact Customer Support for information on replacement or individual parts available. Appropriate tools and fasteners will be provided with extra and replacement parts.

Información de Contact

The Comfort Companies, LLC
1960 South Calhoun Road
New Berlin, Wisconsin 53151 EE. UU.

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