
Operation Manual

Permobil Direct Mount Hardware

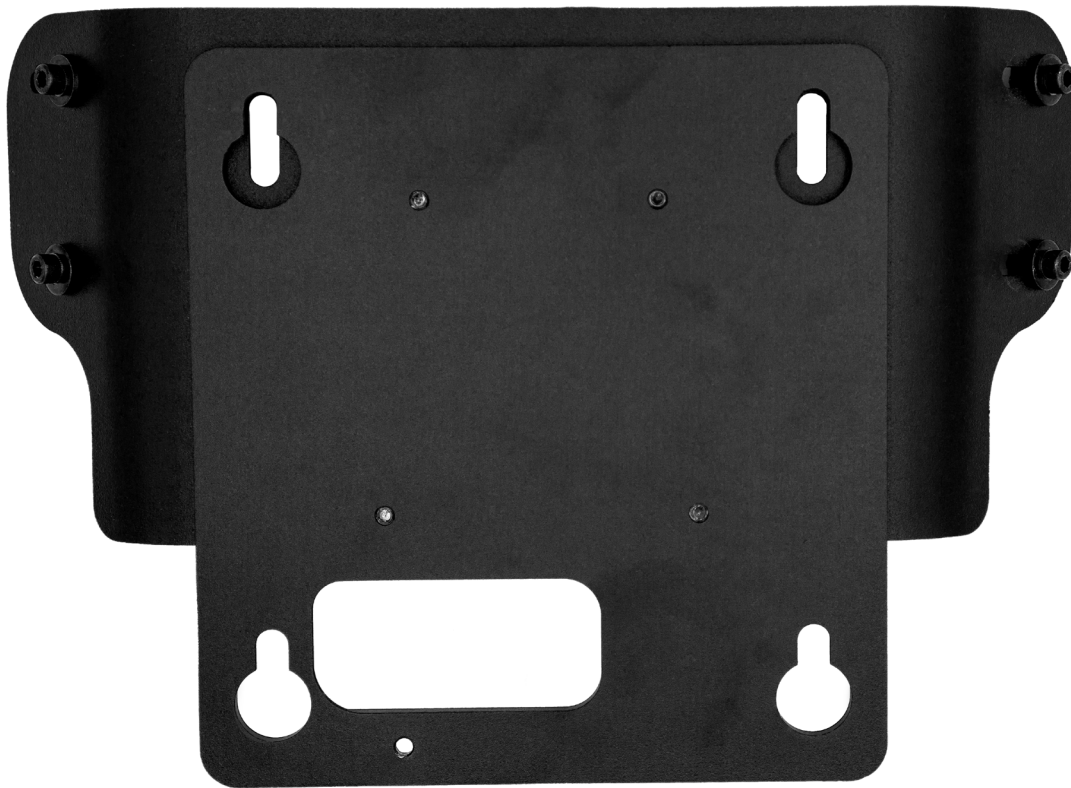permobil

Contact Information

 The Comfort Companies, LLC, 1960 South Calhoun Road
New Berlin, Wisconsin 53151 USA

Phone:

U.S.A.: +1 800-736-0925

Outside the USA: +1 618-277-9150

Fax: +1 406-522-8563

Customer Support:

orders.comfort@permobil.com

permobil.com

⚠ WARNING

Indicates that not following the specified procedure could lead to potentially hazardous conditions resulting in serious injury.

⚠ CAUTION

Indicates that not following the specified procedure could lead to potentially hazardous conditions resulting in minor to moderate injury or damage to the equipment or other property.

Incident Reporting

In the event of a serious incident or injury, contact Customer Support.

Before Use

  This manual references the package contents shown in "Parts Detail." For additional parts information, including information about the back shell, cover and cover components, refer to the instructions provided with the back support.

The images in this manual may not show the specific wheelchair model. The images are for illustration purposes only. For additional information, contact the equipment provider.

Contents

Intended Use	2
Parts Detail and Product Specifications	
Important Safety Information	2-3
Installation On the Wheelchair	4
Height Adjustment	5
Removal from the Wheelchair	
Separately Purchased Back Support:	6
Installation of Direct Mount	
Direct Mount Removal	
Troubleshooting	
Maintenance, Cleaning and Disinfecting	7
Storage, Transport, Disposal and Recycling	
Warranty	
Español	8
Français	15
Português	22

Symbols Glossary

 Manufacturer	 Consult instructions for use	 Date of manufacture	 Model number
 Batch code	 Medical Device	 Quantity	 Unique Device Identifier

Please Note

Supplier: This manual must be given to the user of this product.

Operator (Individual or Caregiver): Before using this product, read the instructions and save for future reference.

Products may be covered by one or more U.S. and foreign patents and trademarks.

Product information is changed as needed; current product information is available at permobil.com. Contact Customer Support for previous versions of product information.

The Comfort Companies, LLC is part of Permobil.

© 2024 Permobil

Rev. 2024-06-11

Intended Use

The Permobil Direct Mount Hardware (Direct Mount) is intended to attach compatible back supports to Permobil power wheelchairs that have compatible Corpus seating. Direct Mount is only compatible with Comfort Unite Back Supports and has a user weight limit of 450 lbs (204 kg).

⚠ Compatibility: The Direct Mount is NOT intended for use with other manufacturer's back supports or for attachment to wheelchairs other than those specified in the intended use statement. The Direct Mount is NOT intended for use with second generation or older Corpus seating.

Wheelchair Configuration: A Permobil-authorized technician and/or a clinician should be consulted to evaluate compatibility, potential limitations or possible interference concerns when using the Direct Mount and the wheelchair configuration.

An experienced clinician should determine whether the product is appropriate for the individual's particular equipment needs. The product's intended use statement alone should not be relied upon to make that determination.

The medical products featured in this document are intended to be only one part of an overall care regimen that includes all therapeutic equipment and interventions. An experienced clinician should identify this care regimen after evaluating the individual's physical needs and overall medical condition.


A clinician should also evaluate for visual, reading, and cognitive impairment to determine the need for caregiver assistance or the need for other assistive technology, such as large-print instructions, to ensure proper product use.

Parts Detail and Product Specifications

Package contents: Direct Mount, components to attach the Direct Mount, UniTrack (for head support attachment; not included with all models)*, operation manual

To order replacement parts, contact the product provider or Customer Support.

* Refer to the instructions provided with the UniTrack. 

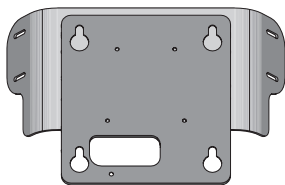Wing bracket / Center bracket

Direct Mount weight range:
3.5 - 4 lb. (1.6 - 1.8 kg)

Expected product lifetime:
5 years

Materials:

Direct Mount wing bracket/center bracket: aluminum

Sizing:

Refer to labeling to confirm that the width of the back shell is compatible with the width of the Direct Mount.

Important Safety Information

⚠ WARNINGS:

- DO NOT use the product without first reading, understanding and following all of the instructions and safety information provided with the product, the Permobil power wheelchair, and any other accessories. Failure to do so could lead to serious injury and damage to the product or other property.
- Skin/soft tissue breakdown can occur due to a number of factors, which vary by individual. Check skin frequently, at least once a day, especially those areas covering bony prominences. Redness, bruising, or darker areas (when compared to normal skin) may indicate superficial or deep tissue injury and needs to be addressed. If there is any discoloration to skin/soft tissue, STOP USE immediately. If the discoloration does not disappear within 30 minutes after disuse, immediately consult a healthcare professional.
- The Direct Mount MUST be correctly sized to the back shell.

⚠ WARNINGS:

- The product should be installed, adjusted and removed by a Permobil-authorized technician or a clinician. If the technician or clinician is unable to perform any tasks in this manual, they should seek assistance by contacting their distributor or Customer Support.
- DO NOT remove or install the back support while the individual is seated in the wheelchair.
- Check the hardware daily for tightness, wear and/or breakage. If the hardware is not securely tightened, the back support may shift up or down during operation, which may reduce or eliminate the product's benefits. DO NOT use the product if any part cannot be tightened or is broken. Periodically, the product should be thoroughly inspected by a Permobil-authorized technician.
- The wheels on the wheelchair MUST be locked before installing, adjusting or removing the back support and/or the Direct Mount from the wheelchair.

Important Safety Information, continued**⚠ WARNINGS:**

- The wheelchair **MUST** be powered off before installing, adjusting or removing the back support and/or the Direct Mount from the wheelchair.
- Dropping, Pinching, Crushing Risk: To avoid injury, hold securely and use extreme caution when installing, adjusting, removing or otherwise handling the back support with the Direct Mount. Ask for assistance if needed. **DO NOT** modify the product or use non-compatible hardware, accessories, or components, which may result in product damage or failure.
- Attaching any object, including the Direct Mount with or without the back support, will add weight that may affect the wheelchair's center of balance. Check with a Permobil-authorized technician or a clinician to determine whether the wheelchair should be modified. Proceed with extreme caution if the back support has been installed: The wheelchair's center of balance can significantly affect the stability, operation and safety of the wheelchair; if the wheelchair becomes unbalanced, injury to the user or damage to the wheelchair may result.
- If the back support is removed from the wheelchair for any reason, it must be reattached at the prescribed height. Before removing the the back support, mark the height of the back shell (using tape or a marking pen) to ensure the product is correctly reinstalled. Follow all instructions in this manual to reattach the back support to the Direct Mount and/or the wheelchair.
- **DO NOT** modify the product or use non-compatible hardware, accessories or components, which may result in product damage or failure.
- The cover should always be used with the back support. Refer to the instructions provided with the Comfort Unite Back Support.
- Metal components may become hot if exposed to high temperatures or cold if exposed to low temperatures. The product's surface will acclimate to the surrounding temperature. Take appropriate precautions, especially when the product will contact unprotected skin.
- **DO NOT** expose the product to high heat, open flames or hot ashes. Testing or certification claims, including for flammability, may no longer apply to this device when it is combined with other products or materials. Check testing and certification claims for all products used in combination with this device.
- Keep small hardware components out of the reach of children and any individuals who have a tendency to place inedible objects in their mouths. Swallowing or inhalation may lead to serious injury or death. **IMMEDIATELY** seek emergency medical assistance.

**⚠ CAUTIONS:**

- Protect the product from blunt force shocks that may cause damage and/or breakage.
- The components of the Direct Mount are not intended to be used separately from a Comfort Unite Back Support or in conjunction with other wheelchair back supports.

⚠ WARNINGS - MOTOR VEHICLE TRANSPORTATION:

- Failure to pay attention to these warnings could result in severe injury to the individual in the wheelchair or to others.
- Whenever possible, transfer out of the wheelchair and into a manufacturer-installed vehicle seat, and use the vehicle's crash-tested occupant restraint system.
- The back support must be correctly and securely installed, used, and maintained according to the instructions and warnings provided by the wheelchair and hardware manufacturers.
- Install the back support on a wheelchair so the individual is seated in an upright position and the wheelchair is facing forward during transportation in a motor vehicle. The back support has only been crash tested on a wheelchair in the forward-facing position per the test information in this operation manual.
- The back support must only be used with wheelchair frames that are equipped with securement points specified in ANSI/RESNA WC-4:2017, Section 19 (WC19) or ISO 7176-19, and that have been successfully crash tested using the methods described in Annex A of WC19 or ISO 7176-19 (respectively).
- The wheelchair must be used with an effective wheelchair securement system and a properly positioned, crash-tested pelvic and shoulder-belt restraint, or Wheelchair Tiedown and Occupant Restraint System (WTORS), following the manufacturer's instructions.
- The top of the back support should be positioned as close as possible to the top of the individual's shoulder to provide good support in the event of a crash. The distance between the top of the individual's shoulder and the top of the wheelchair back should not exceed 6½ in. (16.5 cm).
- **ALWAYS** safely store during transport to avoid damage.
- If the product has been involved in an accident during transport, see the "Maintenance" section.
- If the Direct Mount is removed from the back support and/or the wheelchair for transport, the height of the back support on the wheelchair must be adjusted during reinstallation. See the height adjustment instructions.
- For more information about safely transporting wheelchair users, refer to: <https://wc-transportation-safety.umtri.umich.edu/ride-safe-brochure/>.

Back support-hardware combinations have been dynamically tested by a third-party testing facility for use in a motor vehicle and meets all applicable criteria for wheelchair seating systems proposed in Section 5.1 of ANSI/RESNA WC-4:2017 Section 20 "Wheelchair Seating Systems for Use in Motor Vehicles" or ISO 16840-4(2009): Wheelchair Seating - Part 4: Seating Systems for Use in Motor Vehicles.

Test results apply to seating systems intended for use by adults and children with a total body mass greater than 23 kg (51 lb, or approximately 6 years of age) who are transported facing forward in their wheelchairs in all types of motor vehicles.


To determine whether motor vehicle third-party testing has been performed with a specific back support hardware and other accessory combination, refer to the operation manual provided with the back support, attachment hardware, and other accessory, or contact the equipment manufacturer.

Installation on the Wheelchair

⚠ WARNING: After installation and adjustment is complete, and before the individual transfers into the wheelchair, make sure all hardware is securely tightened. If a torque wrench is available, tighten to 9.6 Nm (85 lb-in).

To prepare for installation:

1. **Remove any existing wheelchair backrests:** Follow the instructions provided with the Permobil Corpus to remove the backrest cushion and plates.
2. **Check the Direct Mount Position:** The Direct Mount **MUST** be in the lowest position on the back support in order to attach it to the wheelchair. If the position of the Direct Mount needs to be adjusted:
 - a. Use the hex key to loosen all four socket head cap screws on the back support.
 - b. Slide down the Direct Mount on the back support as far as it will go.
 - c. Use the hex key to tighten all four socket head cap screws until snug.

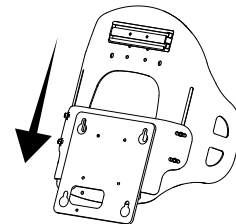To install:

3. Attach the Direct Mount to the wheelchair:

- a. Hold the back support firmly and securely and fit the four keyholes of the Direct Mount over the four pins on the wheelchair.
- b. Slide the back support straight down until it is securely attached.

IMPORTANT! All four keyholes **MUST** engage with all four pins.

4. Secure the Direct Mount and adjust the height of the back support:



- a. **Place the threaded knob in the opening and tighten.** After all pins are engaged with all four keyholes, the threaded knob will keep the Direct Mount attached to the wheelchair. The opening for the threaded knob is located beneath the lower front of the cover.

Note: If unable to access the threaded knob, install only the Direct Mount (without the back support) on the wheelchair and secure it with the threaded knob. Then install the back support on the Direct Mount.

- ⚠ WARNING:** Tightening the threaded knob does NOT necessarily mean that the Direct Mount is securely installed. The threaded knob can be tightened while the pins are NOT fully engaged in the keyholes. After tightening the threaded knob, test for a secure installation by pushing forward on the back support (toward the front of the wheelchair).
- b. **IMPORTANT! Test to make sure the Direct Mount is securely installed:** Push forward on the back support (toward the front of the wheelchair). It should NOT move forward if all four pins are engaged with all four keyholes. If it moves, repeat the installation instructions and test again.

- ⚠ WARNING:** If the Direct Mount is not securely attached to the wheelchair, the product is NOT safe to use.

- c. **IMPORTANT! Adjust the height of the back support to fit the individual:** To ensure a correct fit for the individual, adjust the height. Refer to the height adjustment instructions and safety information.

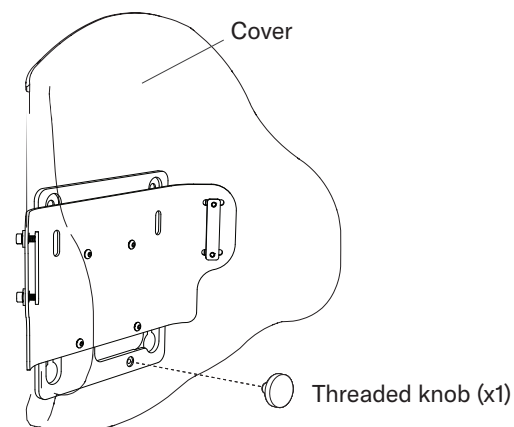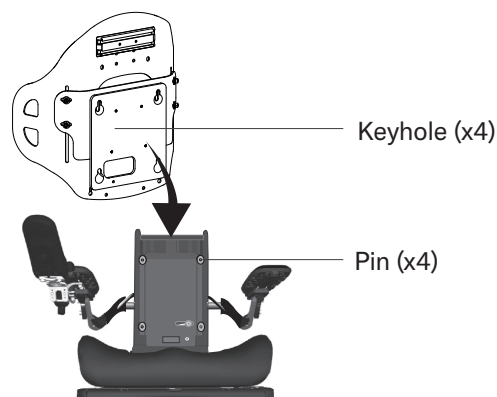The opening for the threaded knob is the same opening used for the "knob" on the wheelchair.


Height Adjustment

⚠️ WARNINGS:

- The appropriate back height for the individual should be determined by a Permobil-authorized technician or a clinician.
- If performed correctly and with caution, height adjustment may be made while the individual is in the wheelchair. If it is too difficult or unsafe to make the adjustment while the individual is in the wheelchair, have the individual transfer safely out of the wheelchair.

Note: Adjust the height only after installing the back support on the wheelchair.

1. On both sides, use the hex key to loosen the socket head cap screw just enough to allow the back support to slide up and down.
2. To select the desired height, slide the back support up or down.

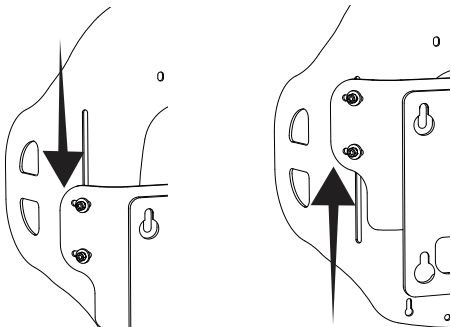
⚠️ WARNING: If a back support is NOT positioned to make contact with the Posterior Superior Iliac Spine (PSIS), pelvic support and trunk control may be reduced.

Make sure the back support height is positioned so that the bottom portion makes contact with the PSIS.

3. Use the hex key to tighten the socket head cap screws.
4. Updates to programming may be required to prevent contact between the wheelchair and the back support / accessories during power recline. Check for interference.


⚠️ WARNING: Interference during recline may result in: the back support slipping or instability; injury to the user; damage to the back support or to accessories attached to the back support; or damage to the wheelchair. Check for interference during recline when: 1) installing the back support for the first time, 2) adding any accessories to the back support, or 3) making any adjustments to the height of the back support. For reprogramming or adjustments, consult with a Permobil-authorized technician.

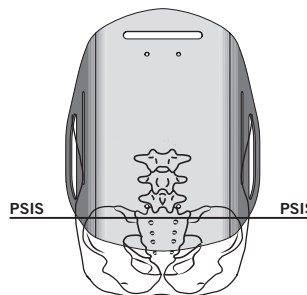Removal from the Wheelchair


To remove the back support and the Direct Mount

1. If a head support is installed, remove it.
2. Use the hex key to loosen all four socket head cap screws, then slide the back support up as far as it will go.
3. To keep the back from sliding up or down during removal, securely tighten one of the socket head cap screws.
4. Remove the back support with the Direct Mount:
 - a. Remove the threaded knob, which is located beneath the lower front of the cover.
 - b. Lift upward and then away from the wheelchair to remove the back support and the Direct Mount. The four keyholes on the Direct Mount will disengage from the four pins on the wheelchair.

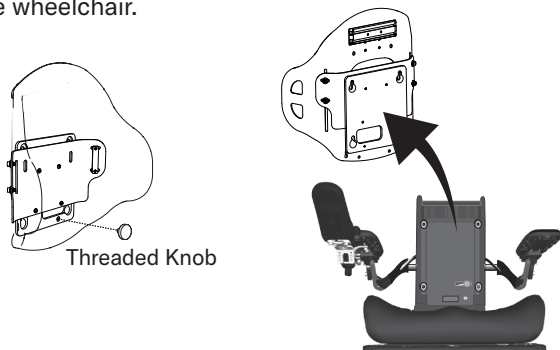
OPTION: If unable to access the threaded knob, remove the back support and Direct Mount separately:

- a. Remove the four socket head cap screws, washers and double nuts.
- b. Remove the back support from the wheelchair, but do not remove the Direct Mount Hardware from the wheelchair.
- c. Replace the four socket head cap screws, washers and double nuts.
- d. Remove the threaded knob.
- e. Lift upward and then away from the wheelchair to remove the Direct Mount from the wheelchair.

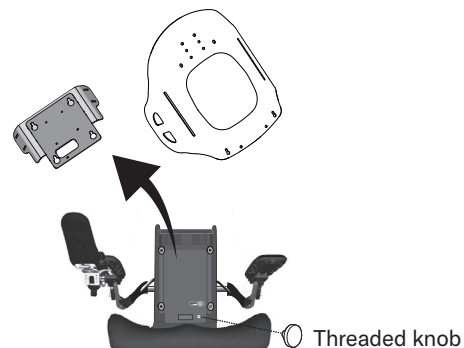Separately Purchased Back Support: Installation of Direct Mount

1. Remove the cover from the back support.
2. If necessary, remove the back support hardware. Set aside and store in a safe location.
3. If assembled, remove the Direct Mount hardware components from the Direct Mount and set aside.
 - a. On the back side of the back support, place a washer, safety washer and socket head cap screw through the opening on the side. Place the double nut on the front side of the back support. Use the hex key to tighten the socket head cap screw. Repeat this process for the second washer, safety washer, socket head cap screw and double nut.
 - b. Repeat step 1 on the other side of the back support. Do not overtighten the socket head cap screws.
 - c. Reattach the cover. Refer to the instructions provided with the Comfort Unite Back Support.
 - d. Place the Direct Mount in the lowest position on the back support. See the "Installation on the Wheelchair" section for instructions on installing on the wheelchair.

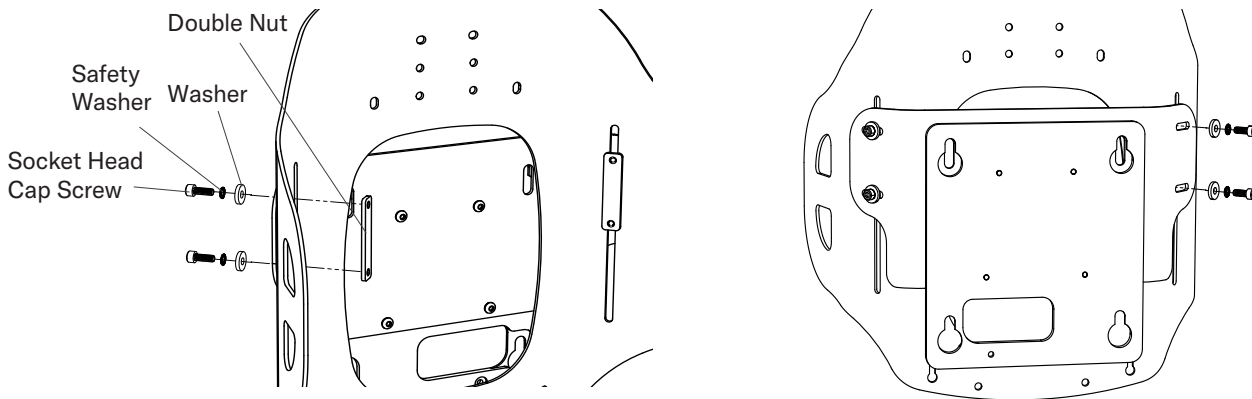Direct Mount Removal

1. Use the hex key to loosen the socket head cap screws and then remove the socket head cap screws, washers and double nuts from the back support. Set hardware components aside.
2. Remove the Direct Mount from the back support. Replace the the socket head cap screws, washers and double nuts onto the Direct Mount. See "Separately Purchased Back Support: Installation of Direct Mount" section for correct placement of hardware components.
3. Store the Direct Mount and hardware components per the instructions in the "Storage, Transport, Disposal and Recycling" section.

Troubleshooting

For additional assistance, contact a Permobil-authorized technician, distributor, or Customer Support.

Not able to access the threaded knob when installing the back support on the wheelchair.	Make sure the Direct Mount is in the lowest possible position on the back shell.
Not able to tighten the threaded knob while installing the back support on the wheelchair.	Make sure that all four keyholes on the Direct Mount are engaged with all four pins on the wheelchair. ⚠ WARNING: Tightening the threaded knob does NOT necessarily mean the Direct Mount is securely installed. The threaded knob can be tightened while the pins are NOT fully engaged in the keyholes. After tightening the threaded knob, test for a secure installation by pushing forward on the back support (toward the front of the wheelchair). See the instructions in the "Installation on the Wheelchair" section.
Not able to attach the Direct Mount to the pins on the wheelchair.	Make sure the wheelchair back support cushion and back support plate(s) are removed. Refer to the wheelchair manufacturer's instructions. Make sure the Direct Mount is not turned upside down. See the images in the installation instructions.
A component is damaged or missing.	See "Warranty" or Contact Customer Support for warranty or replacement.

Maintenance, Cleaning and Disinfecting

⚠ WARNING: DO NOT continue to use the product if any part is broken, cannot be tightened, or is not functioning properly. Periodically, the product should be thoroughly inspected by a Permobil-authorized technician.

Check daily to confirm the tightness of screws and hardware. Re-tighten loose screws, following the instructions in this manual. Check components daily for wear or damage and replace if necessary.

⚠ WARNINGS:

- Cleaning and disinfecting are separate processes. Cleaning must precede disinfection. Before use by a different individual: clean, disinfect, and check product for proper functioning.
- Clean the product regularly and when it becomes soiled.
- Follow safety guidelines on the chemical containers.

⚠ CAUTIONS:

- DO NOT use abrasives such as steel wool or scouring pads on the product.
- DO NOT use cleaning products containing petroleum or organic solvents [i.e. acetone, toluene, Methyl Ethyl Ketone (MEK), naphtha, dry cleaning fluids, adhesive removers] on the product.

To Clean: Hand wash with a sponge, using soap and water. Rinse with clean water. Wipe dry with a clean cloth and then allow to air dry completely.

To Disinfect: Wipe gently with a cloth dampened with a household disinfectant. Rinse with clean water. Wipe dry with a clean cloth and then allow to air dry completely.

Recommendation: When washing by hand, use room-temperature water.

Care Symbols

Wipe with damp cloth.

Dry flat.

Hand wash.

Storage, Transport, Disposal and Recycling

Storage or Transport: Clean and disinfect the product. Place in a container that protects the product from moisture, contaminants and damage. After storing or transporting, inspect the product to ensure it was not damaged (see "Troubleshooting"), and follow instructions to prepare the product for use.

Disposal: The components of the products in this manual are associated with no known environmental hazards when properly used and when disposed of in accordance with all local/regional regulations. At the end of its useful life, treat as healthcare waste and dispose of according to local healthcare-waste guidelines. Incineration must be performed by a qualified, licensed waste management facility.

Recycling: Contact a local recycling agency to determine recycling options for the product.

Warranty

This product is covered by a 5 year warranty against manufacturer defects.

If a manufacturing defect should occur, discontinue use immediately; note that the product will be replaced at no cost to the buyer. Wear from regular use, misuse, damage by an accident or damage from "an act of Nature" (e.g., flood, tornado, earthquake and/or fire) is not considered a manufacturing defect. Replacement parts are available by contacting your supplier, distributor, or retailer.

The following will void the warranty: modifications made in attempt to interface with other products (where compatibility is limited, as described in product information); alterations made to the product; or defects caused by irregular use.

Warranty Replacement Procedure

Warranty claims should be initiated by the original purchaser, who should contact the Durable Medical Equipment (DME) supplier, distributor or retailer from whom the product was purchased, if they are available, or any authorized Permobil DME supplier, distributor or retailer.

Whenever possible, the DME supplier, distributor or retailer with a warranty issue should contact Customer Support and provide the original purchase order number, sales order number, or invoice number.

Before a warranty replacement order or quote may be processed, Customer Support must issue a Return Merchandise Authorization (RMA). Repeated unsubstantiated warranty claims may void the warranty of the product.

Información de contacto

 The Comfort Companies, LLC, 1960 South Calhoun Road
New Berlin, Wisconsin 53151 EE. UU.

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