


Permobil tough

Finally, armrest and legrests designed to withstand the rigors of daily repositioning and transfers. Unlike traditional, wobbly rehab arms, the M1 locking flip back armrests are durable and secure. The center mount and multi-angle adjustable swingaway legrests are built to last with a full range of options to address your lower extremity positioning needs.



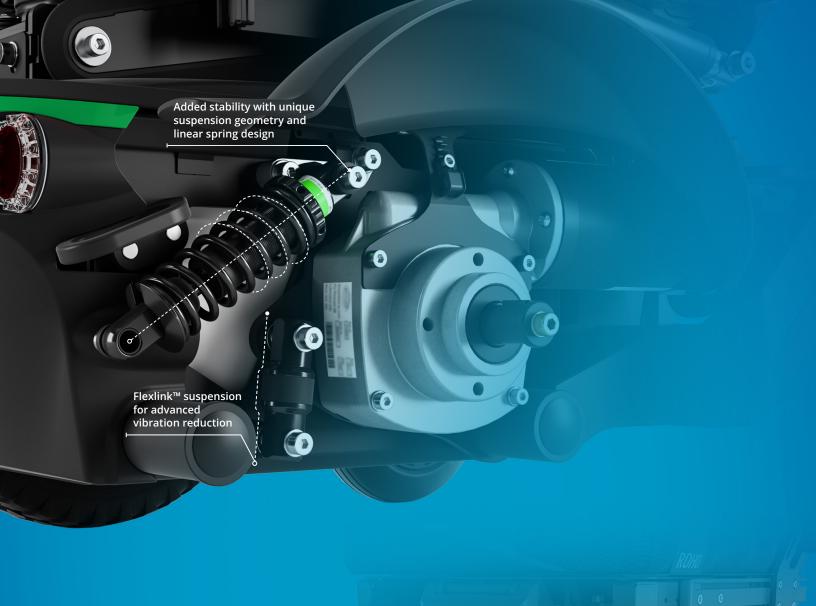

permobil

Proven technology

Built on the modular chassis of the Permobil M3

Corpus®, the M1 features the same reliable drive units,
patented FlexLink™ suspension and precision engineered

CasterGlide™ technology to help ensure smooth, confident
operation in a variety of driving scenarios.

7" performance caster design + adjustment free flutter dampening



Specifications

Maximum user weight	300 lbs ^{††}
Maximum speed	5 mph
Range	16 miles*
Base width	24"
Base length	35.5"
Minimum turning radius	20.5"
Ground clearance	3"
Recommended battery type	Group 34
Weight including batteries	319 lbs**
Drive electronics	R-Net 90A, VR2 PM 90 (Advanced) or PM 70 (Advanced and Basic)
Seat-to-floor height	16", 17", 18"
Tilt options	0° – 50°

Recline options	85° – 110° (manual)
Legrest angle	Center Mount 80° – 140° (power) [>19" seat depth] 90° – 150° (power) [≤19" seat depth] 80° – 135° (manual)
	Swing-away 65° – 80°
Backrest height	18" – 29" (by 1" incr.)
Seat depth	14" – 22" (by 1" incr.)
Seat width	14" – 22" (by 2" incr.)
Armpad options	2"x10", 2"x14", 4"x10", 4"x13"

[†]Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions.

Color options



Radiant Red



Sunburst Orange



Power Pink



Galactic Green



Ocean Blue



Midnight Black

Hydrographic options



Carbon

Fiber



CDG



Mossy Oak®

Breakup®

. Country™





American Liberty

Midnight Patriot



[#]The maximum user weight for a seat depth of 14" cannot exceed 165 lbs.

^{*} Actual driving range will vary based on driving and battery conditions.

^{**}Includes weight of base, curved backrest and batteries.