

ActiveHeight & ActiveReach packages



Up to 45° of ActiveReach and 14" of ActiveHeight

Power wheelchairs

Corpus® Seating System




	M3	M5	MVS	F3	F5	F5VS
↑ ActiveHeight	12"	14"	12"	12"	14"	14"
↶ ActiveReach	5°, 10° & 20°	5°, 10° & 20°	Standing	5°, 10°, 20° & 30°	5°, 10°, 20° & 45°	Standing

*Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions.



Created for you

Corpus is the most ergonomic and capable wheelchair seating system in the world. It is elegantly designed and configured to empower functional movements while supporting your optimal body position.

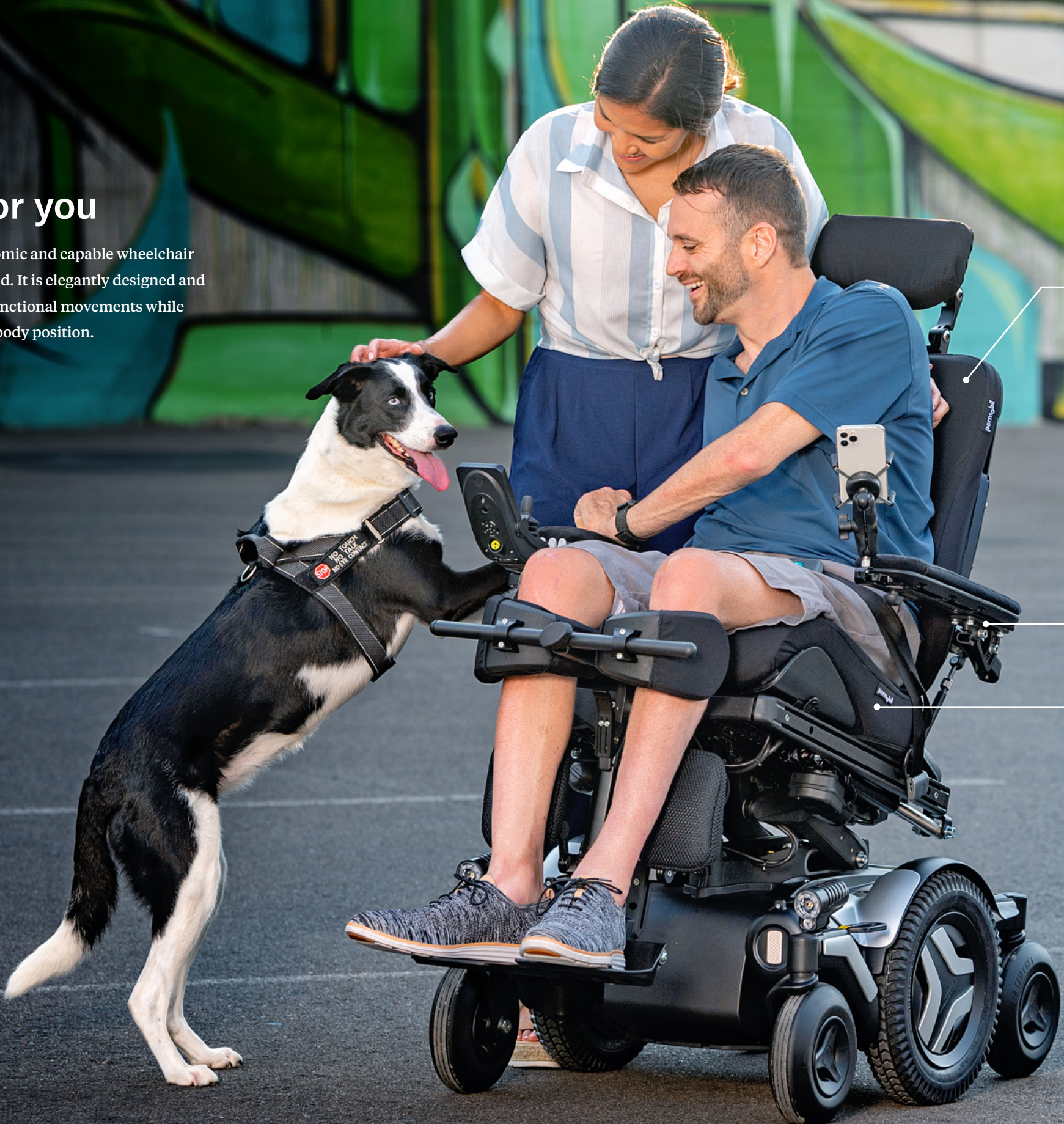

Dual density back support foam for greater comfort and surface area contact



Leatherette



STRETCH-AIR®

Multi-axis armrests allow independent height, angle and optional medial swivel adjustment

Redesigned pelvic well insert for improved pressure redistribution



Meeting your positioning needs

Manage your seating goals with multiple memory seat positions and sequences such as Independent Repositioning Mode (IRM). IRM enables automatic sequencing of tilt, leg elevation, and recline promoting maximum benefit of pressure relief protocols without loss of positioning. Additional positions and sequences allow complete customization of everyday positioning needs.

Access memory positions and sequences at the touch of a button

Up to 180° power legrest elevation

Seamless shear reduction and optimized pivot points keep you in position

Up to 180° power recline and 50° posterior tilt

Programming the Corpus seating system with your wheelchair team is simplified through the QuickConfig app



3.5" rearward shift during elevate optimizes your center of gravity and reduces chair footprint

12" & 14" of ActiveHeight
Drive up to 3.5 mph with full suspension at any elevated height

3 points of contact for added stability



ActiveHeight

Dynamic seat elevation

ActiveHeight raises your seat up to 14". It is unique in that while elevating, the seating system repositions 3.5" rearward to maintain forward stability and create a more compact footprint enabling you to get closer to the important things in your life. ActiveHeight empowers you to have face-to-face conversations, order food at a high counter, and accomplish many other daily tasks with ease.

ActiveHeight + ActiveReach

Seat elevation and forward tilt combine seamlessly to give you the ultimate advantage in accessing your environment.

Swing-away knee support option, for use on Corpus VS legrests, are 75% lighter than single-post knee supports



Flexible on-chair storage makes it easy to add or remove knee supports as needed



ActiveReach Functional forward tilt

Face daily tasks head on with ActiveReach. By combining Corpus seat functions with up to 45 degrees of forward tilt, you can achieve improved ergonomics for overhead reaching and table top activities. ActiveReach can also help improve posture for social engagement and body mechanics during transfers.

10° of no charge ActiveReach with ActiveHeight





Adjustable chest support

Swing-away or single-post knee support options

Articulating power footplate to support standing


Corpus VS

Superior power standing

Fully programmable standing sequences allow user position to be optimized using supine, semi-reclined, or sit-to-stand techniques. Further customize standing positions with adjustable chest support, multiple knee support options, and power articulating footplates. From face-to-face interactions to improved health outcomes, the Corpus VS seating system empowers you to stand with confidence.



Stand while driving

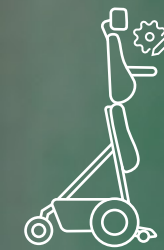

Customizable standing sequences



Full range power functions