

BodiLink<sup>®</sup>  
Accessories

permobil






# BodiLink Lateral Pelvic/Thigh Supports

Compatible on virtually any power or manual wheelchair

## Hardware

- Telescoping Taper (TT) and Gear Tooth (GT) hardware options
- Fixed or removable with high-visibility release lever
- Extra links (GT) and long extension arm (TT) allow for greater medial reach where needed

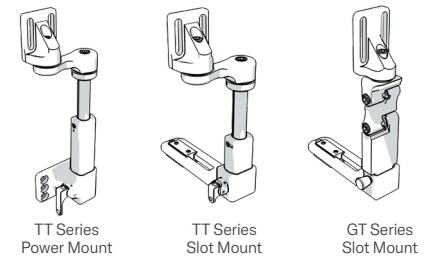

### TT Series

- Mounts to Permobil, Quantum Rehab®, Quickie®, Rovi® and Avid power seating systems and manual chairs with solid seats
- Precision telescoping height and angle adjustment
- Small, medium, and large hardware options accommodate for a range of cushion thicknesses
- Adjustable release lever position and direction with one-handed operation is great for patient ADLs




Medium Size TT Slot Mount Hardware Pictured

## Pads

### Premium

- New molded foam pads in 4 sizes
- 3 cover choices: Comfort-Tek® with foam, Stretch-Air® with foam and GlideWear® with foam

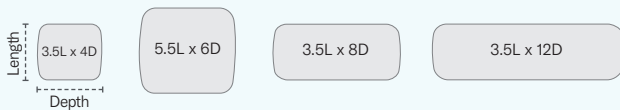


Premium pad covers enclose attachment points and allow easy removal for washing.

### Basic

- 10 pad sizes offered in zero elevation
- 3 cover choices: Comfort-Tek, Stretch-Air and GlideWear

#### Zero elevation

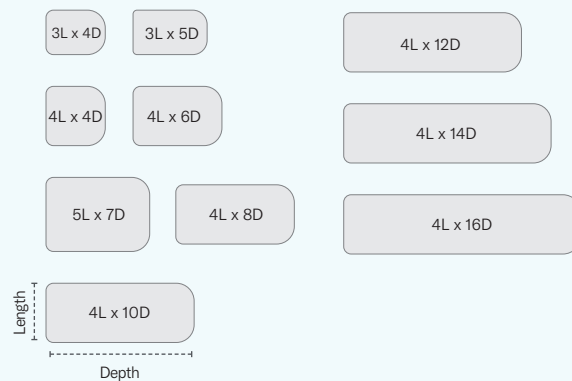


# BodiLink Head Supports

Advanced proprietary design is durable and reliable


## Pads

### Premium

- New piano hinges on 10W and wider pads offer more reliable and precise lateral adjustment
- 3 cover choices: Comfort-Tek, Stretch-Air and GlideWear®
- 6 sizes: soft and extra soft firmness levels
- Pad shape and adjustment capabilities offer precise alignment, greater surface contact, stable positioning, and clear line of vision for the user

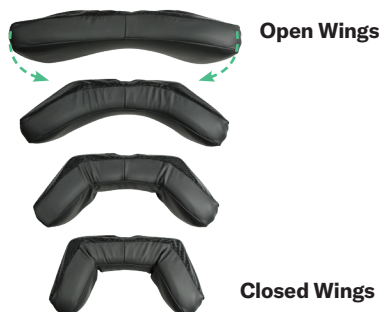

### Premium

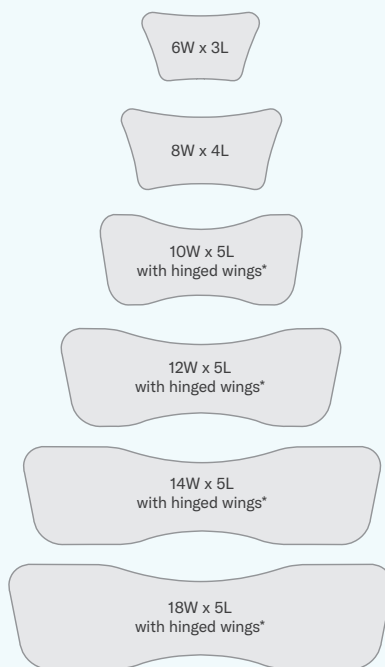

\*Hinged wings allow for lateral contouring and support on 10W, 12W, 14W, & 18W pads.

### Basic



Basic pads are moldable for contoured support.


## Hardware

- 7 Styles: 5 Piano Taper (PT) and 2 Gear Tooth (GT) options
- Dedicated mid-line styles when lateral positioning is unnecessary
- Mounting options for Acta-Series, Corpus®, AGILITY® back supports and more
- Easy installation and high-visibility release lever


### PT Series

- Durable PT hardware comes in 5 styles ranging in functionality and reach
- Precision link adjustment for ideal positioning
- Tool-less wing lock option for on-the-spot adjustments
- Optional flip-down release allows the head support pad to move away from the user's head when transferring, transporting, or performing other ADLs

### Mounting Blocks

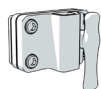

Posterior -  
Acta-Series &  
Comfort VISCO

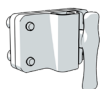

Anterior -  
ROHO® AGILITY® &  
Tru-Shape™



Permobil A-  
Corpus® backs 23"-28" H  
Uni-Track compatible



Permobil B-  
Corpus® backs 20" H



Quantum-  
TRU-Comfort 2 back  
supports




### PT Hardware Options

Choice of 7" or 3" stem length

Tool-Less Wing Lock

2", 3", and 4" Extra Links

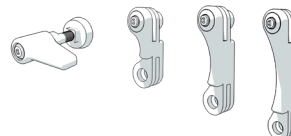

# BodiLink Lateral Trunk Supports

Designed with the body's shape in mind with rounded edges to ensure comfort

## Hardware

Durable hardware ensures pads stay firmly in place longer

### GT2 Series

- Spring-loaded joints for easy set-up and positioning
- Unique wrist joint allows pad rotation to match user shape
- Modular design for customizable link configuration
- Mounting options are durable and a positive fit for Acta-Series, Corpus, AGILITY, and TRU-Comfort 2 back supports
- Winter/Summer Lateral Trunk Support mounting hardware available for ACTA-EMBRACE® back support





GT2 Swing-Away Hardware Pictured

## Pads

- 6 pad sizes offered in each premium & basic, moldable & rigid, supportive & active, giving 48 options to fit any user

### Premium

- 7 cover choices: bare, Comfort-Tek, Stretch-Air and GlideWear® for skin protection, and each fabric available with added foam for a non-zeroing, soft feel
- Anatomically designed pad offers more natural positioning and greater surface contact
- Molded foam is durable, reducing risk of foam breakdown

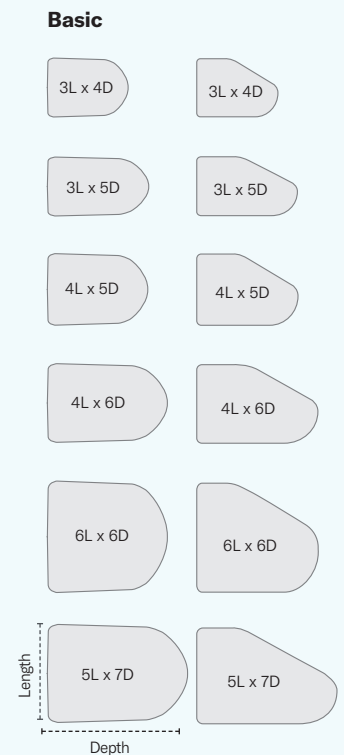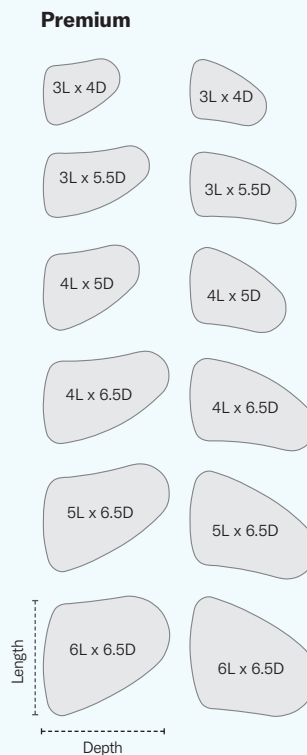

# BodiLink Medial Knee/Thigh Supports

Let the accessory do it! Medial Knee/Thigh Supports decrease the need for dramatic positioning elements within a cushion

## Hardware


- Mounts to manual wheelchairs with Comfort Removable Solid Seat Pan
- 3 height options available to fit cushions from 1" - 7" thick
- Flip-down mechanism makes transferring easier by swinging the pad and hardware underneath the wheelchair
- Height adjustment is made through a vertical slider
- Depth adjustable through a slotted track underneath the pad




## Pads

- 2 styles: wedge and oval
- 3 sizes in each style
- 3 cover choices: Comfort-Tek, Stretch-Air and GlideWear®
- Pad is depth adjustable

### Wedge

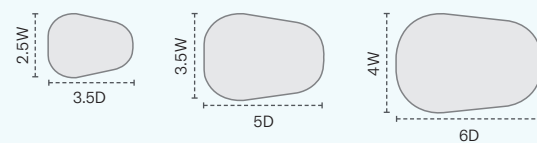


### Oval

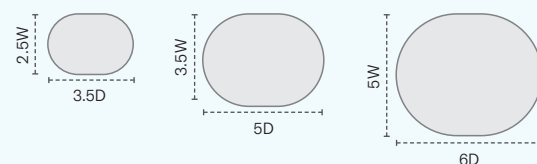

# GlideWear® Covers

Premium shear reduction included

Get the positioning and support you need with tissue protection in mind. Experience the benefits of up to 88% reduction in shear forces\* by choosing from our variety of exclusive GlideWear cover options.

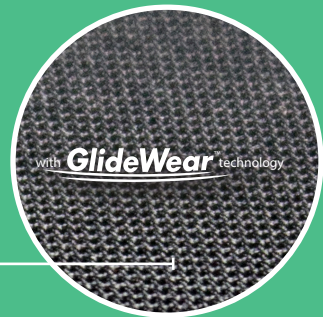

Patented 2-ply, low-friction fabric transfers shear forces from the user to the fabric.

## GlideWear cover options are currently available on:

BodiLink Head Support	Premium Pads Only
BodiLink Lateral Trunk Supports	Premium Pads Only
BodiLink Lateral Pelvic/Thigh Supports	All Support Pads
BodiLink Medial Knee/Thigh Supports	All Support Pads

*\*Results of modified ISO 16840-2 testing of a foam cushion comparing a typical cover to a GlideWear cover.*

*Decisions about how to code a product are the responsibility of the provider/supplier. Since Medicare coding is subject to change, the provider should always confirm the HCPCS code and coverage criteria as part of the client assessment process*

*Lifetime guarantee on all BodiLink Accessories. Should you experience performance issues related to the PT Head Support joints, we will refund service costs. Some restrictions apply.*

*GlideWear is a registered trademark of MIPS AB CORPORATION SWEDEN.*

*Quantum Rehab is a registered trademark of Pride Mobility Products Corporation.*

*Quickie is a registered trademark of Sunrise Medical US LLC.*

*ROVI is a registered trademark of Shoprider Mobility Products, Inc.*