

Power wheelchairs


# Corpus® Seating System



permobil

# Created for you

Corpus is the most ergonomic and capable wheelchair seating system in the world. It is elegantly designed and configured to empower functional movements while supporting your optimal body position.

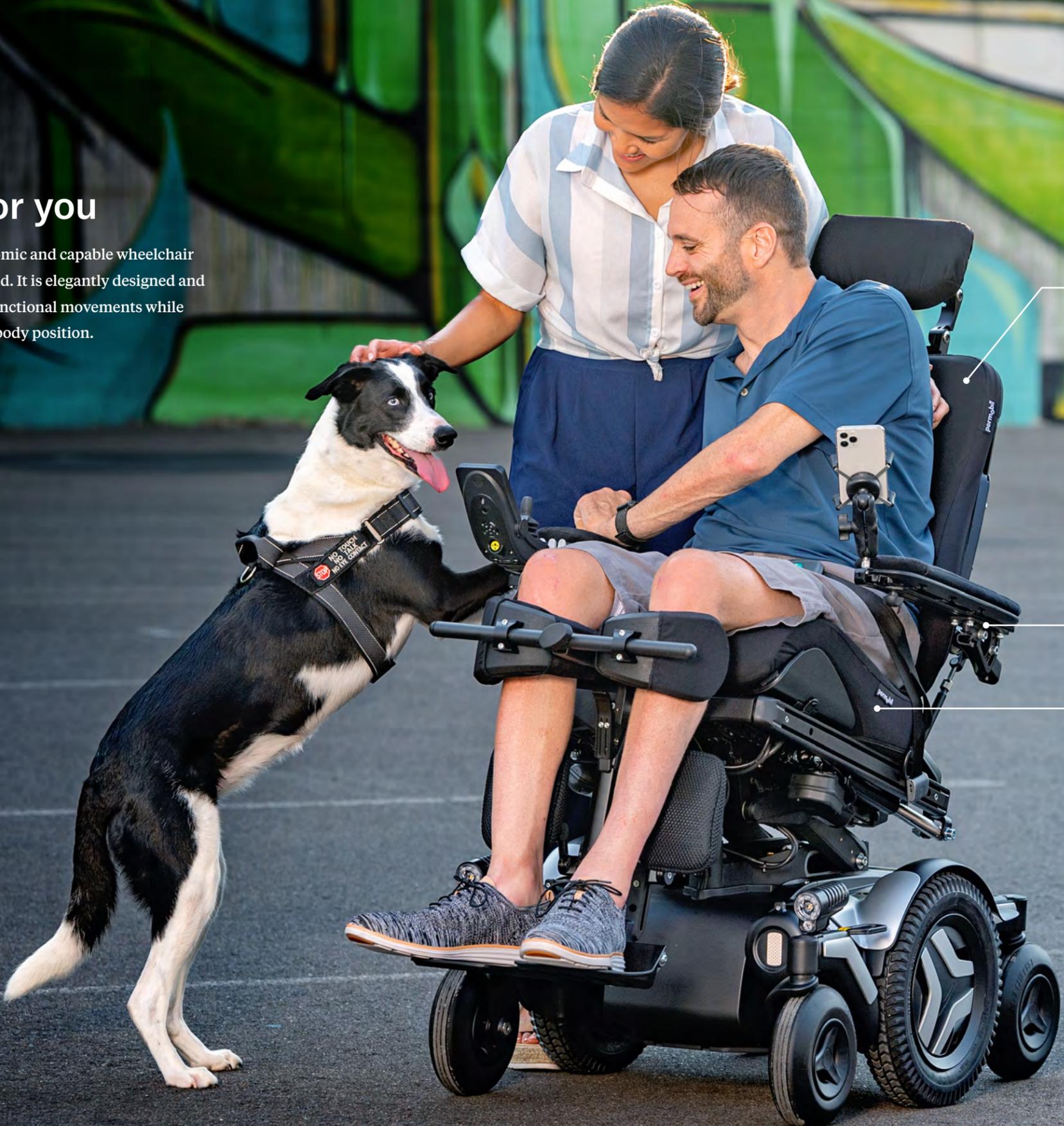

Dual density back support foam for greater comfort and surface area contact



Leatherette



STRETCH-AIR®

Multi-axis armrests allow independent height, angle and optional medial swivel adjustment

Redesigned pelvic well insert for improved pressure redistribution



Access Independent  
Repositioning Mode at  
the touch of a button

Up to 180° power  
legrest elevation

Seamless shear reduction  
and optimized pivot points  
keeps you in position

Up to 180° power recline  
and 50° posterior tilt

## Independent repositioning

Manage your seating goals with multiple memory seat positions and Independent Repositioning Mode (IRM). IRM enables automatic sequencing of tilt, leg elevation and recline promoting maximum benefit of pressure relief protocols without loss of positioning.



3.5" rearward shift during elevate optimizes your center of gravity and reduces footprint

12" & 14" of ActiveHeight  
Drive up to 3.5 mph with full suspension at any elevated height

3 points of contact for added stability



## ActiveHeight

Dynamic seat elevation

ActiveHeight raises your seat up to 14". It is unique in that while elevating, the seating system repositions 3.5" rearward to maintain forward stability and create a more compact footprint enabling you to get closer to the important things in your life. ActiveHeight empowers you to have face-to-face conversations, order food at a high counter, and accomplish many other daily tasks with ease.

## ActiveHeight + ActiveReach

Seat elevation and forward tilt combine seamlessly to give you the ultimate advantage in accessing your environment.



## ActiveReach

Functional forward tilt

Face daily tasks head on with ActiveReach. By combining Corpus seat functions with up to 45 degrees of forward tilt, you can achieve improved ergonomics for overhead reaching and table top activities. ActiveReach can also help improve posture for social engagement and body mechanics during transfers.

10° of no charge ActiveReach with ActiveHeight





Adjustable chest support

Moldable knee block

Articulating power footplate

Intelligent support wheels deploy automatically for stability in standing

## F5 Corpus VS

### Superior power standing

Choose how you stand from three customizable standing sequences. The F5 Corpus VS offers sit-to-stand, semi-reclined or lay-to-stand techniques that are fully programmable to your individual needs. From face-to-face interactions to improved health outcomes, the F5 Corpus VS empowers you to experience the vast array of standing benefits.



Stand while driving



Customizable standing sequences



Full range power functions

# ActiveHeight & ActiveReach packages



	M3	M5	F3	F5	F5VS
↑ <b>ActiveHeight</b>	12"	14"	12"	14"	14"
↶ <b>ActiveReach</b>	5°, 10° & 20°	5°, 10° & 20°	5°, 10°, 20° & 30°	5°, 10°, 20° & 45°	Standing

\*Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions.