


Power wheelchairs

# M-Series




permobil



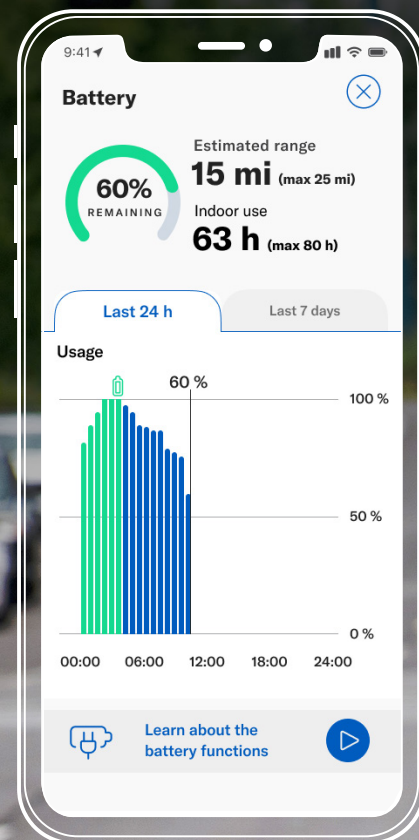

Battery status and travel range estimates help you know how far you can go

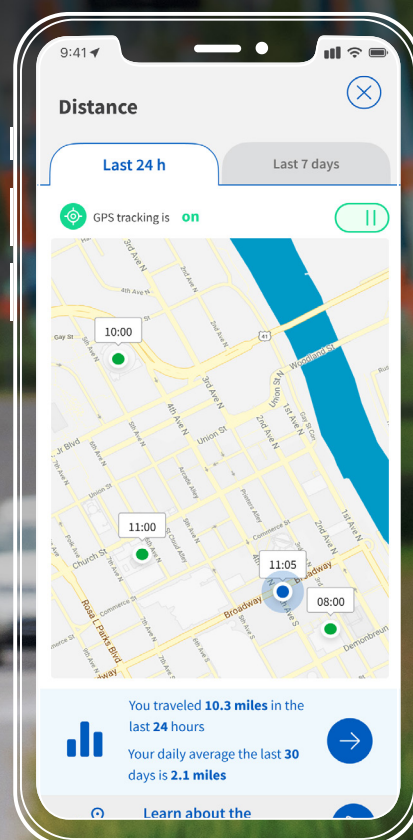

Integrated map illustrates daily and historic travel behavior



MyPermobil

MyPermobil connectivity is an available option on all new M-Series models.



## More adventure

The Permobil battery and VoltPro charger work together to provide you with up to 25% greater driving distance, reduced charging times, and longer battery life. Use the MyPermobil app to get even more from every adventure with seat activity tracking, battery range estimates, and an integrated map.



Driving range  
up to 20 miles



Reduced  
charging times



Longer  
battery life



# More freedom

Corpus®, our most capable power wheelchair seating system, now provides you with even more positioning options. Enhanced stability on the M-Series enables you to optimize your repositioning, functional reach, and power drive positions for even greater freedom of movement.



Full suspension at any elevation, seat position, or drive speed

# Reach for more

ActiveHeight™ and ActiveReach™ combine seamlessly to give you the ultimate advantage in accessing your environment. Seat elevation with integrated forward tilt creates the most innovative, functional, and stable seating solution on the market today.

	M3	M5	MVS
↑ <b>ActiveHeight</b>	12"	14"	12"
↶ <b>ActiveReach</b>	20°	20°	Standing

10° of no charge ActiveReach with ActiveHeight



Up to 3.2 mph elevated drive speed

3.5" rearward shift during elevate for added stability and reduced footprint

12"

3 points of contact for added stability



# Conquer more

FlexLink™ and FlexLink Pro™ suspension balances stability, traction, and ride comfort for optimal performance and advanced vibration reduction. DualLink suspension on the new M Corpus VS optimizes stability in standing without sacrificing active suspension in any position. The result is an unbeatable stand and drive experience. Factory optimized and adjustable suspension adapts to your needs and environment. Be ready to handle any adventure with high-powered LED lights to guide your way.



FlexLink  
M3 Corpus



FlexLink Pro  
M5 Corpus

2X more powerful  
LED lights for true  
pathway lighting and  
enhanced visibility





# M3 Corpus

Meticulously engineered to help you maintain better positioning, feel more secure, and reduce driving fatigue so you can conquer daily activities with ease. Its compact base offers maneuverability in tight spaces without sacrificing performance. Featuring the innovative FlexLink suspension geometry with independent caster swing-arms helping to keep you stable and in control as you traverse the obstacles in your world.




# M5 Corpus

Take control with high-performance motors enabling a top speed of 7.5 mph and exceptional low-end torque that lets you pivot like a pro. Experience unsurpassed stability and ride comfort with FlexLink Pro independent drive wheel suspension. The M5 Corpus will take you to new heights with increased stability and full suspension even in 14" of seat elevation.



# M Corpus VS

The new M Corpus VS brings best-in-class power standing to the maneuverability of a mid-wheel drive base. Fully programmable standing sequences allow for user position to be optimized using supine, semi-reclined, or sit-to-stand techniques. DualLink suspension balances stability, traction, and ride comfort for a smooth power stand and drive experience at home or on the go.

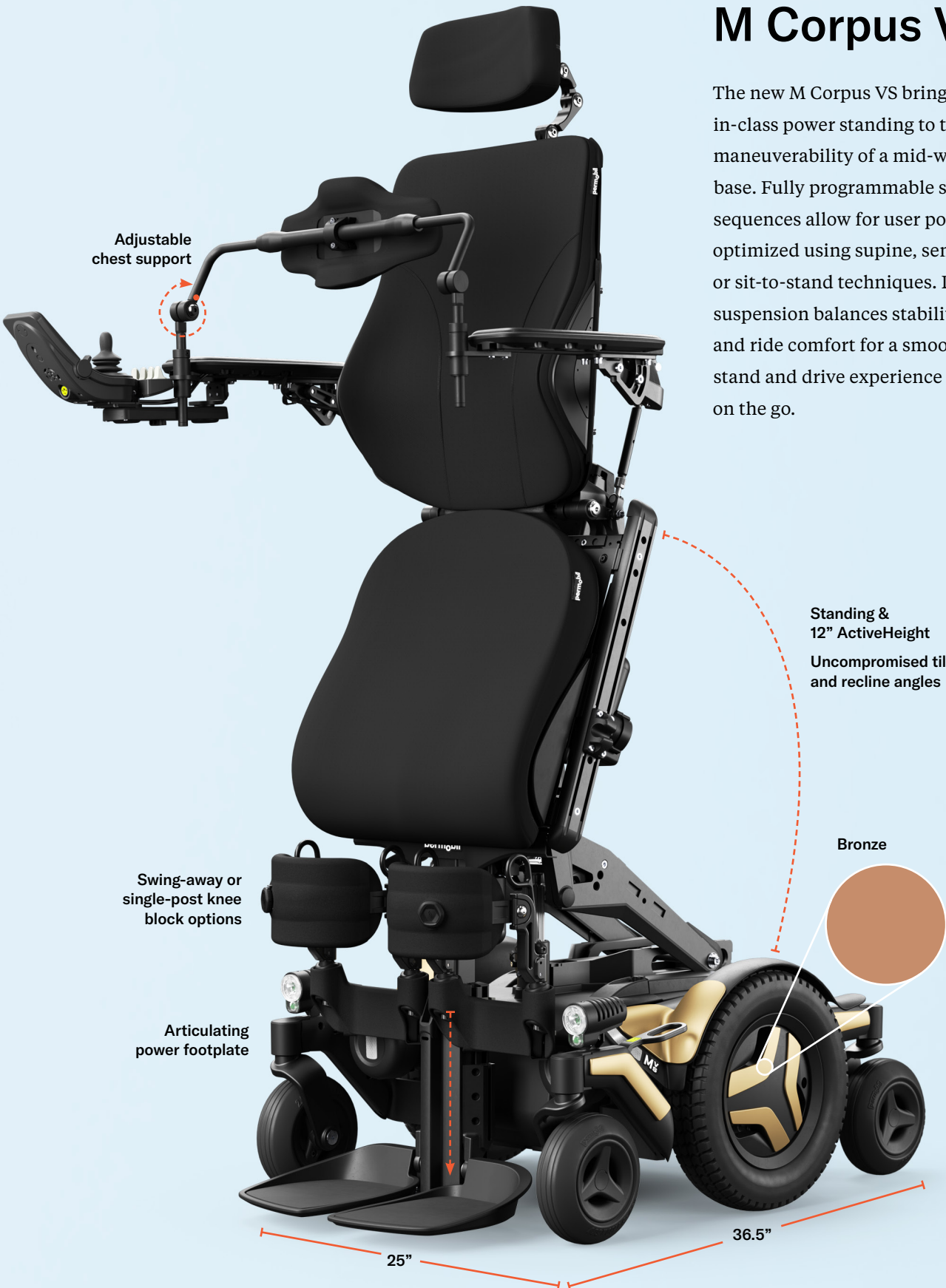



# Swing-away knee supports

The new swing-away knee support, for use on Corpus VS models, are 75% lighter than previous single-post knee supports and offer support through the standing sequence. Thanks to flexible on-chair storage, they are easy to add or remove as needed.



Release mechanism designed for a wide variety of hand function

Seat-rail mounted parking spot for convenient on-chair storage

Integrated BOA® Fit System for a precision fit

BOA is a trademark of BOA.

# BodiLink on Corpus

Specially designed for the Corpus seating system, BodiLink accessories feature best-in-class hardware. Offering true anatomical support with the highest level of adjustability and comfort to meet your individual positioning needs.



**Head Support—**  
Easy to adjust piano hinge and double-taper joints for maximum strength

**Lateral Trunk Support—**  
Designed for enhanced comfort and trunk alignment with unique wrist joint rotation adjustability allowing pad to match body contours

**Lateral Pelvic/Thigh Support—**  
Easily position along the UniTrack rail with 360° pad adjustment and user-friendly quick release lever


# Accent options




Midnight Patriot



American Liberty

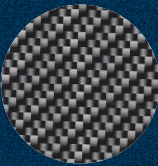

Carbon Fiber

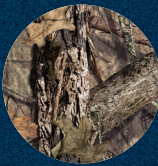

Mossy Oak® Breakup® Country™



Water Lily



Liquid Fusion\*



Midnight

Arctic

Violet

Sky

Aqua

Spring

Canary

Sunset

Chili

Flamingo

Platinum

Bronze

Gold

Please note: Actual product colors may vary from colors shown in this brochure

\*Liquid Fusion is currently only available on the M Corpus VS





Live seating angles and activity tracking in MyPermobil

## Specifications

	M3 Corpus	M5 Corpus	M Corpus VS
<b>Maximum user weight</b>	330 lb	330 lb	265 lb
<b>Maximum speed</b>	6 mph	7.5 mph	6 mph
<b>Range</b>	16 mi (Grp 34)   20 mi (Grp 24)*	20 mi*	20 mi*
<b>Base width</b>	24"	25"	25"
<b>Base length</b>	35.5"	37.5"	36.5"
<b>Base turning radius</b>	20.5"	21.5"	22"
<b>Ground clearance</b>	3"	3"	3"
<b>Recommended battery type</b>	Group 34 or 24	Group 24	Group 24
<b>Weight including batteries</b>	385 lb (Grp 34)   405 lb (Grp 24)	423 lb	424 lb
<b>Drive electronics</b>	R-net 120A	R-net 120A	R-net Power Platform 120A
<b>Seat-to-floor height</b>	17.5", 18.5", 19.5"	17.5", 18.5", 19.5"	17.5", 18.5", 19.5"
<b>Seat elevation</b>	Optional 12"	Optional 14"	12"
<b>Drive wheel</b>	14"	14"	14"
<b>Casters</b>	7"	8"	7"
<b>Tilt options</b>	Posterior: 0° – 50° ActiveReach: 5°, 10° & 20° (up to 300 lb)	Posterior: 0° – 50° ActiveReach: 5°, 10° & 20° (up to 300 lb)	Posterior: 0° – 50° Standing
<b>Legrest angle - power</b>	85° – 180°	85° – 180°	85° – 180°
<b>Backrest angle - power recline</b>	85° – 180°	85° – 180°	85° – 180°
<b>Backrest height</b>	20", 23" – 28" (1" incr.)	20", 23" – 28" (1" incr.)	23" – 28" (1" incr.)
<b>Seat depth</b>	14" – 22" (1" incr.)	14" – 22" (1" incr.)	14" – 22" (1" incr.)
<b>Seat width</b>	17" – 23" (2" incr.)	17" – 23" (2" incr.)	17" – 23" (2" incr.)

\* Actual driving range will vary based on driving and battery conditions.

Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions. Chair weights are approximate and include base, batteries, power seat functions, and seating system without accessories.