


ROHO® Hybrid Select

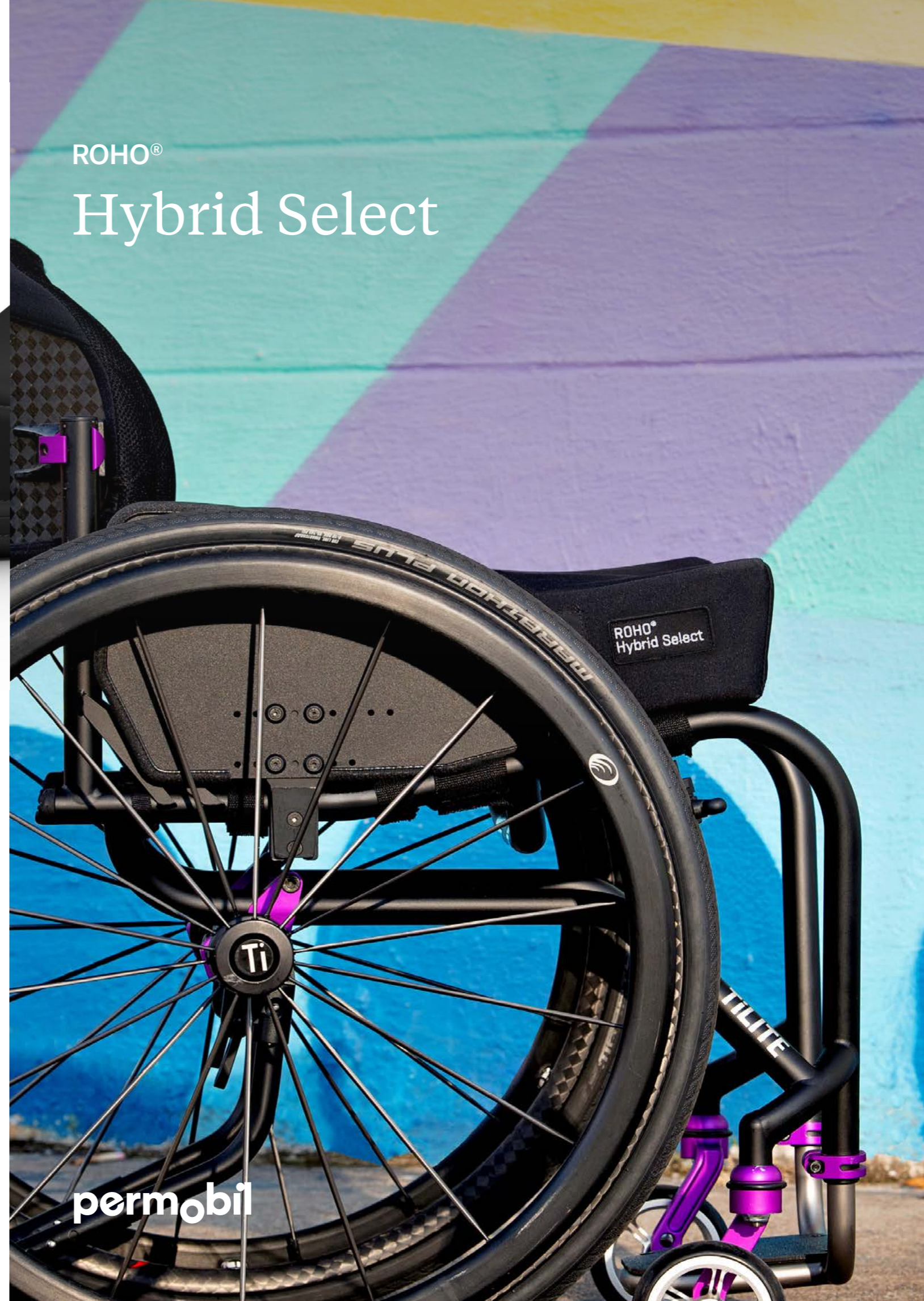

Specifications

Cover type	Stretch-Air®
Product weight	1.45 kg (3.2 lbs)*
Weight limit	227 kg (500 lbs)
Cushion warranty	36 month limited
Cover warranty	6 month limited

*Product weight is based on a size 16" x 16" cushion. Product weight varies depending on cushion size.

Size options

14" x 14" (35.5 cm x 35.5 cm)	16" x 18" (40.5 cm x 45.5 cm)	18" x 18" (45.5 cm x 45.5 cm)
14" x 16" (35.5 cm x 40.5 cm)	16" x 20" (40.5 cm x 51.0 cm)	18" x 20" (45.5 cm x 51.0 cm)
15" x 15" (38.0 cm x 38.0 cm)	17" x 16" (43.0 cm x 40.5 cm)	20" x 16" (51.0 cm x 40.5 cm)
15" x 16" (38.0 cm x 40.5 cm)	17" x 17" (43.0 cm x 43.0 cm)	20" x 17" (51.0 cm x 43.0 cm)
15" x 17" (38.0 cm x 43.0 cm)	17" x 18" (43.0 cm x 45.5 cm)	20" x 18" (51.0 cm x 45.5 cm)
16" x 15" (40.5 cm x 38.0 cm)	17" x 20" (43.0 cm x 51.0 cm)	20" x 20" (51.0 cm x 51.0 cm)
16" x 16" (40.5 cm x 40.5 cm)	18" x 16" (45.5 cm x 40.5 cm)	
16" x 17" (40.5 cm x 43.0 cm)	18" x 17" (45.5 cm x 43.0 cm)	

The following are trademarks and registered trademarks of Permobil: ROHO®, ISOFLO®, Stretch-Air® and DRY FLOATATION®.

Advanced hybrid seating with superior positioning

The Hybrid Select pelvic loading area adjusts to your unique shape, weight, and positioning changes. The ROHO overlay and IT insert utilises the ISOFLO® Memory Control that allows you to easily find and set the right position each and every time.



3-Zone ISOFLO Memory Control for adjustable positioning

New unique ROHO cell design for skin protection and immersion



Coccyx Relief Cutout
Provides relief to the
coccyx/tailbone

Contoured Foam
Leg Troughs

Pelvic well using patented
ROHO DRY FLOATATION®
technology for hydrostatic
distribution and 1.25" (3.2 cm)
immersion.

Unique lateral stability

The innovative design of the Hybrid Select uses the positioning capability of a contoured foam base with high resiliency (HR) foam and leg troughs to offer you support, lateral stability, and provide relief to the coccyx/tailbone.

Proven skin protection

With Hybrid Select, experience truly functional stability combined with great skin and soft tissue protection. The contoured design of the Hybrid Select allows your pelvis to be immersed in air which shapes to you, moves with you, and protects your skin and soft tissue.

Optional offloading comes standard

The Hybrid Select provides flexibility and forgiveness that other offloading cushions do not, making it easier to be in the right position each and every time. The Hybrid Select offers an exclusive optional offloading technique by removing the ischial tuberosities (IT) air insert.



Offloading
Removable IT air insert
for optional offloading



“The Hybrid Select gives me the positioning I need and makes it easier for me to transfer.

— Eric, Permobil User