

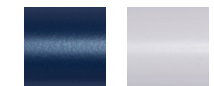
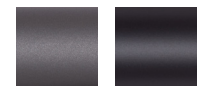

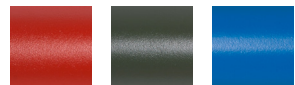




## Tattoo options

- 1 Mosaic
- 2 Flower garden
- 3 Need for speed



## Frame colors

<b>High gloss paint</b>  TiLite red Neon pink High gloss black	<b>Metallic paint</b>  Tangerine metallic Acid green metallic Ocean blue metallic	<b>Pearlescent paint</b>  Midnight blue pearl White pearl	<b>Matte paint</b>  Titanium grey Black lava
<b>Candy paint</b>  Ruby red Azure blue Electric plum	<b>Matte textured paint</b>  Flame orange OD green Hyper blue	<b>Natural Titanium</b>  Satin Polished (CMD 87)	
<b>Component colors</b> <b>Anodized</b>  Red Orange Blue Purple Green Gold Black			

## Specifications

### TiLite Aero T

<b>Frame features</b>	Aluminum Dual-Tube Rigid, TiFit, full adjustability	<b>Folding seat back</b>	Adj. height 9.5" – 22" Adj. angle 80° – 101°
<b>Total chair width*</b>	23"	<b>Front seat width (seat taper)</b>	0", 1", 2", 3"
<b>Seat width</b>	12" – 20"	<b>Ergonomic seat</b>	6", 7", 8"
<b>Seat depth</b>	12" – 20"	<b>Transit options</b>	Q-Strait lock down approved
<b>Front seat height</b>	15" – 21"	<b>Transport weight*</b>	12 lb.
<b>Rear seat height</b>	13" – 21"	<b>Weight limit</b>	265 lb Standard 350 lb Heavy duty
<b>Front frame angle</b>	70°, 75°, 80°, 85°, or 90°	<b>HCPCS Code</b>	K0005, E1235
<b>Camber</b>	0°, 2°, 4°, or 6°	<b>Warranty</b>	Five years on frame One year on TiLite components
<b>Center of gravity adjustment</b>	5"		

### TiLite Aero Z

<b>Frame features</b>	Aluminum Mono-Tube Rigid, TiFit, full adjustability	<b>Folding seat back</b>	Adj. height 9.5" – 22" Adj. angle 80° – 101°
<b>Total chair width*</b>	23"	<b>Front seat width (seat taper)</b>	0", 1.5", 2.5"
<b>Seat width</b>	10" – 20"	<b>Ergonomic seat</b>	6", 7", 8"
<b>Seat depth</b>	10" – 20"	<b>Transit options</b>	Transit Tie-Down
<b>Front seat height</b>	16" – 21"	<b>Transport weight*</b>	12.4 lb
<b>Rear seat height</b>	12.5" – 20"	<b>Weight limit</b>	265 lb
<b>Front frame angle</b>	70°, 75°, 80°, 85°, or 90°	<b>HCPCS Code</b>	K0005, E1235
<b>Camber</b>	0°, 2°, 4°, or 6°	<b>Warranty</b>	Five years on frame One year on TiLite components
<b>Center of gravity adjustment</b>	6"		

\*Weights and widths based on 16" x 16" chair with 0° camber, minimum rear wheel spacing and all standard features. Rear wheels are not included in transport weight. For reference, a pair of standard 24" rear wheels weighs 7.2 lb.



Aluminum adjustable wheelchairs

TiLite® Aero T | TiLite® Aero Z



# TiLite Aero T

Our most advanced aluminum chair, the Aero T's dual-tube frame minimizes weight while maximizing performance.

Custom designed and built to your specifications, with available ergo seating

Find your perfect ride balance with our adjustable center of gravity

Adjustable seating angle to maximize comfort

Dual-tube geometry provides better load distribution and less frame flex

Adjustable caster fork angle to ensure a precise ride



# TiLite Aero Z

Distinguished by a clean, minimal look, the mono-tube frame's 1.25" tubing gives the Aero Z strength and durability.

Custom designed and built to your specifications, ergo seating available

Adjustable center of gravity to find that perfect ride balance

Adjustable seating angle to maximize comfort

1.25" diameter tubing to maximize performance and minimize visual impact

Adjustable caster fork angle to ensure a precise ride

