


Tattoo options



Frame colors

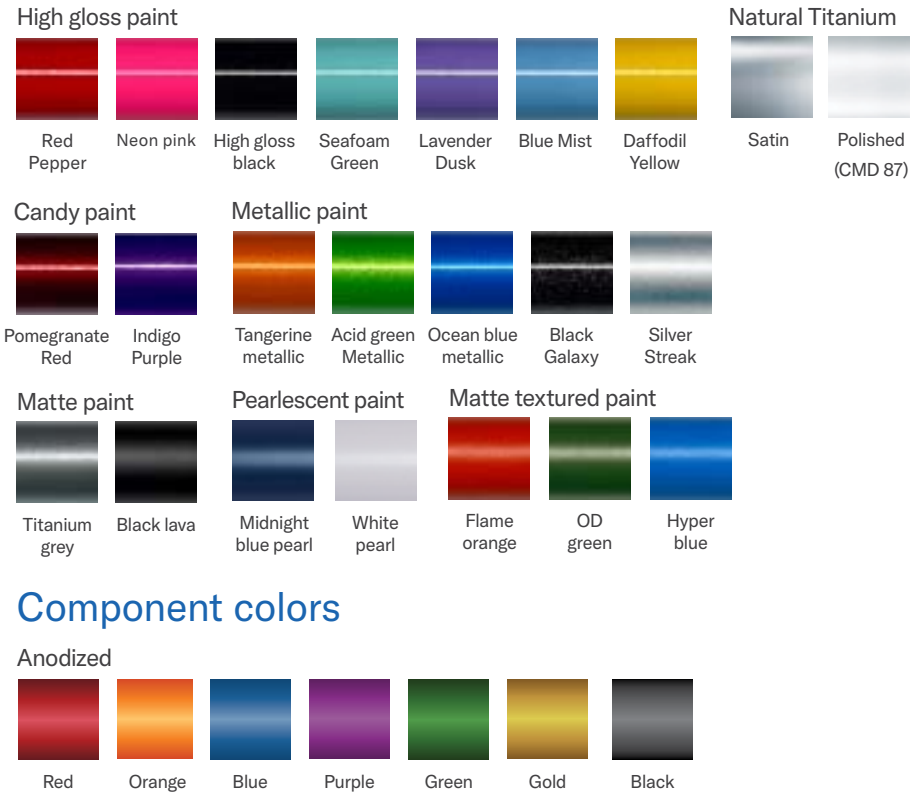

Specifications

TiLite ZR

Frame features	Titanium Mono-Tube Rigid, TiFit, Limited adjustability	Center of gravity adjustment	1.6"
Total chair width*	23"	Folding seat back	Adj. height 9.5" – 22" Adj. angle 80° – 101°
Seat width	10" – 20"	Fixed seat back	Fixed height 8" – 20" Fixed angle 80° – 101°
Front seat width (seat taper)	0", 1.5", 2.5"	Transit options	Transport bracket
Seat depth Ergonomic seat	10" – 20" 5" – 9"	Transport weight*	10.9 lb
Front seat height	16" – 21"	Weight limit	265 lb Standard 300 lb Heavy duty
Rear seat height	12.5" – 20"	Warranty	Five years on frame One year on TiLite components
Front frame angle	70°, 75°, 80°, 85°, or 90°		
Camber	0°, 2°, 4°, or 6°		

TiLite ZRA

Frame features	Titanium Mono-Tube Rigid, TiFit, Full adjustability	Front seat width (seat taper)	0", 1.5", 2.5"
Total chair width*	23"	Folding seat back	Adj. height 9.5" – 22" Adj. angle 80° – 101°
Seat width	10" – 20"	Ergonomic seat	6", 7", 8"
Seat depth	10" – 20"	Transit options	Transit tie-down
Front seat height	16" – 21"	Transport weight*	11.2 lb
Rear seat height	12.5" – 20"	Weight limit	265 lb Standard 300 lb Heavy duty
Front frame angle	70°, 75°, 80°, 85°, or 90°	Warranty	Five years on frame One year on TiLite components
Camber	0°, 2°, 4°, 6°, 8° or 12°		
Center of gravity adjustment	6"		

*Weights and widths based on 16" x 16" chair with 0° camber, wheel locks, folding backrest, minimum rear wheel spacing and all standard features. Rear wheels are not included in transport weight. For reference, a pair of standard 24" rear wheels weighs 7.2 lb.



Titanium rigid wheelchairs

TiLite® ZR | TiLite® ZRA



TiLite ZR

A fixed-frame mono-tube wheelchair for individuals who want a solid, streamlined ride with a minimalist aesthetic. Configured with TiFit, the TiLite ZR offers the highest level of configuration while keeping moving pieces to only those that are truly necessary for experienced individuals.

Custom designed and built to your specifications

Fixed or ergo seating position is custom fabricated to the user's specifications

Adjustable center of gravity allows the individual to fine-tune the balance and feel of their chair.



Straight caster arms and a redesigned caster bullet that transitions to a titanium fork stem saves weight and improves the chair's responsiveness

TiLite ZRA

A fixed-frame mono-tube wheelchair for individuals who want a TiFit configuration with unmatched strength and weight while leaving some room for adjustments, allowing for a wheelchair that adapts to their changing seating and mobility needs.

Custom designed and built to your specifications, ergo seating available

Adjustable seating angle to maximize comfort

Adjustable center of gravity to find your perfect ride balance

1.25" diameter tubing to maximize performance and minimize visual impact

Adjustable caster fork angle to ensure a precision ride

