

ActiveHeight & ActiveReach packages



	M3	M5	MVS	F3	F5	F5VS
↑ ActiveHeight	12"	14"	12"	12"	14"	14"
↶ ActiveReach	5°, 10° & 20°	5°, 10° & 20°	Standing	5°, 10°, 20° & 30°	5°, 10°, 20° & 45°	Standing

**Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions.*



Permobil New Zealand | 0800 115 222 | sales.nz@permobil.com | www.permobil.com/en-nz

BR-NZ-PCS-001.0


Power wheelchairs

Corpus® Seating System



Created for you

Corpus is the most ergonomic and capable wheelchair seating system in the world. It is elegantly designed and configured to empower functional movements while supporting your optimal body position.

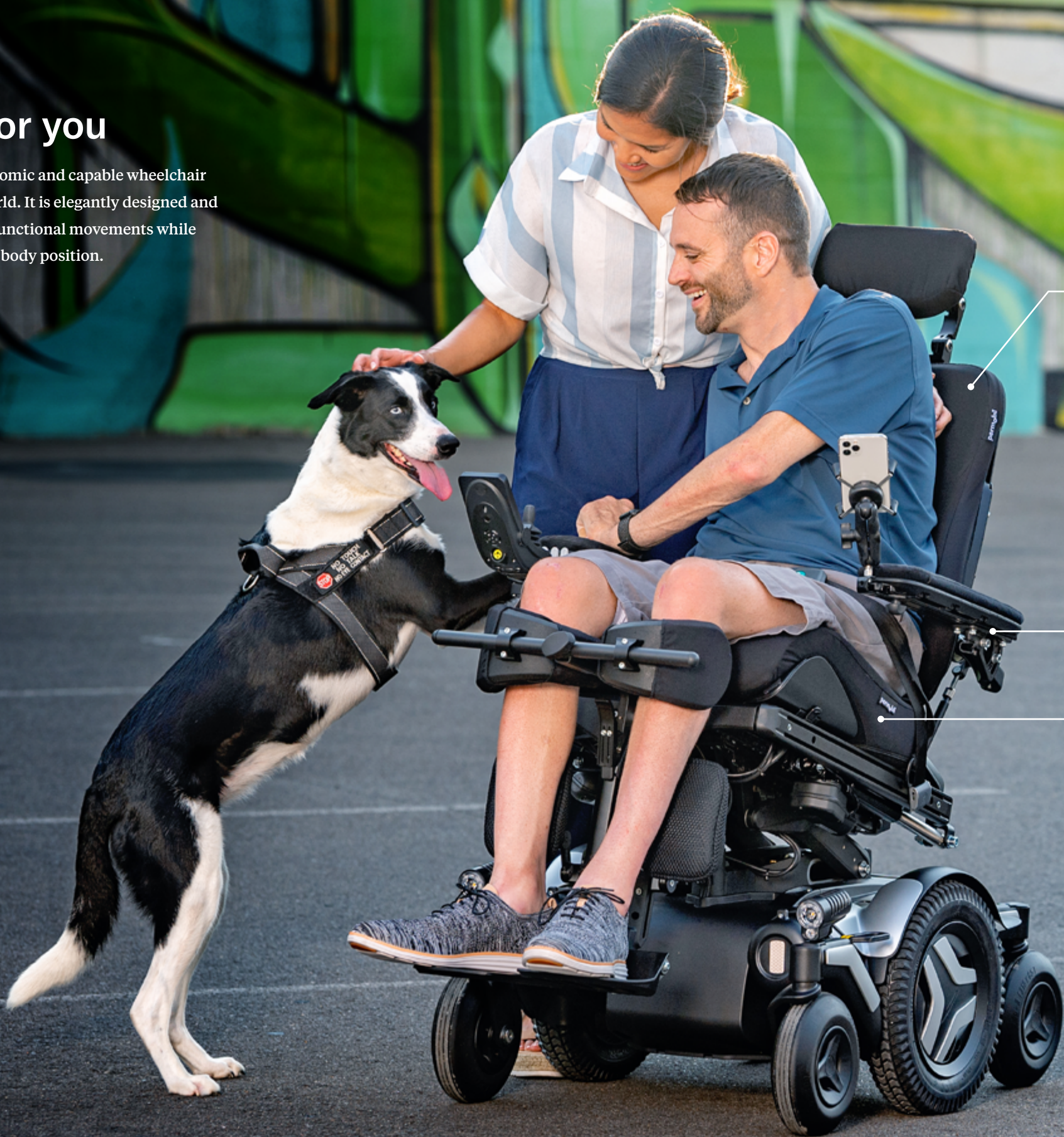

Dual density back support foam for greater comfort and surface area contact



Leatherette



STRETCH-AIR®

Multi-axis armrests allow independent height, angle and optional medial swivel adjustment

Redesigned pelvic well insert for improved pressure redistribution



Independent repositioning

Manage your seating goals with multiple memory seat positions and Independent Repositioning Mode (IRM). IRM enables automatic sequencing of tilt, leg elevation and recline promoting maximum benefit of pressure relief protocols without loss of positioning.

Access Independent Repositioning Mode at the touch of a button

Up to **180°** power legrest elevation

Seamless shear reduction and optimised pivot points keeps you in position

Up to **180°** power recline and **50°** posterior tilt

F5



ActiveHeight

Dynamic seat elevation

ActiveHeight raises your seat up to 14". It is unique in that while elevating, the seating system repositions 3.5" rearward to maintain forward stability and create a more compact footprint enabling you to get closer to the important things in your life. ActiveHeight empowers you to have face-to-face conversations, order food at a high counter, and accomplish many other daily tasks with ease.

ActiveHeight + ActiveReach

Seat elevation and forward tilt combine seamlessly to give you the ultimate advantage in accessing your environment.



ActiveReach

Functional forward tilt

Face daily tasks head on with ActiveReach. By combining Corpus seat functions with up to 45 degrees of forward tilt, you can achieve improved ergonomics for overhead reaching and table top activities. ActiveReach can also help improve posture for social engagement and body mechanics during transfers.



10° of no charge ActiveReach with ActiveHeight



Adjustable
chest support

Swing-away or
single-post knee
block options

Articulating power
footplate

Corpus VS

Superior power standing

Fully programmable standing sequences allow user position to be optimised using supine, semi-reclined, or sit-to-stand techniques. Further customise standing positions with adjustable chest support, multiple knee support options, and power articulating footplates. From face-to-face interactions to improved health outcomes, the Corpus VS seating system empowers you to stand with confidence.



Stand while
driving



Customisable
standing sequences



Full range power
functions